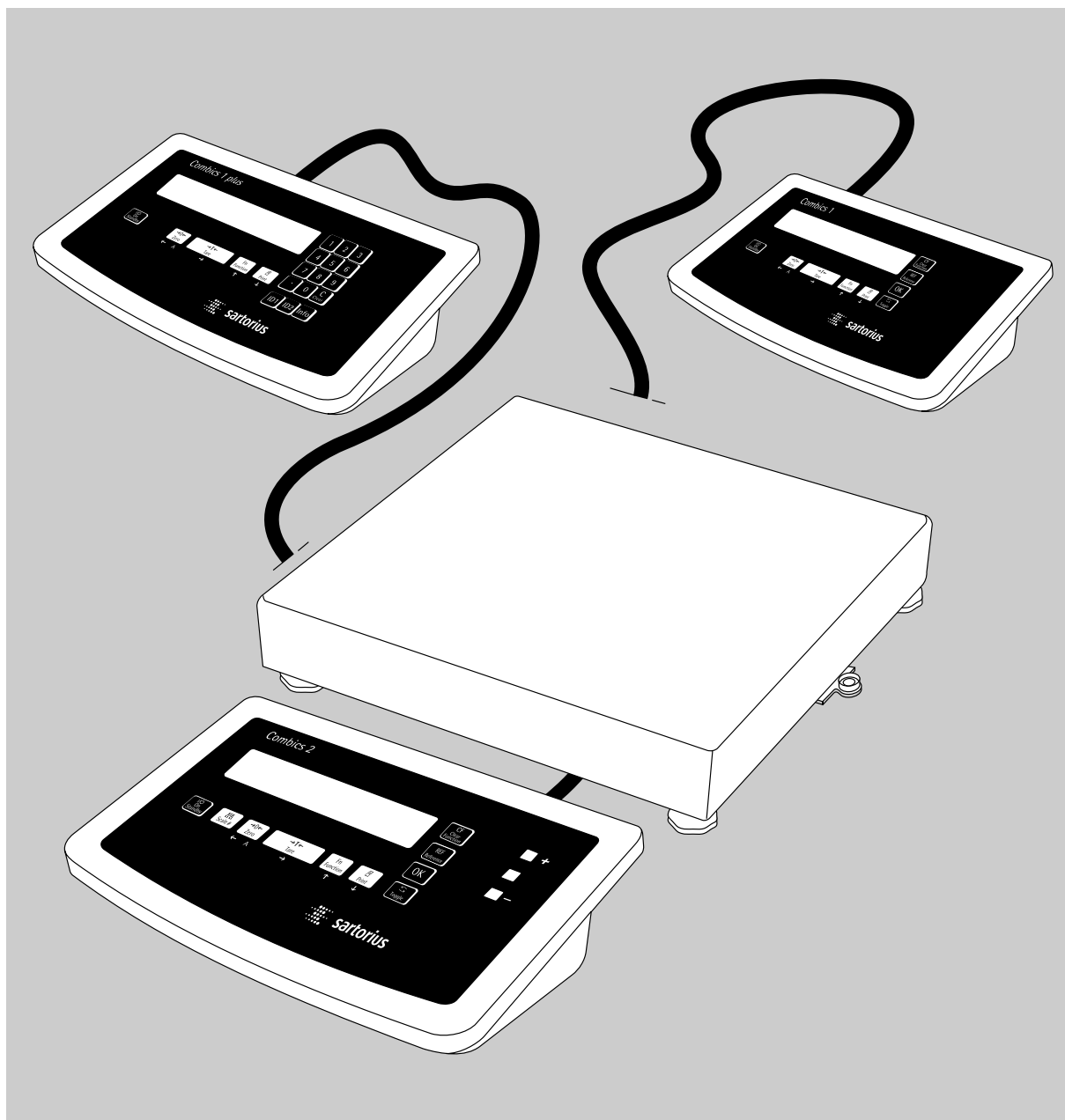


Operating Instructions

Sartorius Combics 1 | Combics 1 plus | Combics 2

Models CW1P | CW1NP | CW2P | CW1S | CW1NS | CW2S

Combics Complete Scales



Intended Use

Combics 1, Combics 1 plus and Combics 2 are rugged display and control units for the complex quality control tasks you perform every day. They meet the highest requirements placed on the accuracy and reliability of weighing results:

- in the food industry
- in the pharmaceutical industry
- in the chemical industry
- in the electronics and metal-working industries

Combics complete scales are:

- Rugged and durable (stainless steel housing)
- Easy to clean and disinfect
- Easy to operate, thanks to the following features:
 - large, backlit display segments
 - large keys with positive click action
- Independent of the weighing instrument location
- Equipped with a range of interfaces for flexible use
- Password-protected from unauthorized changes in parameters (optional)

Combics 1 plus speeds up your routine procedures with:

- Input functions for tare values through numeric keypad
- Option for 2 alphanumeric lines to identify samples
- Connectivity for bar code scanner to enter tare values or ID codes

Combics 2 indicators have the following features:

- Built-in application programs:
 - Counting
 - Neutral measurement
 - Weighing in percent
 - Averaging
 - Checkweighing
 - Classification
 - Net-total formulation
 - Totalizing
- Automatic initialization when you switch on the Combics
- Automatic taring when a load is placed on the weighing instrument
- Optional remote control using an external computer

Symbols

The following symbols are used in these instructions:

- indicates required steps
- describes what happens after you have performed a certain step

⚠ indicates a hazard

Hotline:

For advice on the use of these applications, just call or fax your local Sartorius office. For the address, please visit our Internet website at: www.sartorius.com

Contents

2	Intended Use	77	Data Interfaces
4	Warnings and Safety Precautions	79	Pin Assignment Charts
5	Getting Started	81	Installing the Interface Cable
6	General View of the Equipment	82	Cabling Diagram
10	Operating Design	83	Synchronization
10	Weighing/Measurement	84	Configuring the Data Interface as a COM Port
12	Configuration (Operating Menu)	84	Data Input Format
13	Operation	85	Data Output Format
13	Weighing $\Delta\bar{\Delta}$	87	Configuring the Data Interface as a Printer Port
21	Individual ID Codes	87	Configuring Printouts
23	Calibration and Adjustment	89	Sample Printouts
25	Counting \clubsuit	92	Error Codes
29	Neutral Measurement $\clubsuit \cap \bar{\cap}$	93	Care and Maintenance
33	Checkweighing $\frac{\pi}{L}$	93	Repairs
38	Classification $\frac{\pi}{L} \bar{L}$	93	Cleaning
43	Weighing in Percent %	93	Safety Inspection
47	Averaging (Animal Weighing) $\frac{\pi}{L}$	93	Recycling
50	Net-total Formulation \bar{L}	94	Overview
53	Totalizing Σ	94	Common Specifications
57	Configuration	94	Model-specific Specifications
57	Operating Menu Overview	95	Type Designation
58	Setting the Language (Example)	97	Dimensions (Scale Drawings)
59	Entering/Changing the Password (Example)	98	Accessories
61	Operating Menu Overview (Parameters)	103	Declaration of Conformity
		105	EC Type-approval Certificate
		107	Plates and Markings
		109	Index
			Appendix:
			General Password

Warning and Safety Precautions

Safety Information:

- Please read these operating instructions carefully before using your scale to prevent damage to the equipment.
- ⚠ Do not use this equipment in hazardous areas.
- ⚠ Use only standard cables that have protective grounding conductors. The protective conductor must not be disconnected for any reason.
- ⚠ Disconnect the scale from power before connecting or disconnecting peripheral devices.
- ⚠ The scale may be opened only by trained service technicians.
- ⚠ If you operate the equipment under ambient conditions subject to higher safety standards, you must comply with the applicable installation regulations.
- ⚠ If there is visible damage to the equipment or power cord, unplug the equipment and make sure it cannot be used for the time being.
- ⚠ If you use electrical equipment in installations and under ambient conditions requiring higher safety standards, you must comply with the provisions as specified in the applicable regulations for installation in your country.

Installation:

- Proceed with extreme caution when using pre-wired RS-232 connecting cables from other manufacturers, as the pin assignments may not be compatible with Sartorius equipment. Check all pin assignments against the cabling diagrams and disconnect any lines that are not assigned.
- ⚠ Always wear gloves, safety boots and protective clothing when lifting the load plate with a vacuum lifting pad. Danger of injury! This work may be carried out only by authorized and properly trained personnel.
- Weighing platforms with dimensions larger than 1 x 1 m are provided with suspension supports. Be careful not to stand under the load when the weighing platform/load plate is being transported or lifted with a crane. Always comply with the applicable accident prevention regulations. Make sure to avoid damaging the terminal box and housing or the load cell during transport.
- Connect only Sartorius accessories and options, as these are optimally designed for use with your Combics indicator. The operator shall be solely responsible for installation and testing of any modifications to Sartorius equipment, including connection of cables or equipment not supplied by Sartorius. Contact Sartorius for detailed operating specifications in accordance with the Standards for immunity to interference.
- Do not expose the indicator to aggressive chemical vapors or to extreme temperatures, moisture, shocks, or vibration.
- Clean your Combics only in accordance with the cleaning instructions (see “Care and Maintenance”).
- If you have any problems with your Combics indicator, contact your local Sartorius customer service center.

IP Rating:

- CISL models are rated to IP44 (with Option L1: IP65); CIS models are rated to IP67.
- The IP65/IP67 protection rating is ensured only if the rubber gasket is installed and all connections are fastened securely (including the caps on unused sockets). Weighing instruments must be installed and tested by a certified technician.
- If you install an interface port after setting up your indicator, keep the protective cap in a safe place for future use. The cap protects the interface connector from vapors, moisture and dust or dirt.

Using the Equipment in Legal Metrology in the EU*:

- When the indicator is connected to a weighing instrument and the resulting weighing instrument is to be verified, make sure to observe the applicable regulations regarding verification. When connecting a Sartorius weighing instrument, please observe the “Guide to Verification” on the enclosed CD and the permitted weighing range as listed in the Declaration of Conformity.
- EU legislation requires that a control seal be affixed to the verified device. The control seal consists of a sticker with the “Sartorius” logo. This seal will be irreparably damaged if you attempt to remove it. If any of the verification seals are damaged, make sure to observe the national regulations and standards applicable in your country in such cases. In some countries, the verification will become null and void and the equipment must be re-verified.

* Including the Signatories of the Agreement on the European Economic Area

Getting Started

The complete scale is available in various versions. If you have ordered special options, the indicator will be supplied with these options premounted at the factory.

Storage and Shipping Conditions

- Allowable storage temperature: -10 ... +40°C (+14°F ... + 104°F)
- Unpackaged equipment may lose its accuracy when exposed to strong vibration. Excessive vibration may compromise the safety of the equipment.
- Do not expose the indicator unnecessarily to extreme temperatures, moisture, shocks, or vibration.

Unpacking

- After unpacking the equipment, please check it immediately for any visible damage.
- If you detect any damage, proceed as directed in the chapter entitled “Care and Maintenance” under “Safety Inspection.”
- If you will need to ship the equipment later, save all parts of the packaging because only the original packaging provides the best protection for shipment.
- Before shipping, be sure to disconnect all cables to prevent damage.

Equipment Supplied

- Indicator
- Weighing platform
- Operating instructions (this manual)
- Special accessories listed on the bill of delivery, if ordered

Installation Instructions

Choose a location that is not subject to the following unfavorable conditions:

- Excessive temperatures (operating temperature range: -10°C to +40°C; +14°F to +104°F)
- Aggressive chemical vapors
- Excessive moisture (depends on IP rating)

Conditioning the Indicator

Moisture in the air can condense on cold surfaces whenever the equipment is moved to a substantially warmer place. To avoid the effects of condensation, condition the indicator for about 2 hours at room temperature, leaving it unplugged from AC power.

Checking the Geographical Data Entered for Use in Legal Metrology

Preparation (see also the “Device Information” menu items listed under “Operating Menu Overview” in the chapter entitled “Settings”).

- Press $\boxed{I/O}$ to turn on the CombiCS
- While all segments are lit, press $\boxed{\rightarrow T \leftarrow}$
 - > *RPPL* is displayed
- Select “Info”: Press \boxed{Fn} repeatedly; press $\boxed{\rightarrow T \leftarrow}$ to confirm
- Select “Device Specific Information”:
 - for WP1 or WP2: Press \boxed{Fn} repeatedly; press $\boxed{\rightarrow T \leftarrow}$ to confirm
 - > Press \boxed{Fn} repeatedly to scroll through the geographical data (data shown depends on input before verification), for example:
 - Latitude (in degrees): 51°4
 - Elevation (in meters): 513⁵
 - or
 - Acceleration of gravity in m/s⁻²: 9.810⁶

The scale can be used anywhere in Germany if the geographical data is as follows:

- Latitude: 51.00 degrees
- Elevation: 513 meters
- This data corresponds to the following value:
 - Acceleration of gravity: 9.810 m/s⁻²

These values are calculated for Germany based on a mean value for the Earth’s acceleration. The greater the precision of the geographical data entered, the greater the precision achieved with the weighing instrument; the tolerance range, however, is also restricted accordingly.

- The tolerance ranges, for example for a scale with 3000 e, are as follows:
- ±100 km for the latitude and
 - ±200 m for the elevation above sea level.

⚠ If used outside the specified zone, the scale must be re-verified for use in legal metrology. Please contact an authorized service technician.


Seal on Indicators Verified for Use in Legal Metrology in the EU*

EU legislation requires that a control seal be affixed to the verified device. The control seal consists of a sticker with the “Sartorius” logo. This seal will be irreparably damaged if you attempt to remove it. If the seal is broken, the validity of the version becomes null and void, and you must have your scale re-verified.

* Including the Signatories of the Agreement on the European Economic Area

General View of the Equipment

Combics 1 | Combics 1 plus

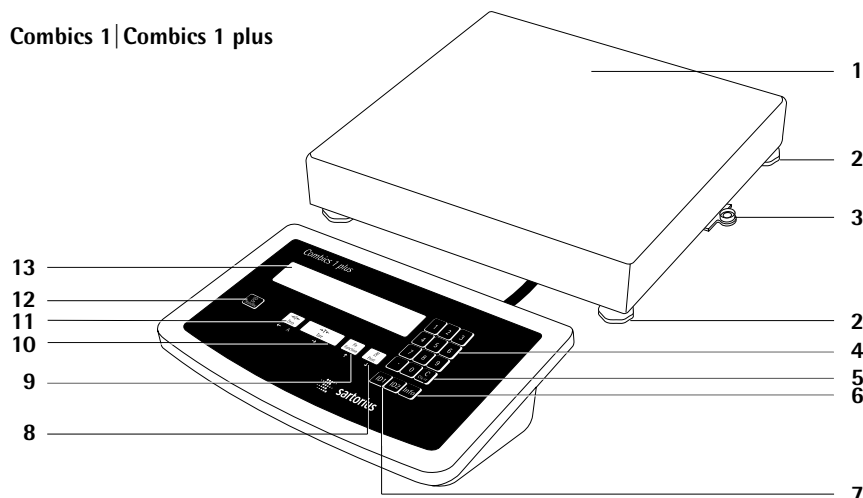


Display and Keypad:


- 1 Load plate
- 2 Leveling feet
- 3 Level indicator
- 4* Keypad for numeric input
- 5* "Clear" key (deletes ID codes and tare input)
- 6* "Info" key (shows ID codes and tare input)
- 7* ID keys (for entering ID codes)
- 8 Data output
- 9 Gross/net; 2nd unit or 10 x higher resolution (depending on the settings)
- 10 Tare
- 11 Zero
- 12 On/off key
- 13 Display (for a detailed view, see chapter entitled "Operating Design")
- 14 Select reference weight (depending on the application)
- 15 "Clear" key
- 16 LEDs (for checkweighing and classification)
- 17 Start application
- 18 Toggle to the application program | application-specific information
- 19 Toggle to different weighing platform

* Combics 1 plus only

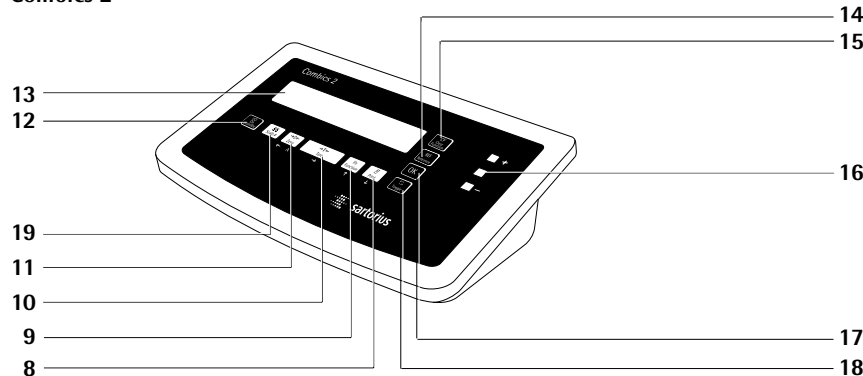Combics 2




Rear view: Model: CW1P | CW1NP | CW2P

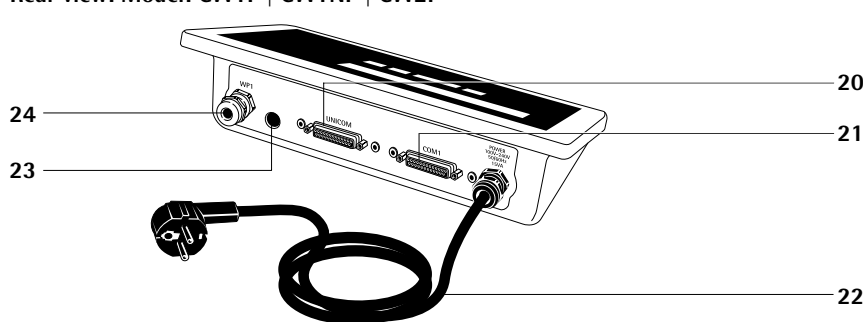


Rear View

- 20 Second „UNICOM“ interface port for bar code scanner or external rechargeable battery pack (Combics 2 only) (additional functions optional)
- 21 RS-232C "COM1" interface port (standard)
- 22 Power cord with country-specific plug
- 23 Menu access switch (standard operating mode or legal metrology mode)
- 24 Connector for weighing platform
- 25 Vent valve

Rear view: Model: CW1S | CW1NS | CW2S

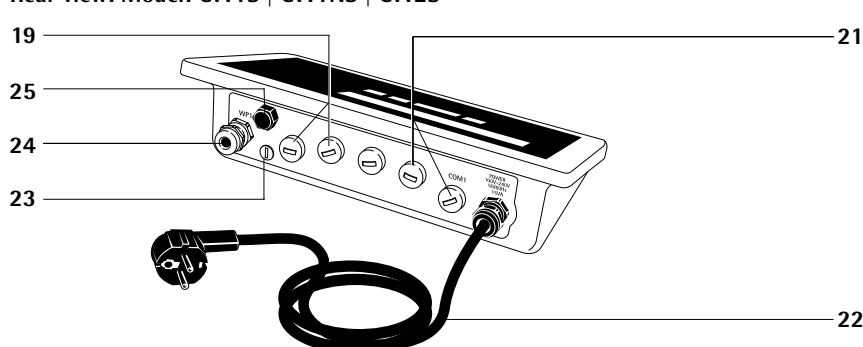

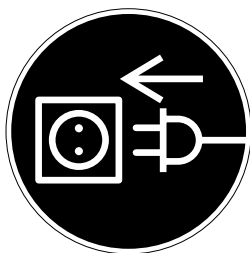

Connecting the Equipment to AC Power

- Check the voltage rating and the plug design.
- The scale is powered through the pre-installed power cord. The power supply is built into the scale, which can be operated with a supply voltage of 100V to 240V. Make sure that the voltage rating printed on the manufacturer's ID label is identical to that of your local line voltage. If the voltage specified on the label or the plug design of the AC adapter does not match the rating or standard you use, please contact your Sartorius office or dealer.
The power connection must be made in accordance with the regulations applicable in your country.
- To power a protective class 1 device, plug the power cord into an electrical outlet (mains supply) that is properly installed with a protective grounding conductor (protective earth = PE).

Safety Precautions

If your local AC output does not have a protective grounding conductor (protective earth), have a certified electrician install equivalent protection according to your country's valid installation requirements. Make sure the protective grounding effect is not neutralized by use of an extension cord that lacks a protective grounding conductor.

Connecting Electronic Peripheral Devices

- Make absolutely sure to unplug the device from AC power before you connect or disconnect a peripheral device (printer or PC) to or from the interface port.

Warmup Time

To deliver exact results, the device must warm up for at least 30 minutes after initial connection to AC power or after a relatively long power outage. Only after this time will the device have reached the required operating temperature.

Using Equipment Verified as Legal Measuring Instruments in the EU*:

- Make sure to allow the equipment to warm up for at least 24 hours after initial connection to AC power or after a relatively long power outage.

Connecting the External Rechargeable Battery Pack (Accessory: YRB10Z)

- △ Disconnect the equipment from AC power (unplug the AC adapter)

● Installation


For model CW1P | CW1NP | CW2P: Connect a 25-pin D-SUB male connector (connecting cable YCC02-RB01) to a second "UNICOM" interface port

For model CW1S | CW1NS | CW2S: see the Section "Pin Assignment Chart" (via connecting cable YCC02-RB02 or as option L2)

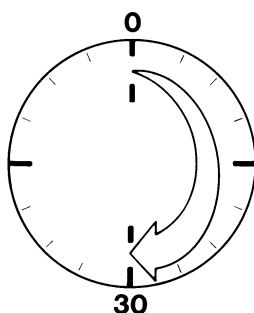
Operating

- Hours of operation: up to 40, depending on the weighing platform connected; without options. The Combics will automatically switch to battery operation whenever there is a power shortage or the power is cut off. Once the mains power supply is reinstated, the Combics will automatically switch back to normal operation.

Battery symbol

Battery fully charged: 

Battery empty: 



* including the Signatories of the Agreement on the European Economic Area

Connecting a Bar Code Scanner (Accessory: YBR02CISL)

⚠ Disconnect the equipment from AC power (unplug the AC adapter)

○ Installation


For model CW1NP | CW2P:

- Connect a 25-pin D-SUB male connector (connecting cable YCC02-BR01) to a second “UNICOM” interface port
- For bar code scanner and external rechargeable battery: please use T connector YTC01.


For model CW1NS | CW2S: see the section entitled “Pin Assignment Chart” (via connecting cable YCC02-RB02 or as option M8)


Installing the Verification Adapter for Use in Legal Metrology (on verifiable models only)

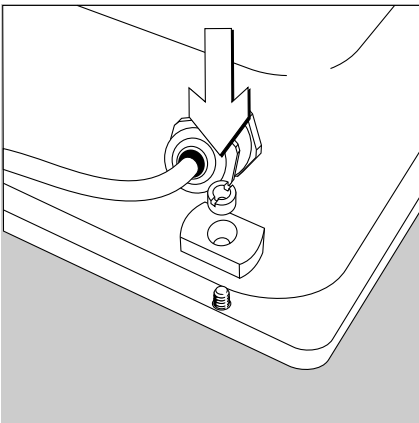- Remove the nut located on the back of the indicator
- Use the slotted screw to install the adapter plate



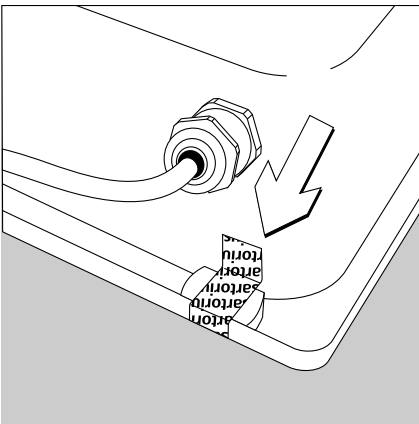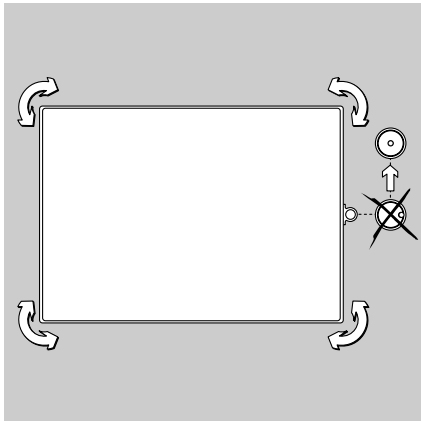- Affix the verification seal over the adapter






Leveling the Weighing Platform


Purpose:

- To compensate for uneven areas at the place of installation
 - Make sure that the equipment is placed in a perfectly horizontal position for consistently reproducible weighing results
- Always re-level the weighing platform after changing the place of installation.

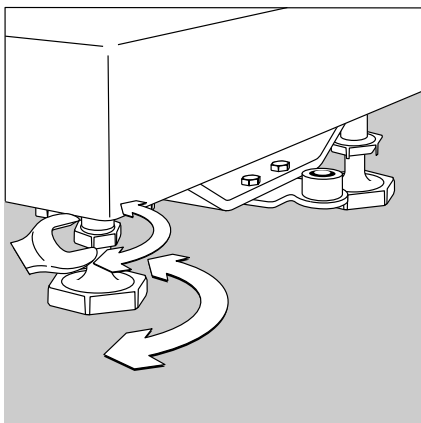- Level the weighing platform using the four adjustable feet. Keep turning the feet until the air bubble is centered in the level indicator circle.



- Make sure that all four feet are touching the ground.
- > The load must be equally distributed over all leveling feet!
- Loosen the lock nuts on the leveling feet with an open-end wrench.
- > Setting the adjustable feet:
Raise the weighing platform by extending the leveling feet (turning to the right).
Lower the weighing platform by retracting the leveling feet (turning to the left).
- After aligning the weighing platform, tighten the lock nuts.
Small platforms (1 weighing cell): against the platform frame,
Large platforms (4 weighing cells): against the platform foot.

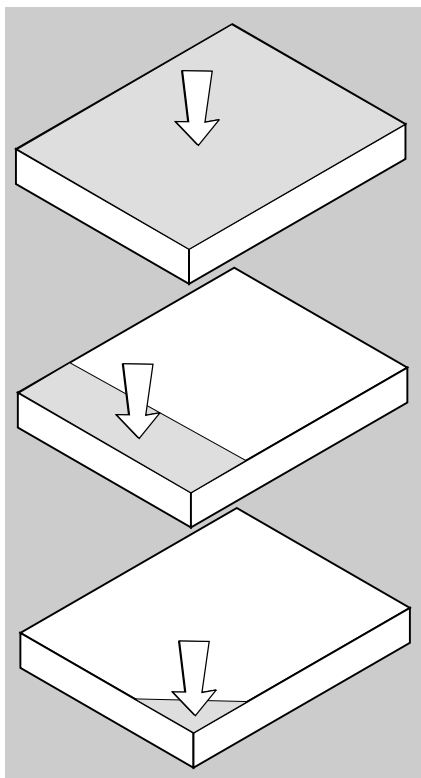

Operating Tolerances

Never exceed the maximum capacity of the weighing platforms.

The maximum loading capacities of the weighing platforms are listed in the table below and depend on the position of the weight loaded on the platform (center, sides, one-sided corner load):

Platform dimensions	Center	Side	Corner
320 x 240	50	35	20
400 x 300	130	85	45
500 x 400	300	200	100
500 x 400 (P*)	600	400	200
650 x 500 (S**)	450	300	150
800 x 600 (P*)	1200	800	400
800 x 600 (S**)	900	600	300
1000 x 1000	4500	3000	1500
1250 x 1000	4500	3000	1500
1500 x 1250	4500	3000	1500
1500 x 1500	4500	3000	1500
2000 x 1500	4500	3000	1500

* Steel

** Stainless steel

Keys

The operation of Combics 1, Combics 1 plus and Combics 2 involves just a few keys. These keys have one function during measurement and another during configuration. Some of the keys have one function when pressed briefly, and another when held for longer than 2 seconds.

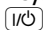
If a key is inactive, this is indicated as follows when it is pressed:


- The error code "-----" is displayed for 2 seconds. The display then returns to the previous screen content.
- An acoustic signal (double-beep) is emitted.


Weighing/Measurement


Input Through the Keypad


Keys below the Display

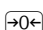
 On/off key
(in standby mode, αFF is displayed).

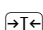 - Press briefly:
Zero the instrument
- Press and hold (> 2 sec):
Show the adjustment/
configuration counter

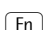 - Tare the instrument
- Save numeric input as tare weight (Combics 1 plus only)
- Press and hold (> 2 sec):
Start calibration/adjustment







 Toggle the display between (depends on configuration):
- first and second weight unit,
- gross and net values, or
- normal and 10-fold increased display resolution

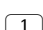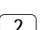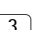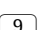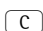

 - Press briefly: Print
- Press and hold (> 2 sec):
Print GMP footer


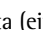
 When two weighing instruments are connected, this key toggles the display between instruments (Combics 2 only).



Keys to the Right of the Display on Combics 1 plus

Keys for entering ID codes and tare values

, , ... , 
Numeric keypad for entering values that are identified by the key subsequently pressed (e.g.,  for tare input or "ID" key for ID codes)


 Press to delete data (either ID code or tare input, depending on subsequent key; e.g.,  for tare)
When pressed during numeric input: deletes the last character entered

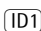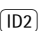
 Press to enter data (either ID code or tare value, depending on subsequent key; e.g.,  for tare)


,  Press to store or view ID codes (user-defined data to identify weight values)

Keys to the Right of the Display on Combics 2

These keys are used for operating applications. Please refer to the individual application descriptions for details.

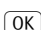 Deletes initialization values or totalizing memory, depending on configuration.

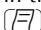

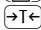
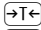
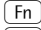
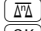
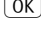
 For modifying reference values.

 Store a value or start an application program.

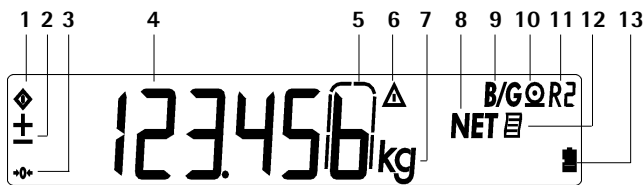
 Toggle between display modes within an application program

Input Through the Digital Input Port


The indicator is equipped with a control input (universal input port). You can connect a hand switch or foot switch to this port, if desired. Assign one of the following functions to this port in the operating menu:

-  key
-  (> 2 sec.)
-  key
-  key (> 2 sec.)
-  key
-  key
-  key

Operating Design



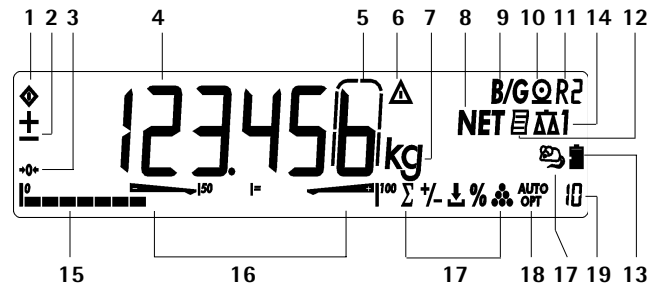Display during weighing/measurement
(Combics 1, Combics 1 plus) (example)



Display during weighing/measurement
(Combics 2) (example)

Display in Weighing Mode


The illustration above shows all of the main display elements and symbols shown during weighing.

- 1 Busy symbol; indicates that an internal process is in progress
- 2 ± sign for the value displayed
- 3 Identifies “zero” as a weight value (after the weighing instrument has been zeroed)
- 4 Weight value or calculated value (main display)
- 5 In legal metrology, on equipment with e = d, the digit shown with a border is not valid
- 6 Identification of calculated value in the main display (value not valid in legal metrology)
- 7 Weight unit of the value displayed
- 8 Net value in the main display (when data is stored in tare memory)
- 9 Gross value in main display (when data is stored in tare memory)
- 10 Printing in progress
- 11 Display of the range on multiple-range instruments
- 12 GMP-compliant printing in progress (Combics 1 plus and Combics 2 only)
- 13 Battery symbol showing status of rechargeable battery (empty outline indicates battery is drained)


- 14 Indicates active weighing instrument; flashes to prompt calibration/adjustment (Combics 2 only)
- 15 Bar graph (Combics 2 only)
 - Shows how much of the available weighing capacity is “used up” by the current load, or
 - Shows the measured value in relation to a target value (with the Checkweighing or Classification application)
- 16 Symbols for Checkweighing and classification (Combics 2 only)
- 17 Application symbols: Totalizing, Checkweighing, Classification, Net-total Formulation, Weighing in Percent, Counting and Neutral Measurement. For details, please refer to the the descriptions of the applications (Combics 2 only).
- 18 Symbols for reference updating (Combics 2 only)
 - Auto: Depending on the weight value, a reaction is triggered in the application
 - Opt: Automatic reference updating has been performed (Counting application)
- 19 Numeric display; e.g., shows reference value (Combics 2 only)

Saving Data in Weighing Mode

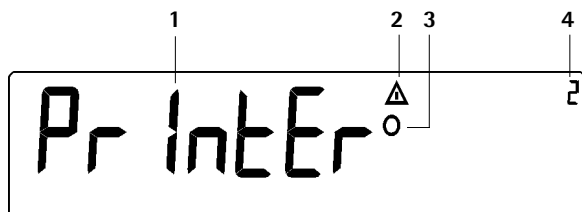- All of the application parameters saved (e.g., reference values) remain stored and are still available when
- you switch the Combics off and then back on again
 - you return to the originally selected application from a second one (for example, when you switch from Averaging back to Counting, all parameters saved for Counting are available)



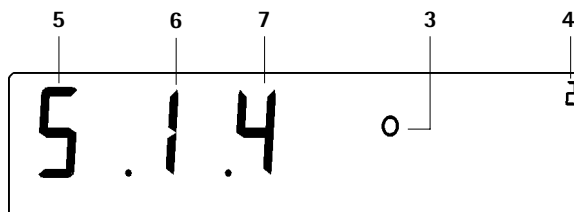Display of menu settings: Text menu (example)



Display of menu settings: Numeric menu (example)

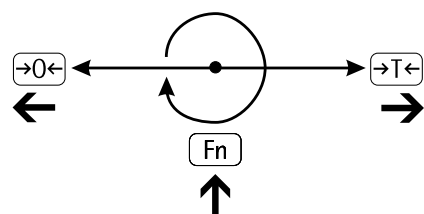
Configuration (Operating Menu)

Navigation and input in the operating menu are implemented using the keys below the display (on Combics 1 plus, numeric values can be entered using the numeric keypad)

Opening the Menu

Press the key to switch the Combics off and then on again; while all segments are displayed, press the key briefly.

Navigating the Menu




Close the active submenu and return to next higher menu level ("back")

– Press briefly (< 2 sec):
Select and store a menu item
– Press and hold (> 2 sec):
Exit the menu

Show the next item on the same menu level (the display scrolls through all items in series)

Print the menu settings starting from the current position, or print Info data

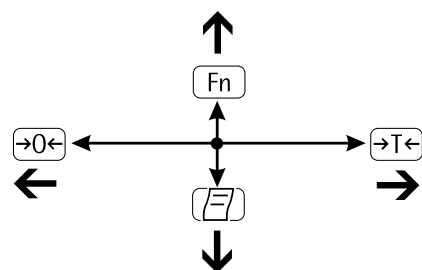Alphanumeric Input in the Menu



– Press briefly:
Activate character to the left of the current character (when first character is active: exit input mode without saving changes)
– Press and hold (> 2 sec):
Exit the input mode without saving changes

– Press briefly (< 2 sec):
Confirm currently active character and move 1 position to the right (after the last character: store input)
– Press and hold (> 2 sec):
Store current input and display the menu item

– Cursor in first position, no characters entered yet: Delete character(s) and enter 0
– Change the displayed character; scroll forward (sequence: 0 to 9, decimal point, minus sign, Z to A, space)

– Cursor in first position, no characters entered yet: Delete entire string and enter a space
– Change the displayed character; scroll backwards (sequence: Space, A to Z, minus sign, decimal point, 9 to 0)

Numeric input in Combics 1 plus operating menu:

Enter values (date and time, etc.) using the 10-key numeric keypad

Display of Menu Settings

The illustration above shows all of the main display elements and symbols shown during menu configuration.

- 1 Selected menu item on text level (e.g., "Printer" for configuring the connected printer)
- 2 Indication that there are additional submenus
- 3 Indication that this is the currently active setting
- 4 Menu history (indicates the highest menu level)
- 5 Highest level in numeric menu
- 6 Second level in numeric menu
- 7 Third level in numeric menu

Saving Data in Configuration Mode

The parameters selected in the operating menu remain stored after you switch off the Combics.

You can prevent unauthorized changes in operating menu settings by requiring password input for menu access.

Operation

Weighing Δ

The basic weighing function is available at all times.

Features:

- Zero the weighing instrument by pressing $\rightarrow 0 \leftarrow$
- Store the weight on the instrument as tare by pressing $\rightarrow T \leftarrow$
- Tare container weight automatically
- Enter tare weight through bar code scanner (Combics 1 plus and Combics 2 only)
- Enter tare weight through numeric keypad (Combics 1 plus only)
- Delete tare values by pressing $\square 0$ and $\rightarrow T \leftarrow$ or $\square C$ and $\rightarrow T \leftarrow$ (Combics 1 plus only)
- Press $\square Fn$ to toggle the display between:
 - Gross and net values,
 - 1st and 2nd weight unit, or
 - normal and 10-fold higher resolution
- Weighing with two weighing instruments (Combics 2 only)
- Individual data ID codes with numeric values for identifying weight values (Combics 1 plus only)
- Print weight value:
 - GMP-compliant printout
 - Automatic printing
 - Automatic data output (see "Data Interfaces")

Automatic Taring (Menu Item 3.7)

When menu item 3.7.2 is active, the first load placed on the weighing instrument that exceeds the specified minimum load is stored, at stability, in the tare memory.

The weighing instrument returns to the initial state when the load is less than 50% of the minimum load.

Minimum Load for Automatic Taring and Automatic Printing (Menu Item 3.5)

You can choose from the following settings for the minimum load:

1 digit (no minimum load)
2 digits
5 digits
10 digits
20 digits
50 digits
100 digits
200 digits
500 digits
1000 digits

The "digits" here refer to the intervals in the connected weighing instrument. If the interval of the connected instrument is 1 g, for example, and 1000 digits are required, the minimum load is 1000 g (=1000 digits).

If the interval of the connected instrument is 5 g and the same number of intervals is required as in the example above, the minimum load is 5000 g.

Once the load on the weighing instrument exceeds the specified minimum, the instrument is tared and/or a printout is generated, if the operating menu is configured for automatic taring (menu item 3.7.2) and/or automatic printing (menu item 7.13.2).

Automatic Printing (Menu Item 7.13)

When menu item 7.13.2 is active, the first weight value that exceeds the specified minimum load is printed. If the menu code for automatic taring is also active, the weighing instrument is only tared when the minimum load is exceeded; the value is not printed.

First Weighing Instrument Displayed (Combics 2 only)

You can define which weighing instrument shows the first weight value when you switch on the Combics, under $\rightarrow E \rightarrow L \rightarrow E$ (menu item 8.9).

Entering Tare Weight using a Bar Code Scanner (Combics 1 plus and Combics 2 only)

You can enter the tare value of a container using a bar code scanner. To do this, the "Store value as tare" ($\rightarrow E \rightarrow E$) menu item must be selected under "Setup > Bar code" in the operating menu. In this case, the value is stored as the tare automatically, without pressing the $\square Tare$ key. The contents of the tare memory are display in Info mode (press and hold $\square \rightarrow$).

Entering ID Codes using a Bar Code Scanner (Combics 1 plus only)

You can use a bar code scanner to enter ID codes. To do this, the "Store value as ID1" ($\rightarrow ID \rightarrow 1$) menu item must be selected under "Setup > Bar code" in the operating menu. In this case, the value is stored as ID1 automatically, without pressing the $\square ID1$ key.

To store the second ID code, the $\square ID2$ key must be pressed.

To view the stored ID codes:

- Press $\square Info$ and $\square ID1$
- Press $\square Info$ and $\square ID2$

Calibration/Configuration Counter on Standard Weighing Instruments

Purpose

These two mutually independent counters automatically keep track of changes made in calibration/adjustment parameters and in the operating menu. Counter values are stored in an EEPROM, and remain stored during the entire service life of this memory chip. To view the current values in the counter, press and hold the $\rightarrow 0 \leftarrow$ key (longer than 2 sec). The readout shows the "configuration counter" value for 3 seconds first (identified by "P"). Then the "calibration counter" value is shown for 3 seconds (identified by "C"). The information display closes automatically after 6 seconds.

Features of the Calibration Counter:

- Limited to a count of 9999
- Counter set to "C 0000" when the hardware is first put into operation
- Counter cannot be reset
- The counter value is updated ("1" is added) automatically following:
 - Successful calibration/adjustment or linearization
 - Changes in the user-defined calibration/adjustment or linearization weight (menu item 1.18)
 - Changes in any of the following parameters:
 - Function of the \square CAL \square key (menu item 1.9)
 - Zero-setting range (menu item 1.11)
 - Initial zero-setting range (menu item 1.12)
 - Resetting of the above parameters to factory settings (menu item 9.1.1)

Features of the Configuration Counter:

- Limited to a count of 9999
- Counter set to "P 0000" when the hardware is first put into operation
- Counter cannot be reset
- The counter value is updated ("1" is added) automatically following:
 - Changes in the following parameters:
 - Place of installation (menu item 1.1)
 - Application filter (menu item 1.2)
 - Stability range (menu item 1.3)
 - Taring (menu item 1.5)
 - Auto zero (menu item 1.6)
 - Weight unit 1 (menu item 1.7)
 - Weight unit 2 (menu item 3.1)
 - Weight unit 3 (menu item 3.3)
 - Resetting of the above parameters to factory settings (menu item 9.1.1)
 - Function of the \square Fn \square key changed to or from 10-fold higher resolution display
 - Activation or deactivation of application-dependent automatic taring (menu item 3.7)
 - Resetting of the application parameters to factory settings (menu code 9.1.1)

Device parameters

Password

You can prevent unauthorized changes in the device settings ("SEtUP") and application settings ("APPt," CombiCS 2 only) by assigning a password under "Setup > Code" (C0dE; see also the chapter entitled "Configuration").

Acoustic Signal

An acoustic signal is emitted when you press a key (active key: single beep; inactive key: double-beep). You can switch the acoustic signal off or on under "Setup > Utilities" (SEtUP, uT IL lE) (menu item 8.2).

Keys

In the Setup menu under "uT IL lE," you can block or release the keypad for input (menu item 8.3).

Automatic Power-off

In the Setup menu under "uT IL lE" you can configure the CombiCS to shut down automatically following a specified interval of no user activity (menu item 8.7).

Display Backlighting

In the Setup menu under "uT IL lE" you can choose from the following settings for the display backlighting:

- On (8.8.1)
- Off (8.8.2)
- Shut off after the specified time period has elapsed (8.8.3)

Timer Mode

In the Setup menu under "uT IL lE" you can set the timer interval to 2, 4 or 10 minutes (menu item 8.9).

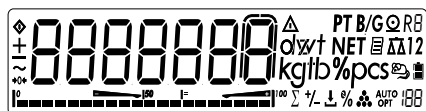
Operating

Example


Switch on the Combics, zero the weighing instrument, tare the container weight, place sample in container, toggle display to gross weight, 2nd weight unit or 10-fold increased resolution



Switch on the Combics



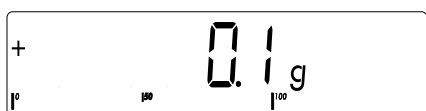
All display segments are shown for approx. 1 second (self-test)




Display with no load on weighing instrument



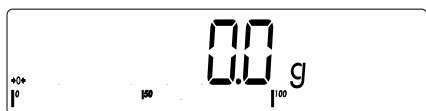
Zero the weighing instrument



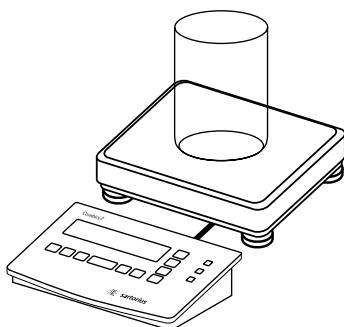Display with no load on weighing instrument




Place container on weighing instrument

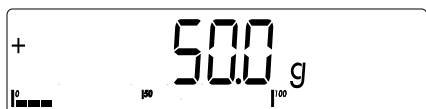


Container weight is displayed



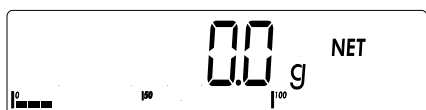Tare the weighing instrument



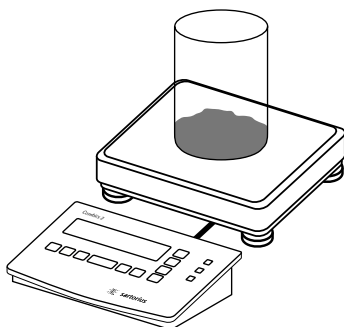Display with tared container on weighing instrument



Place sample in container (in this example, 120.2 g)



Display with tared weighing instrument and sample in container



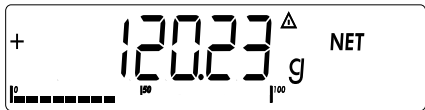
Toggle display; readout depends on your settings:



gross weight
(in this example, 50 g for container + 120.2 g substrate) or



display in 2nd weight unit (in this example, kg) or



display with 10-fold increased resolution



Return to previous readout
(if 10-fold increased resolution is shown, display returns to previous readout automatically after 10 seconds)



Print results

```

ACE HARDWARE
GOETTINGEN
24.02.2002    15:10
-----
G#   +   170.2 g
T    +   50.0 g
N    +   120.2 g
-----
  
```

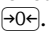

Operating the Combics 1 plus


Example:

Combics 1 plus: Weighing with tare value entered using the numeric keypad; generate printout of results



Switch on the Combics 1 plus

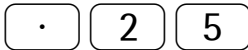
The automatic self-test runs. Once a readout is shown, Combics 1 plus is automatically zeroed and ready to operate. With no load on the weighing instrument, you can zero the instrument at any time by pressing .



Enter the tare weight in the current weight unit using the keypad (in this example, 0.25 kg).



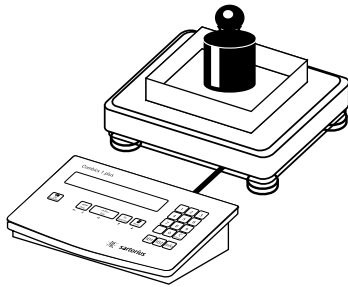Tare the weighing instrument




Place container with sample on the weighing instrument



Read the result



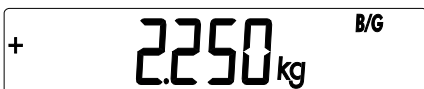Toggle the display from net to gross weight values. The display shows



the gross weight (in this example, 0.250 kg for the container plus 2.000 kg for the sample)

Fn

Toggle to display of net value

[=]

Print the results

Start of GMP header (only if GMP-compliant printout is configured)

05.04.2004 10:09
Model CW1NP1-30ED-LOCE
Ser.no. 12345678
Vers. 1.0001.04.4
BVers. 01-33-01

ACE HARDWARE
GOETTINGEN
Batch no. 123456
Customer 6.789
05.04.2004 10:09

End of GMP header
Header lines

ID code 1
ID code 2

G# + 2.250 kg
T + 0.000 kg
PT2 + 0.250 kg
N + 2.000 kg

05.04.2004 10:10
Name :

GMP footer (only if GMP-compliant printout is configured)

End of GMP footer

C + **→T←**

Delete tare weight entered through the keypad

or

0 + **→T←**

Operating the Combics 1 plus

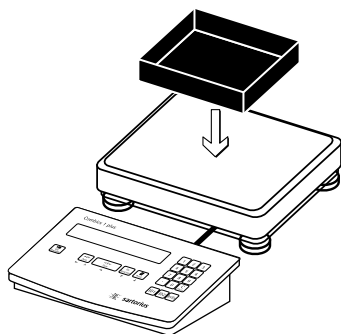
Example:

Combics 1 plus: Weighing with varying tare values; generate printout of results, delete tare values



Switch on the Combics 1 plus

The automatic self-test runs. Once a readout is shown, the Combics is automatically zeroed and ready to operate. With no load on the weighing instrument, you can zero the instrument at any time by pressing



Place empty container on the weighing instrument

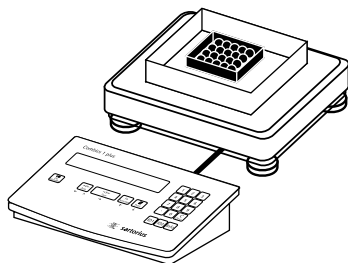


Tare the weighing instrument


Note: With the automatic tare function enabled, you do not need to press to tare the weighing instrument; the tare weight is stored automatically when you place the container on the weighing instrument.



Wait until a zero value is displayed together with the NET symbol.



Place packaged sample in the container

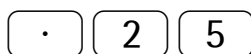

Enter the tare weight of the packaging in the current weight unit using the keypad (in this example, 0.25 kg).



Save the package weight. The package tare is added to the container tare.



Read off net weight



Print the results

G#	+	6.433	kg
T	+	4.183	kg
PT2	+	0.250	kg
N	+	2.000	kg



Enter a zero ("0") using the keypad



Save the value entered. This deletes tare values; the display shows the gross value

⬆
+

6.433 kg



Print the results

G#	+	6.433	kg
T	+	0.000	kg
N	+	6.433	kg

Operating the Combics 1 plus

Individual ID Codes (Identifiers)

You can assign ID codes (such as product name, batch number, etc.) while weighing, for identification of measured values.

Features

- Assign up to two ID codes.
- Assign both a name and a value for each ID code.
- The name is left-justified and the value is right-justified on the printout. If the entire code is too long for one line, additional lines are printed.
- Enter ID code names in Setup under: *SEtUP: PrtProt* (printouts):
7.4.3 (ID1)
and
7.4.4 (ID2)
- Enter up to 20 characters for the ID code name. No more than 11 characters are displayed during input; all 20 characters are printed.
- Maximum length for values: 21 characters.
- Enter numeric values for ID codes using the numeric keypad and press (ID1) or (ID2) to save.
- To delete the last character entered in the ID code value: Press (C)
To delete the entire ID code: Press (C) and then (ID1) or (ID2)
- If both the name and value fields are empty, no ID code is printed.

- Print ID codes: Refer to the Setup menu for menu item numbers.
Data is output to COM1 menu item 7.7.x) or UniCOM (menu item 7.8.x).
- Show ID codes:
Press (Info) and then (ID1) or (ID2)
- To store values entered using a bar code scanner: Scan the value for ID1


Function Keys

(ID1) Store the input as
or value for first or
(ID2) second ID code.

(C) Delete the selected ID code value

(Info) Display ID codes

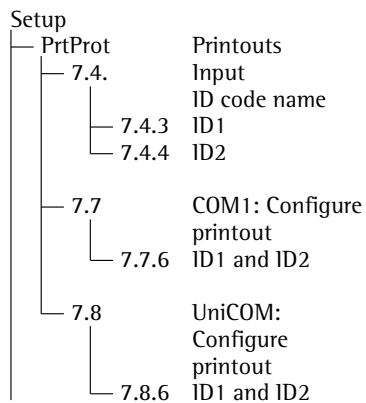Application Parameters: ID Codes



Example:
Entering ID code values.
Enter “123.456” and “678.9” as values for ID codes 1 and 2.
See “Entering the Password” in the chapter entitled “Configuration” for details on how to enter ID code names.

0000kg

1

2

3

.

4

5

6

Enter value for ID code 1 (in this example, 123.456)

123456.

ID1

Confirm value for the first ID code

0000kg

6

7

8

.

9

Enter value for ID code 2 (in this example, 678.9)

6789.

ID2

Confirm the value for the second ID code

6.433kg

Place container with sample on the weighing instrument

[Print]

Print weight value (perform further weighing operations as desired, and print)

ID1		123.456
ID2		6.789
24.02.2003		10:09

Ser.no.		12345678
G#	+	6.433 kg
T	+	0.000 kg
N	+	6.433 kg

ID code 1
ID code 2

C

ID1

Delete ID code:
You can delete each ID code individually when the weighing series has been completed.

C

ID2

Operation

Calibration and Adjustment

Purpose

The accuracy of weighing results must be carefully controlled. This is achieved through calibration and adjustment.

Perform calibration to determine the difference between the value displayed and the actual weight on the weighing instrument. Calibration does not entail making any changes within the weighing instrument.

The adjustment procedure actually eliminates the difference between the readout and the actual weight, or reduces it to a level within the permissible tolerance limits.

Features

Which of the following features are available depends on the weighing instrument:

- External calibration/adjustment with the default weight value or standard weight (1.9.1) (not available on verified instruments)
- External calibration/adjustment with a user-defined weight (1.9.3) (not available on verified instruments)
- Internal calibration/adjustment for IS platform (1.9.4), WP2 only
- Block the $\rightarrow T \leftarrow$ key to prevent use of the two functions described above (1.9.10)
- Calibrate first; then adjust automatically (1.10.1) (not available on verified instruments)
- Calibrate, then prompt for manual input of adjustment command (1.10.2)
- Calibration prompt: flashing $\Delta \Delta$ symbol (1.15.2). If more than one weighing instrument is connected, the instrument number is also displayed.
- Block external calibration/adjustment (1.16.2)

- Display altitude and geographical latitude or acceleration of gravity after $\Delta \Delta$ is shown at the beginning of the calibration procedure (menu item 8.12.2). These values are shown only if they have been entered in the service menu and activated.
For each of these parameters, the term is displayed first (ALT or $GRAV$) for 1 second, and then the corresponding value is displayed continuously until you press $\rightarrow T \leftarrow$.

Note

On verified weighing instruments, the external calibration/adjustment function is available only when the menu access switch is in the “open” position, which entails breaking the verification seal (refer to the chapter entitled “Service”). The equipment must be re-verified after the seal has been broken.

Example

External calibration and manual adjustment with default weights

Setup menu settings:
1.9.1; 1.10.2



Unload and zero the weighing instrument



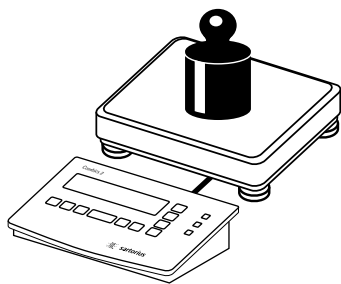
Start calibration (e.g., when calibration prompt is flashing: $\Delta \Delta$)




The following is displayed for 2 seconds



You are prompted to place the required weight on the weighing instrument (e.g., 10 kg)



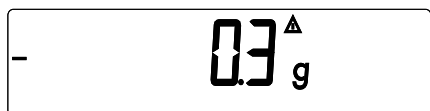Position the calibration weight on the weighing instrument



The difference between the weight value and the true mass is displayed, with \pm sign.

```

External calibration
Nom. + 10000.0 g
Diff. - 0.3 g
-----
  
```

Calibration record is printed, if adjustment was not performed and the process was stopped by pressing $\rightarrow 0 \leftarrow$



Activate calibration/adjustment manually (press the $\rightarrow 0 \leftarrow$ key to cancel)

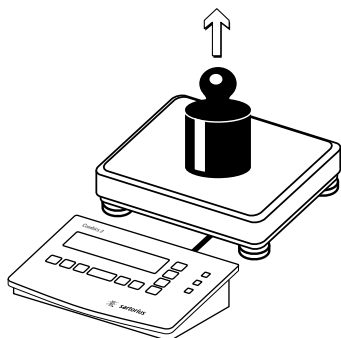


The calibration weight is displayed at the conclusion of calibration

```

-----
14.01.2002    13:00
Model        CISL2
Ser.no.      12345678
Vers.        1.1007.12.1
BVers.       01-25-01
-----
External calibration
Nom. + 10000.0 g
Diff. - 0.3 g
External adjustment
Diff. + 0.0 g
-----
14.01.2002    13:02
Name:
-----
  
```

A GMP-compliant printout is generated





Unload the weighing instrument

Operating the Combics 2


Counting ☼


With the Counting application, you can determine the number of parts that each have approximately equal weight.

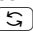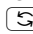Features:


- Store the weight on the weighing instrument as reference weight
- Enter reference weight using a bar code scanner
- Enter tare weight using a bar code scanner
- Automatic reference sample updating (user-definable)
- Counting with two weighing instruments
- Toggle the display between piece count and weight by pressing 
- Info mode for display of average piece weight and reference sample quantity by pressing  (> 2 sec)

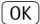
Before the quantity on the weighing instrument can be calculated, the average piece weight must be entered in the application. There are 2 ways to do this with the Combics:

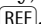- By placing the number of parts defined as the reference sample quantity on the weighing instrument and pressing  to store the average piece weight.

The reference sample quantity is shown in the numeric display, and can be changed by pressing .

How the reference weight is calculated depends on the menu setting for resolution (3.9). Either the value is rounded off in accordance with the display resolution, or the display resolution is increased 10-fold (+1 decimal place) or 100-fold (+ 2 decimal places), or maximum internal resolution is applied.

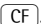- By entering the reference piece weight (i.e., the weight of one piece) using a bar code scanner (menu setting: "Setup> Bar code > Store value as reference (REF)"). In this case, the value is stored as a reference automatically, without pressing the  key.

This value remains active in the reference memory until you delete it by pressing , overwrite it or until you select a different application. It also remains in memory when you switch to a different application program, or switch off the Combics.

Application Parameters: Counting

3. 5.	Minimum Load for Automatic Taring and Automatic Printing
3. 5. 1 *	1 digit
3. 5. 2	2 digits
3. 5. 3	5 digits
3. 5. 4	10 digits
3. 5. 5	20 digits
3. 5. 6	50 digits
3. 5. 7	100 digits
3. 5. 8	200 digits
3. 5. 9	500 digits
3. 5. 10	1000 digits
3. 6.	Minimum Load for Initialization
3. 6. 1 *	1 digit
3. 6. 2	2 digits
3. 6. 3	5 digits
3. 6. 4	10 digits
3. 6. 5	20 digits
3. 6. 6	50 digits
3. 6. 7	100 digits
3. 6. 8	200 digits
3. 6. 9	500 digits
3. 6. 10	1000 digits
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3. 8.	Start Application with Most Recent Application Data when Combics is Switched On
3. 8. 1	Automatic (on)
3. 8. 2 *	Manual (off)
3. 9.	Resolution for Calculation of Reference Value
3. 9. 1 *	Display resolution
3. 9. 2	Display resolution + 1 decimal place
3. 9. 3	Display resolution + 2 decimal places
3. 9. 4	Internal resolution
3.11	Storage Parameter
3.11. 1*	At stability
3.11. 2	At increased stability
3.12.	Reference Sample Updating
3.12. 1	Off
3.12. 3*	Automatic
3.13.	Reference Weighing Instrument
3.13. 1*	No reference instrument selected
3.13. 2	WP1
3.13. 3	WP2

Storage Parameter

The reference weight is saved when the weighing instrument has stabilized. "Stability" is defined as the point at which fluctuation of a measured value lies within a defined tolerance range. The narrower the tolerance range, the more stable the weighing instrument is at "stability." This setting is also applied when you zero the weighing instrument.

Under menu item 3.11 you can determine whether the value is saved "At stability" (normal tolerance range) or "At increased stability" (narrower tolerance range).

If you select "At increased stability," the average piece weight stored will be more accurate and the results more reproducible, but the response time of the weighing instrument might be longer.

Minimum Load for Initialization

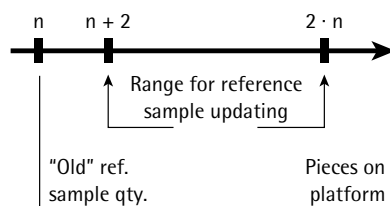
The minimum load for initialization is configured under menu item 3.6. If the load exceeds this limit, the weighing instrument can be initialized. If the load is too light, the following will occur when you try to save a value:

- The error code *lnF 29* is displayed
- A warning signal is emitted (double-beep)
- The weighing instrument is not initialized
- The preset reference sample quantity is stored

Reference Sample Updating

In the Application settings under 3.12, you can define whether the reference sample is updated automatically. Reference sample updating is performed automatically only when the following 6 criteria are met:

1. Menu item 3.12.3 is active
2. The current piece count exceeds the original piece count by at least two
3. The current piece count is less than twice the original piece count (does not apply for the first update if the piece count is entered using a bar code scanner).



4. The current piece count is less than 100.
5. The internally calculated piece count (such as 17.24 pcs) differs by less than ± 0.3 pcs from the whole number (17 pcs in this example).
6. The weighing instrument is stable in accordance with the defined stability parameter.

When automatic reference sample updating is active, the AUTO symbol is displayed in addition to the Counting symbol (⚖). When the reference sample has been updated, OPT is displayed below AUTO. During an updating operation, OPT and the updated piece count are displayed briefly in the main display.

At the conclusion of reference sample updating, a beep is sounded and the new reference weight and reference sample quantity are stored. Activate the "Info" mode to view the reference values (press and hold $\text{[S]} > 2$ sec).

Counting with Two Weighing Instruments

You can use two weighing instruments simultaneously with the Counting application. When using two instruments, you can choose from the following operating modes:

- Counting with two equivalent weighing instruments
- Counting with one reference weighing instrument and one counting platform

Counting with Two Equivalent Weighing Instruments

Use this mode when samples of widely varying weight are counted at one workstation. Count the lighter-weight pieces on one weighing instrument and the heavier pieces on another. When you press [AB] to toggle from one weighing instrument to the other, the application is re-initialized.

You can define which of the two weighing instruments is active in the display when the Combics is switched on (menu item 8.11). This is the first weighing instrument active when you switch on the Combics, regardless of the setting for automatic initialization of the Counting application.

Counting with One Reference Weighing Instrument and One Counting Platform

In this mode, a high-resolution weighing instrument with a relatively low maximum capacity is used as a reference weighing instrument. The weighing platform is used for weighing heavier samples, and has a high capacity with a relatively low resolution.

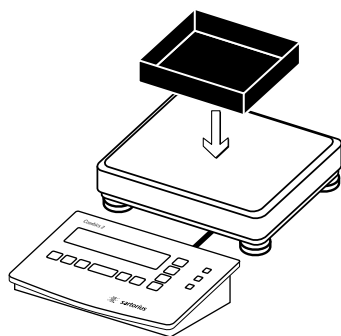
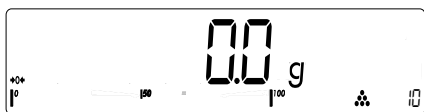
This way, you can both determine the reference sample quantity very precisely and count large amounts of parts, without requiring an expensive high-resolution, high-capacity weighing platform.

Specify which weighing instrument is the reference instrument under menu item 3.13. The system toggles automatically to the reference weighing instrument for initialization (*rEF* is displayed). Following initialization, the system toggles to the counting platform.

The system does not toggle automatically for automatic reference sample updating; the update is based on whichever instrument is active.

Operating the Combics 2

Example:



Place empty container on the weighing instrument

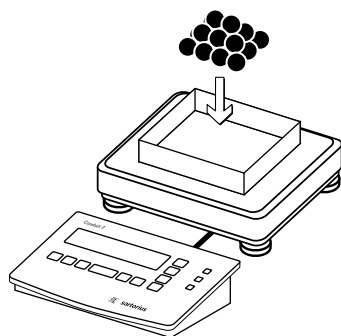


Tare the weighing instrument



(repeatedly, if necessary)


Set the desired reference sample quantity (in this example, 20)




Place the corresponding number of pieces (20) in the container




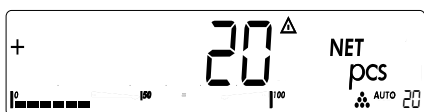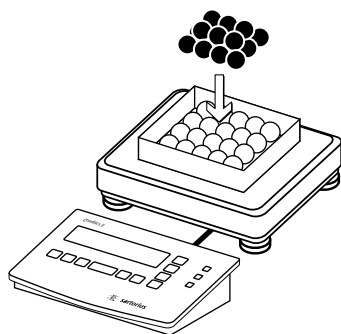Confirm reference sample weight





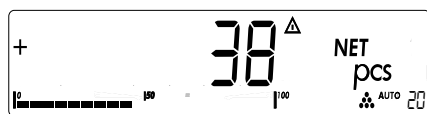
Add more parts to the container (in this example, 18 pcs)



Δ PCt is displayed if automatic reference sample updating is enabled




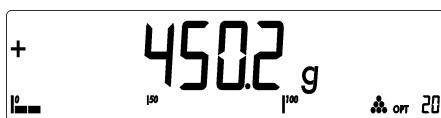Toggle the display from piece count to weight



< 2 sec.

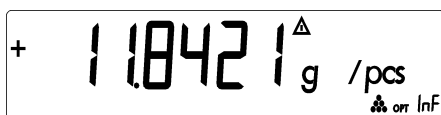
Toggle to Info mode



Piece count is displayed



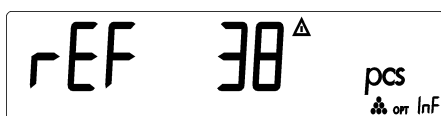Toggle to the next display mode



Reference sample quantity is displayed



(> 2 sec)

Return to weighing mode



Print results



G#	+	610.0 g
T	+	200.0 g
N	+	410.0 g
Qnt		38 pcs

Operating the Combics 2


Neutral Measurement


When the Neutral Measurement application is selected, you can use your weighing instrument to measure the length, surface and volume of parts that have roughly the same specific weight. The *o* symbol is displayed as the weight unit.

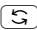Features:


- Store the weight on the weighing instrument as reference weight
- Enter reference weight using a bar code scanner
- Enter tare weight using a bar code scanner
- Toggle the display between weight and calculated value by pressing 
- Info mode for display of reference weight by pressing  (> 2 sec)


To calculate the result correctly, the average weight of a reference quantity must be known (in the example below, the reference is 1 meter of electrical cable). There are 2 ways to do this with the Combics:

- By placing the number of parts defined for the reference value on the weighing instrument and pressing  to store the average piece weight.

The reference value is shown in the numeric display, and can be changed by pressing .

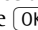How the reference weight is calculated depends on the menu setting for resolution (3.9). Either the value is rounded off in accordance with the display resolution, or the display resolution is increased 10-fold (+1 decimal place) or 100-fold (+ 2 decimal places), or maximum internal resolution is applied.

- By entering the reference weight (i.e., the weight of one piece) using a bar code scanner (menu setting: "Setup> Bar code > Store value as reference (REF)"). In this case, the value is stored as a reference automatically, without pressing the  key.

This value remains active in the reference memory until you delete it by pressing , overwrite it or until you select a different application. It also remains stored after you switch off the Combics.

Application Parameters: Neutral Measurement

- 3. 5. Minimum Load for Automatic Taring and Automatic Printing
 - 3. 5. 1 * 1 digit
 - 3. 5. 2 2 digits
 - 3. 5. 3 5 digits
 - 3. 5. 4 10 digits
 - 3. 5. 5 20 digits
 - 3. 5. 6 50 digits
 - 3. 5. 7 100 digits
 - 3. 5. 8 200 digits
 - 3. 5. 9 500 digits
 - 3. 5. 10 1000 digits
- 3. 6. Minimum Load for Initialization
 - 3. 6. 1 * 1 digit
 - 3. 6. 2 2 digits
 - 3. 6. 3 5 digits
 - 3. 6. 4 10 digits
 - 3. 6. 5 20 digits
 - 3. 6. 6 50 digits
 - 3. 6. 7 100 digits
 - 3. 6. 8 200 digits
 - 3. 6. 9 500 digits
 - 3. 6. 10 1000 digits
- 3. 7. Automatic Taring: 1st Weight Tared
 - 3. 7. 1 * Off
 - 3. 7. 2 On
- 3. 8. Start Application with Most Recent Application Data when Combics is Switched On
 - 3. 8. 1 Automatic (on)
 - 3. 8. 2 * Manual (off)
- 3. 9. Resolution for Calculation of Reference Value
 - 3. 9. 1 * Display resolution
 - 3. 9. 2 Display resolution + 1 decimal place
 - 3. 9. 3 Display resolution + 2 decimal places
 - 3. 9. 4 Internal resolution
- 3. 10. Decimal Places for Display of Results
 - 3. 10. 1 * None
 - 3. 10. 2 1 decimal place
 - 3. 10. 3 2 decimal places
 - 3. 10. 4 3 decimal places
- 3. 11. Storage Parameter
 - 3. 11. 1 * At stability
 - 3. 11. 2 At increased stability
- 3. 13. Reference Weighing Instrument
 - 3. 13. 1 * Off
 - 3. 13. 2 WP1
 - 3. 13. 3 WP2

Storage Parameter

The reference weight is saved when the weighing instrument has stabilized. "Stability" is defined as the point at which fluctuation of a measured value lies within a defined tolerance range. The narrower the tolerance range, the more stable the weighing instrument is at "stability." This setting is also applied when you zero the weighing instrument.

Under menu item 3.11 you can determine whether the value is saved "At stability" (normal tolerance range) or "At increased stability" (narrower tolerance range).

If you select "At increased stability," the average piece weight stored will be more accurate and the results more reproducible, but the response time of the weighing instrument might be longer.

Decimal Places for Display of Results

In neutral measurement, not only whole numbers but also decimal numbers (for example, 1.25 g electrical cabling) can be displayed. The number of decimal places displayed in neutral measurement is configured under menu item 3.10. The measured result can be displayed with 0, 1, 2 or 3 decimal places.

Minimum Load for Initialization

The minimum load for initialization is configured under menu item 3.6.

If the load exceeds this limit, the weighing instrument can be initialized. If the load is too light, the following will occur when you try to save a value:


- The error code *inF 29* is displayed
- A warning signal is emitted (double-beep)
- The weighing instrument is not initialized
- The preset reference value is stored

Neutral Measurement with Two Weighing Instruments

You can use two weighing instruments simultaneously with the Neutral Measurement application. When using two instruments, you can choose from the following operating modes:

- Neutral measurement with two equivalent weighing instruments
- Neutral measurement with one reference weighing instrument and one measurement platform

Neutral Measurement with Two Equivalent Weighing Instruments

Use this mode when samples of widely varying weight are measured at one workstation. Measure the lighter-weight pieces on one weighing instrument and the heavier pieces on another. When you press  to toggle from one weighing instrument to the other, the application is re-initialized.

You can define which of the two weighing instruments is active in the display when the Combics is switched on (menu item 8.11). This is the first weighing instrument active when you switch on the Combics, regardless of the setting for automatic initialization of the Neutral Measurement application.

Neutral Measurement with One Reference Weighing Instrument and One Measuring Platform

In this mode, a high-resolution weighing instrument with a relatively low maximum capacity is used as a reference weighing instrument. The measuring platform has a high capacity, but a relatively low resolution.

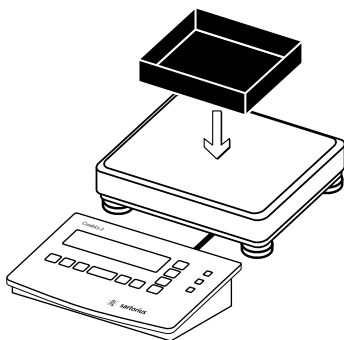
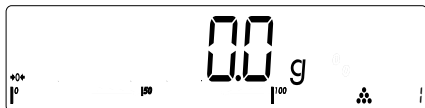
This allows you to both determine the reference value with high resolution; i.e., very precisely, and to measure large samples, without requiring an expensive high-resolution, high-capacity weighing platform.

Specify which weighing instrument is the reference instrument under menu item 3.13. The system toggles automatically to the reference instrument for initialization (*rEF* is shown in the main display). Following initialization, the system toggles to the measuring platform.

Operating the Combics 2

Example:


Determine the length of an amount of electrical cable after weighing in the defined reference unit value



Place empty container on the weighing instrument





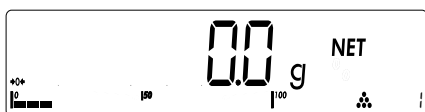Tare the weighing instrument



(repeatedly, if necessary)

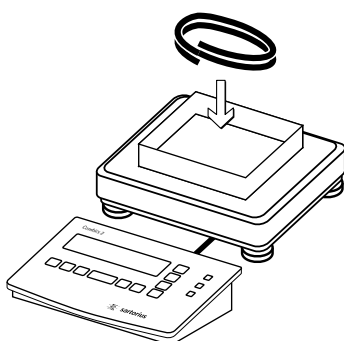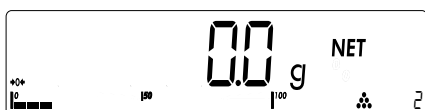
Set the desired reference value (in this example, 2)




Place a sample corresponding to the reference quantity in the container
(in this example, 2 meters of electrical cable)



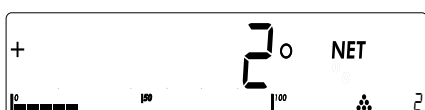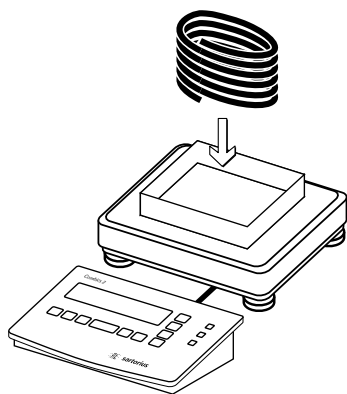Confirm reference sample weight





Remove the reference material and place the sample to be measured on the weighing instrument (in this example, 8 meters of electrical cable)



Print results

G#	+	734.1	g
T	+	200.0	g
N	+	534.1	g
Qnt		8	o

Operating the Combics 2

Checkweighing %

When the Checkweighing application is selected, you can check whether sample weights correspond to a specified target weight; i.e., whether the weight on the weighing instrument is within a given tolerance range. The tolerance range is defined by upper and lower limits. The result is displayed in the main indicator, in the bar graph and by color-coded LEDs.

Features:

- The target value can be taken over as a weighed value from a weighing instrument, and the tolerance limits are defined as a perceptual deviation from the target value. The following percentages can be selected as the deviation: 0.1%, 0.2%, 0.5%, 1%, 1.5%, 2%, 3%, 5% or 10%.
- The target value, lower tolerance limit (minimum) and upper tolerance limit (maximum) can be taken over as weighed values from the weighing instrument.
- Target and tolerance limits checked during input; values must conform to: Upper limit \geq Target \geq Lower limit \geq 1 digit
- Checkweighing range: either 30% to 170% of the target, or from 10% to infinity
- Application started automatically with most recent application data when Combics switched on
- Automatic taring
- Automatic printing
- Toggle the display between weight and tolerances limits by pressing $\left[\begin{smallmatrix} \text{S} \end{smallmatrix} \right]$. When tolerances are displayed, weights exceeding the tolerance limits are shown with "LL" (too low) or "HH" (too high).
- Digital input/output interface
- Info mode for display of tolerance limits by pressing $\left[\begin{smallmatrix} \text{S} \end{smallmatrix} \right]$ (> 2 sec)

Checkweighing entails comparing the current weight value to a defined target. The target value has a tolerance range which can be entered either as an absolute value or percentage (menu item 4.5).

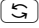
- Entering the tolerance range as an absolute value (weighed value) (menu item 4.5.1):
Start initialization by pressing $\left[\text{OK} \right]$; the middle segments of the bar graph flash to prompt the placement of the weight on the weighing instrument. Place the weight on the instrument and press $\left[\text{OK} \right]$ to store.
The bar graph segment for the lower limit flashes to prompt the weight for lower limit. Place the weight on the weighing instrument and press $\left[\text{OK} \right]$ to store.
The bar graph segment for the upper limit flashes to prompt the weight for upper limit. Place the weight on the weighing instrument and press $\left[\text{OK} \right]$ to store.
- Entering the tolerance range as a percentage (menu item 4.5.2):
A value for the percentage is shown in the numeric display (lower right-hand corner) together with the "%" sign. Press $\left[\text{REF} \right]$ to change the percentage value (0.1%, 0.2%, 0.5%, 1%, 1.5%, 2%, 3%, 5%, 10%).
Start initialization by pressing $\left[\text{OK} \right]$; the middle segments of the bar graph flash to prompt the placement of the weight on the weighing instrument. Place the weight on the instrument and press $\left[\text{OK} \right]$ to store. To can change the percentage value for the tolerance limits again before the target value is stored, press $\left[\text{REF} \right]$.

Before the new initialization, the previous initialization values must be deleted by pressing $\left[\text{CF} \right]$.

Application Parameters: Checkweighing

- 3. 5. Minimum Load for Automatic Taring and Automatic Printing
 - 3. 5. 1 * 1 digit
 - 3. 5. 2 2 digits
 - 3. 5. 3 5 digits
 - 3. 5. 4 10 digits
 - 3. 5. 5 20 digits
 - 3. 5. 6 50 digits
 - 3. 5. 7 100 digits
 - 3. 5. 8 200 digits
 - 3. 5. 9 500 digits
 - 3. 5. 10 1000 digits
- 3. 7. Automatic Taring: 1st Weight Tared
 - 3. 7. 1 * Off
 - 3. 7. 2 On
- 3. 8. Start Application with Most Recent Application Data when Combics is Switched On
 - 3. 8. 1 Automatic (on)
 - 3. 8. 2 * Manual (off)
- 4. 2. Checkweighing Range
 - 4. 2. 1 * 30% to 170%
 - 4. 2. 2 10% to infinity
- 4. 3. Activate Control Line for "Set" as:
 - 4. 3. 1 * "Set" output
 - 4. 3. 2 Ready to operate
- 4. 4. Activation of Outputs
 - 4. 4. 1 Off
 - 4. 4. 2 Always active
 - 4. 4. 3 Active at stability
 - 4. 4. 4 * Active within checkweighing range
 - 4. 4. 5 Active at stability within the checkweighing range
- 4. 5. Parameter Input
 - 4. 5. 1 * Min, max, target
 - 4. 5. 2 Only target with percent limits
- 4. 6. Automatic Printing
 - 4. 6. 1 * Off
 - 4. 6. 2 On
 - 4. 6. 3 Only values within tolerance
 - 4. 6. 4 Only values outside tolerance

Display

The result of a measurement is shown either as a weight value or in relation to the target. You can toggle between these two display modes by pressing .

- **Weight display mode:**
The main display always shows the weight value, even if the value is outside the tolerance range. The bar graph is displayed with symbols indicating lower limit, target and upper limit. It shows a logarithmic display of the current load if the weight is anywhere from 0 to the minimum load, and a linear display for weights beyond that range. The LEDs indicate the following:
Yellow: weight value > upper tolerance limit
Green: weight value is within tolerance
Red: weight value < lower tolerance limit
If no LED lights up:
– the application is not completely initialized, or
– the weight value is outside the checkweighing range (see menu item 4.2)
– the weighing instrument has not stabilized
- **Tolerance limit display mode:**
As “Weight display mode” above, with the exception that:
– *LL* is shown on the main display if the weight value is lower than the target, and
– *HH* is shown on the main display if the weight value is higher than the target

Digital Input/Output Interface

The Checkweighing application supports the digital input/output interface.

There are 4 control lines, or outputs, which are activated as follows (also refer to the diagram):

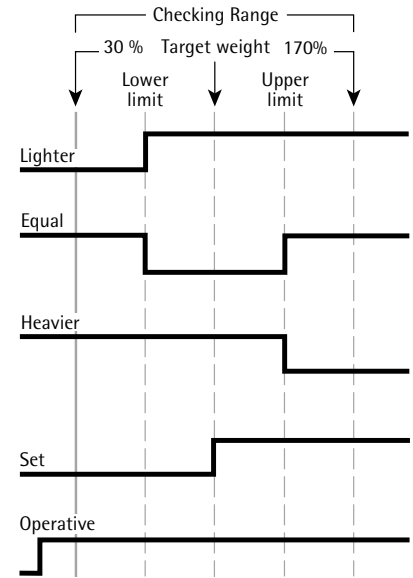
- **Lighter**
- **Equal**
- **Heavier**
- **Set**
The “SET” output normally changes its voltage level when the load is near the target weight. Alternatively, you can assign the “Ready-for-use” function to this port (menu item 4.3.2).

Under menu item 4.4, you can define whether these control ports are

- inactive (4.4.1)
 - always active (4.4.2)
 - active at stability (4.4.3)
 - active within the check range, or
 - active at stability within the check range
- For example, you can use this function to show the weighed or measured result on a simple external indicator, similar to the 3 LEDs on the Combics 2.


All data output ports have a high voltage level when:

- the application has not been initialized,
- the weighing instrument is not at stability and one of the “at stability...” parameters, (4.4.3 or 4.4.5) is selected, or
- the weight is not within the check range (4.4.4).



Digital I/O Interface

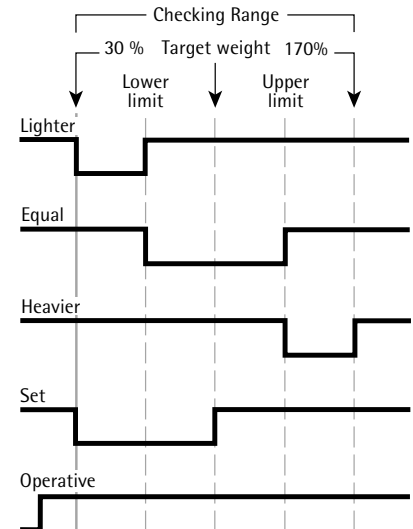“SET” control line: set and control lines:
Always active/Active at stability



Digital I/O Interface

“SET” control line: set and control lines:
active within checkweighing range/
Active within checkweighing range at stability

Output line specifications:

- In the inactive state, the levels are set to “high:” >3.7 V/+4 mA
- In the active state, the levels are set to “low:” <0.4 V/-4 mA

⚠ The output lines are not short-circuit proof!

Operating the Combics 2

Example:
Initialize the Checkweighing application by taking over the weighed value as target; select percentage to define tolerance range (menu item 4.5.2)

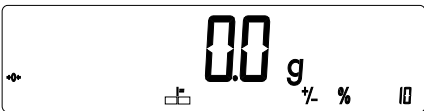
CF

Delete any existing initialization data




OK

Begin initialization



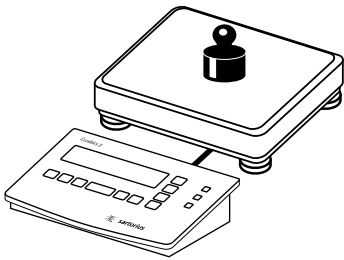The middle segments of the bar graph flash to prompt the placement of the weight on the weighing instrument




Place the weight for the target on the weighing instrument (in this example, 100.0 g)




REF repeatedly, if necessary

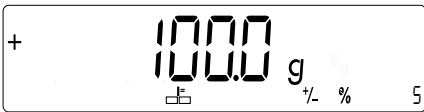Select the desired percentage for tolerance limits (in this example, 5%)



OK

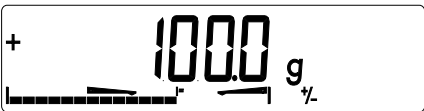Store target and calculate tolerance range based on the selected percentage



Optional:  < 2 sec.

Toggle to Info mode to view stored values



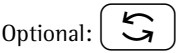
Target value is displayed



Toggle to the next display mode



Percentage for limits is displayed



Toggle to the next display mode



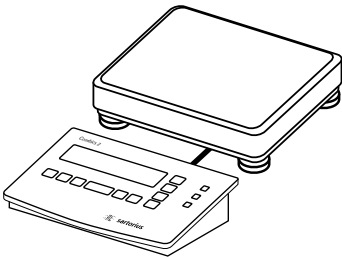
Weight value for the lower limit (minimum) is displayed



Toggle to next display mode (weight value for the upper limit is displayed)



Exit the Info mode




Unload weighing instrument




Operating the Combics 2

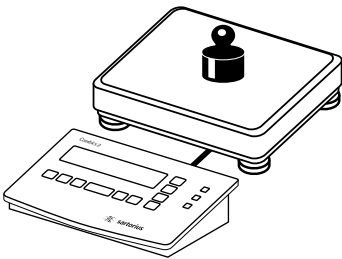Example:
Check the weight of a sample; use Tolerance Limit display mode



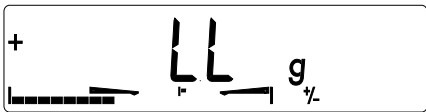
Place a load of unknown weight on the weighing instrument



If the weight is under the lower tolerance limit, LL is shown in the tolerance limit display (load is too light) (the weight display shows the measured weight value)



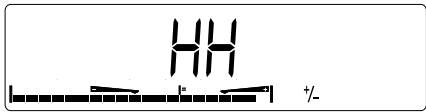If the weight is within the tolerance range (in this example, 103.2 g), the weight value is displayed



If the weight is over the upper tolerance limit, HH is shown in the tolerance limit display (load is too heavy) (the weight display shows the measured weight value)



Print results

ACE HARDWARE		
GOETTINGEN		
19.03.2002		15:43

Setp	+	100.0 g
Min	+	95.0 g
Max	+	105.0 g
G#	+	103.2 g
T	+	0.0 g
N	+	103.2 g
Lim	+	3.20 %
W.Diff	+	3.2 g

Printout

Target value
Minimum
Maximum

Gross weight
Tare weight
Net weight

Percentage of deviation from target*
Absolute difference from target

* In Tolerance Limit display mode:
If the weight is lighter than the target, the display shows: **S t a t** **LL**
If the weight is heavier than the target, the display shows: **S t a t** **HH**

Operating the Combics 2

Classification $\frac{\%}{\square}$ \square

With the Classification application, you can determine whether the weight of a given sample lies within the limits of a defined weight class.

Features:

- Configure 3 or 5 classes (menu item 4.8)
- Define contiguous classes
- Define classes that cover the entire weighing range of the weighing instrument
- Range below the specified minimum load is designated "Class 0"
- Define the upper limit of a given class by storing weight on weighing instrument or by entering a weight value and a percentage
- Show the current weight in the main display as a weight value or as belonging to a certain class
- Class of current weight also indicated by 1 LED (when using 3 classes) or 2 LEDs (when using 5 classes) 2 LEDs (when using 5 classes)
- Toggle the display between weight and class by pressing \square .
- Digital input/output interface
- Info mode for display of class limits by pressing \square (> 2 sec)

To use the Classification application, you need to enter the delimiters that separate one class from another.

- Start the initialization by pressing \square . The lower limit of Class 1 is defined by the preset minimum load. Weight values below minimum load are classified as belonging to Class 0, in part to prevent an unloaded weighing instrument from displaying a value that designates a class.

You can enter the other class delimiters as absolute values or percentages (menu item 4.9).

- Entering the class delimiters as absolute values (weighed values) (menu item 4.9.1):
For each upper limit (except that of the highest class), place a load on the weighing instrument and press \square to store it as the upper limit of a class. The bar graph and the numeric display show the limit currently being configured.
Each time you store a value for an upper limit, a check is run to make sure the value is equal to or greater than the upper limit of the previous class. If the value is invalid, an acoustic signal is emitted; the limit must be entered again.

- Entering delimiters as percentages (menu item 4.9.2):
Define the upper limit of Class 1 by placing the corresponding weight on the weighing instrument and pressing \square . The value is equal to 100 %. The upper limit of Class 2 is equal to 100 % plus the percentage you select by pressing \square (1, 2, 5, 10, 15, 30, 50, 70, 100, 150, 199 %).
Example: A load of 100 g is stored from the weighing instrument as the upper limit of Class 1. The value 15% is then entered, defining the upper limit of Class 2 as 115 g.
In this case, the weight ranges when working with 5 classes are:
Class 0: up to the minimum load
Class 1: Minimum load – 100 g
Class 2: 100 g – 115 g
Class 3: 115 g – 130 g
Class 4: 130 g – 145 g
Class 5: 145 g – maximum capacity

Before a new initialization, the previous initialization values must be deleted by pressing \square .

Application Parameters: Classification

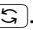
- 3. 5. Minimum Load for Automatic Taring and Automatic Printing
 - 3. 5. 1 * 1 digit
 - 3. 5. 2 2 digits
 - 3. 5. 3 5 digits
 - 3. 5. 4 10 digits
 - 3. 5. 5 20 digits
 - 3. 5. 6 50 digits
 - 3. 5. 7 100 digits
 - 3. 5. 8 200 digits
 - 3. 5. 9 500 digits
 - 3. 5.10 1000 digits
- 3. 6. Minimum Load for Initialization and Defining the Class 1 Lower Limit
 - 3. 6. 1 * 1 digit
 - 3. 6. 2 2 digits
 - 3. 6. 3 5 digits
 - 3. 6. 4 10 digits
 - 3. 6. 5 20 digits
 - 3. 6. 6 50 digits
 - 3. 6. 7 100 digits
 - 3. 6. 8 200 digits
 - 3. 6. 9 500 digits
 - 3. 6.10 1000 digits
- 3. 7. Automatic Taring: 1st Weight Tared
 - 3. 7. 1 * Off
 - 3. 7. 2 On
- 3. 8. Start Application with Most Recent Application Data when Combics is Switched On
 - 3. 8. 1 Automatic (on)
 - 3. 8. 2 * Manual (off)
- 4. 3. Activate Control Line for "Set" as:
 - 4. 3. 1 * "Set" output
 - 4. 3. 2 Ready to operate (for process control systems)
- 4. 7. Activation of Outputs
 - 4. 7. 1 Off
 - 4. 7. 2 Always active
 - 4. 7. 3 * Active at stability
- 4. 8. Number of Classes
 - 4. 8. 1 * 3 classes
 - 4. 8. 2 5 classes
- 4. 9. Parameter Input
 - 4. 9. 1 * Weight values
 - 4. 9. 2 Percentage
- 4.10. Automatic Printing
 - 4.10. 1* Off
 - 4.10. 2 On











Operating the Combics 2

Minimum Load for Initialization and Defining the Class 1 Lower Limit

The minimum load defines the lower limit for Class 1 (i.e., the delimiter between Classes 0 and 1). The application can only be initialized (to define the upper limit of Class 1) when the load on the weighing instrument exceeds the minimum load.

Display

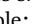The result of a given measurement is shown as either a weight value or the class number. You can toggle between these two display modes by pressing .

- Weight display; 3 classes:
The current weight is shown in the display; for example: 108.7 g
Class 1 Bar graph:  Numeric display: 1
LED: red is lit
Class 2 Bar graph:  Numeric display: 2
LED: green is lit
Class 3 Bar graph:  Numeric display: 3
LED: yellow is lit
- Weight display; 5 classes
The current weight is shown in the display; for example: 108.7 g
Class 1 Bar graph:  Numeric display: 1
LED: red is lit
Class 2 Bar graph:  Numeric display: 2
LED: red and green are lit
Class 3 Bar graph:  Numeric display: 3
LED: green is lit
Class 4 Bar graph:  Numeric display: 4
LED: green and yellow are lit
Class 5 Bar graph:  Numeric display: 5
LED: yellow is lit
- Class display; 3 classes
The current class is shown in the display; for example:  1 and the bar graph indicates the current weight. The LEDs are lit as for the weight display described above.
- Class display; 5 classes
The current class is shown in the display; for example:  1, and the bar graph indicates the current weight. The LEDs are lit as for the weight display described above.

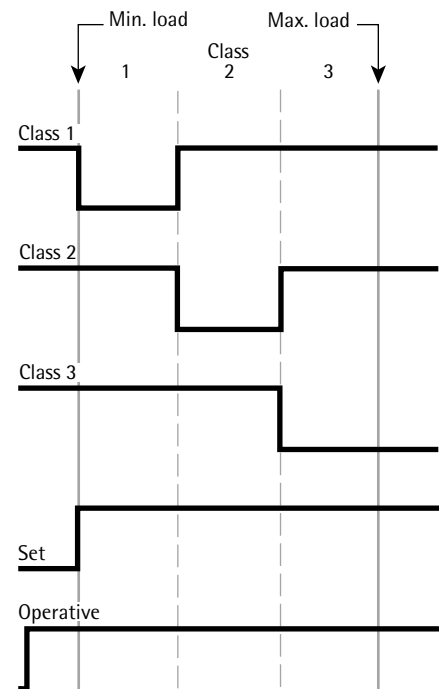
Digital Input/Output Interface

The Classification application supports the digital input/output-interface. There are 4 control lines, or outputs, which are activated as follows (also refer to the diagram):


- Lighter
 - Equal
 - Heavier
 - Set
- The "Set" line usually indicates that the minimum load is exceeded. Alternatively, you can assign the "Ready-for-use" function to this port.

Under menu item 4.7, you can define whether these control ports are

- inactive (4.7.1),
 - always active (4.7.2), or
 - active at stability (4.7.3).
- This makes it possible, for example, to connect a simple indicator for weighing or calculation results, similar to the 3-segment checkweighing display on the Combics, or control an automatic process.



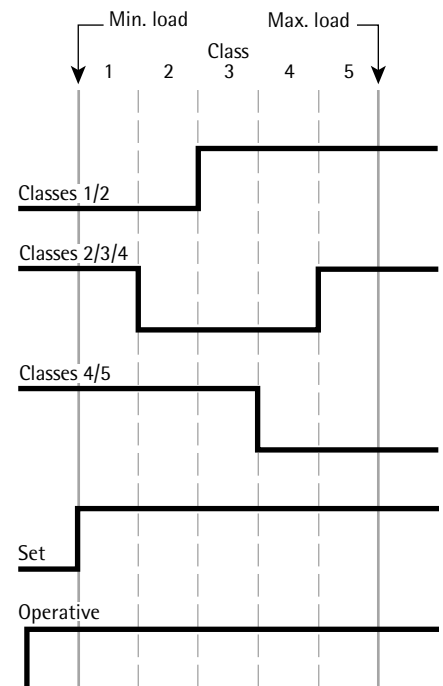Digital I/O Interface
Control lines when working with 3 classes



Digital I/O Interface
Control lines when working with 5 classes

Example:
Initialize the Classification application by taking over weighed values (menu item 4.9.1) to define 3 classes (menu item 4.81.)

CF

0.0 g

100 %

Delete any existing initialization data

OK

0.0 g

1 %

Begin initialization

0.0 g

1 %

The minus sign of the bar graph flashes to prompt the weight for the upper limit of the first class. A “1” is shown in the numeric display.

10.0 g

1 %

Place the weight for the upper limit of the first class on the weighing instrument (in this example, 10.0 g)

OK

10.0 g

2 %

Store the weight value as the upper limit for the first class

10.0 g

2 %

The middle segments of the bar graph flash to prompt the placement of the weight for the upper limit of the second class on the weighing instrument. A “2” is shown in the numeric display.

20.0 g

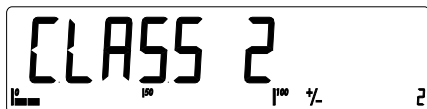
2 %

Place the weight for the upper limit of the second class on the weighing instrument (in this example, 20.0 g)

Operating the Combics 2

OK

Store the weight value as the upper limit for the second class




When initialization is complete, the display shows the class to which the weight on the weighing instrument belongs.

Optional:  < 2 sec.

Toggle to Info mode to view stored values



Upper limit of Class 1 is displayed

Optional: 

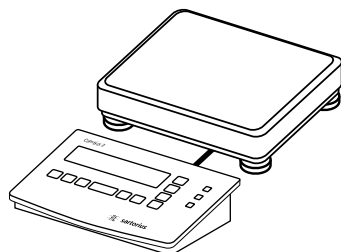
Toggle to the next display mode




Upper limit of Class 2 is displayed




Exit the Info mode



Unload weighing instrument

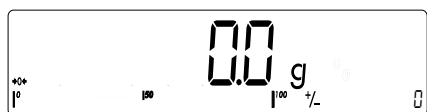

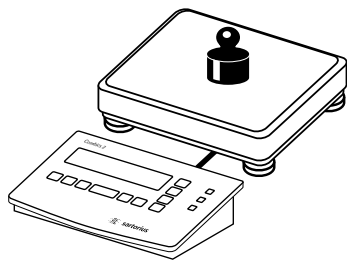Example:
Determine the class of a given weight



Place a load of unknown weight on the weighing instrument



The display shows which class the weight falls into; in this example, Class 2 (10 g – 20 g)



Toggle display to weight readout



The weight of the load is shown (in this example, 17.3 g)



Print results

PET STORE		
GOETTINGEN		
05.03.2002		09:43

Lim1	+	10.0 g
Lim2	+	20.0 g
G#	+	17.3 g
T	+	0.0 g
N	+	17.3 g
Class		2

Printout

Upper limit, Class 1
Upper limit, Class 2

Gross weight
Tare weight
Net weight

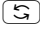
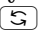
Weight class

Operating the Combics 2


Weighing in Percent %


With the Weighing in Percent application, you can use your weighing instrument to obtain weight readouts in percent which are in proportion to a reference weight. The % symbol is displayed as the weight unit.

Features:


- Store the weight on the weighing instrument as reference weight
- Enter reference weight using a bar code scanner
- Enter tare weight using a bar code scanner
- Toggle the display between percentage and weight by pressing 
- Show percentage as loss or residue
- Display up to 3 decimal places
- Info mode for display of reference weight by pressing  (> 2 sec)


To determine the weight of a sample relative to a reference weight, you must enter the reference weight in the weighing instrument. There are 2 ways to do this with the Combics:

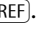- By placing the amount of reference material on the weighing instrument specified by the reference percentage and pressing  to calculate the average and store the value.

The reference percentage is shown in the numeric display, and can be changed by pressing .

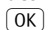How the reference weight is calculated depends on the menu setting for resolution (3.9). Either the value is rounded off in accordance with the display resolution, or the display resolution is increased 10-fold (+1 decimal place) or 100-fold (+2 decimal places), or maximum internal resolution is applied.

- By entering the reference weight (i.e., the weight of 100%) using a bar code scanner (menu setting: "Setup> Bar code > Store value as reference (REF)"). In this case, the value is stored as a reference automatically, without pressing the  key.

This value remains active in the reference memory until you delete it by pressing , overwrite it or until you select a different application. It also remains stored after you switch off the Combics.

Application Parameters: Weighing in Percent

3. 5.	Minimum Load for Automatic Taring and Automatic Printing
3. 5. 1 *	1 digit
3. 5. 2	2 digits
3. 5. 3	5 digits
3. 5. 4	10 digits
3. 5. 5	20 digits
3. 5. 6	50 digits
3. 5. 7	100 digits
3. 5. 8	200 digits
3. 5. 9	500 digits
3. 5. 10	1000 digits
3. 6.	Minimum Load for Initialization
3. 6. 1 *	1 digit
3. 6. 2	2 digits
3. 6. 3	5 digits
3. 6. 4	10 digits
3. 6. 5	20 digits
3. 6. 6	50 digits
3. 6. 7	100 digits
3. 6. 8	200 digits
3. 6. 9	500 digits
3. 6. 10	1000 digits
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3. 8.	Start Application with Most Recent Application Data when Combics is Switched On
3. 8. 1	Automatic (on)
3. 8. 2 *	Manual (off)
3. 9.	Resolution for Calculation of Reference Value
3. 9. 1 *	Display resolution
3. 9. 2	Display resolution +1 decimal place
3. 9. 3	Display resolution +2 decimal places
3. 9. 4	Internal resolution
3. 10.	Decimal Places for Display of Results
3. 10. 1 *	None
3. 10. 2	1 decimal place
3. 10. 3	2 decimal places
3. 10. 4	3 decimal places
3. 11	Storage Parameter
3. 11. 1 *	At stability
3. 11. 2	At increased stability
3. 13.	Reference Weighing Instrument
3. 13. 1 *	No reference instrument selected
3. 13. 2	WP1
3. 13. 3	WP2
3. 15.	Display of calculated value
3. 15. 1 *	Residue
3. 15. 2	Loss

Storage Parameter

The reference weight is saved when the weighing instrument has stabilized. "Stability" is defined as the point at which fluctuation of a measured value lies within a defined tolerance range. The narrower the tolerance range, the more stable the weighing instrument is at stability. This setting is also applied when you zero the weighing instrument.

Under menu item 3.11 you can determine whether the value is saved "At stability" (normal tolerance range) or "At increased stability" (narrower tolerance range).

If you select "At increased stability," the reference weight stored will be more accurate and the results more reproducible, but the response time of the weighing instrument might be longer.

Minimum Load for Initialization

The minimum load for initialization is configured under menu item 3.6.

If the load exceeds this limit, the weighing instrument can be initialized.

If the load is too light, the following will occur when you try to save a value:

- The error code $\text{rEF } 29$ is displayed
- A warning signal is emitted (double-beep)
- The weighing instrument is not initialized
- The preset reference percentage is stored

Weighing in Percent with Two Weighing Instruments

You can use two weighing instruments simultaneously with the Weighing in Percent application. When using two instruments, you can choose from the following operating modes:

- Weighing in Percent with two equivalent weighing instruments
- Weighing in Percent with one reference weighing instrument and one weighing platform

Weighing in Percent with Two Equivalent Weighing Instruments

Use this mode when samples of widely varying weight are measured at one workstation. Measure the lighter-weight pieces on one weighing instrument and the heavier pieces on another.

When you press ALT to toggle from one weighing instrument to the other, the application is re-initialized.

You can define which of the two weighing instruments is active in the display when the Combi is switched on (menu item 8.11). This is the first weighing instrument active when you switch on the Combi, regardless of the setting for automatic initialization of the Weighing in Percent application.

Weighing in Percent with One Reference Weighing Instrument and One Weighing Platform

In this mode, a high resolution weighing instrument with a relatively low maximum capacity is used as a reference weighing instrument. The measuring platform has a high capacity, but a relatively low resolution.

This allows you to both determine the reference value with high resolution; i.e., very precisely, and to measure large samples, without requiring an expensive high-resolution, high-capacity weighing platform.

Specify which weighing instrument is the reference instrument under menu item 3.13. The system toggles automatically to the reference instrument for initialization (rEF is shown in the main display). Following initialization, the system toggles to the measuring platform.

Display

With the Weighing in Percent application, the result can be displayed as a remainder or loss. This is configured in Setup under menu item 3.15.

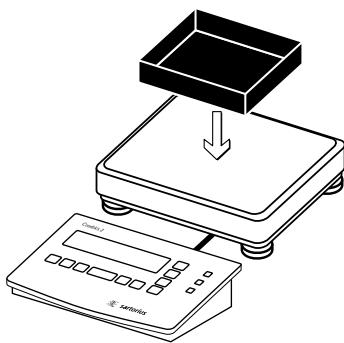
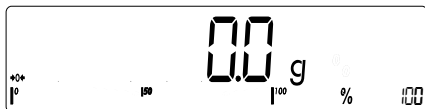
Equations:

$$\text{Residue} = \frac{\text{Current weight}}{100\% \text{ weight}} \cdot 100$$

$$\text{Loss} = \frac{(\text{Current weight} - 100\% \text{ weight})}{100\% \text{ weight}} \cdot 100$$

Operating the Combics 2

Example:



Place empty container on the weighing instrument

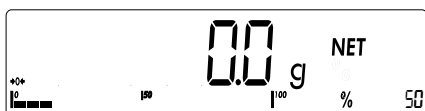


Tare the weighing instrument



REF (repeatedly, if necessary)


Set the desired percentage (in this example, 50%)




Place reference material corresponding to the reference percentage in the container



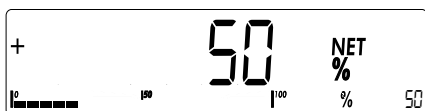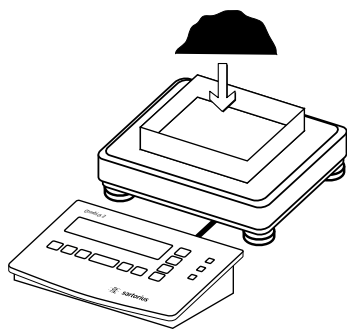Confirm the reference weight





Remove the reference material and place the sample to be measured on the weighing instrument



The percentage of the reference weight equal to the sample is displayed



Print results

G#	+	641.2 g
T	+	200.0 g
N	+	441.2 g
Prc	+	21.00 %

Residue

G#	+	641.2 g
T	+	200.0 g
N	+	441.2 g
D	-	79.00 %

Loss

Operating the Combics 2

Averaging (Animal Weighing)

With the Averaging application, you can use your weighing instrument for calculating weights as the average of a number of individual weighing operations. These individual operations are also known as “subweighing operations.”

This function is used to determine weights under unstable ambient conditions or for weighing unstable samples (such as live animals).

Features:

- Averaging started manually or automatically
- Press the **[REF]** key to select the desired number of subweighing operations
- Toggle the display between last result and current weight by pressing **[S]**

A number of subweighing operations are required to form the basis for calculation of an average weight. These are performed consecutively as soon as the weight of one sample differs from that of the previous sample by less than the preset deviation.

The value for this deviation is entered in the Application settings (3.19) as a percentage. This value and the minimum load must be reached to start the averaging routine.

Whether averaging starts automatically or manually is defined under menu item 3.18.

There are four modes for calculating average weights:

- Manual start with preset number of subweighing operations
To use this mode, place the sample or animal on the weighing instrument and press **[OK]**.
- Manual start with the number of subweighing operations entered manually
To use this mode, place the sample/animal on the weighing instrument, press **[REF]** to enter the number of subweighing operations (10, 20, 30, 40) and then press **[OK]**.


- Automatic start with preset number of subweighing operations
The subweighing operations begin when the weight of the sample/animal on the weighing instrument exceeds the minimum load.
- Automatic start with the number of subweighing operations entered manually
To use this mode, press **[REF]** and enter the number of subweighing operations before placing the sample or animal on the weighing instrument. The subweighing operations begin when the weight of the sample/animal on the weighing instrument exceeds the minimum load.

If the “automatic tare” function is configured, the weight of the first load is stored as the tare value, and averaging begins only when the second load is placed on the weighing instrument (provided this weight exceeds the preset minimum load).

Minimum Load for Automatic Start

The minimum load for automatic start is configured under menu item 3.6. If the load exceeds this limit, the averaging routine can begin.

Display


A calculated average value is shown continuously on the main display. The symbol  is displayed. Press **[S]** to toggle the display between this result and the current weight value.

If the menu setting “Display is static until unload threshold reached” (3.21.1) is active, the program returns to the weight display automatically when the weighing instrument is unloaded; i.e., when the load is less than half the minimum load. The result of the most recent averaging operation is not stored. If the menu setting “Display is static until **[CF]** is pressed” (3.21.2) is active, the average weight remains on the weight display after the weighing instrument is unloaded, until the **[CF]** key is pressed or until a new measurement series is started, whether manually or automatically.

Application Parameters: Averaging (Animal Weighing)

3. 5.	Minimum Load for Automatic Taring and Automatic Printing
3. 5. 1 *	1 digit
3. 5. 2	2 digits
3. 5. 3	5 digits
3. 5. 4	10 digits
3. 5. 5	20 digits
3. 5. 6	50 digits
3. 5. 7	100 digits
3. 5. 8	200 digits
3. 5. 9	500 digits
3. 5. 10	1000 digits
3. 6.	Minimum Load
3. 6. 1 *	1 digit
3. 6. 2	2 digits
3. 6. 3	5 digits
3. 6. 4	10 digits
3. 6. 5	20 digits
3. 6. 6	50 digits
3. 6. 7	100 digits
3. 6. 8	200 digits
3. 6. 9	500 digits
3. 6. 10	1000 digits
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3. 8.	Start Application with Most Recent Application Data when Combics is Switched On
3. 8. 1	Automatic (on)
3. 8. 2 *	Manual (off)
3. 18.	Start of Averaging Routine
3. 18. 1 *	Manual
3. 18. 2	Automatic
3. 19.	Animal Activity
3. 19. 1	0.1 % of the animal/object
3. 19. 2 *	0.2 % of the animal/object
3. 19. 3	0.5 % of the animal/object
3. 19. 4	1 % of the animal/object
3. 19. 5	2 % of the animal/object
3. 19. 6	5 % of the animal/object
3. 19. 7	10 % of the animal/object
3. 19. 8	20 % of the animal/object
3. 19. 9	50 % of the animal/object
3. 19. 10	100 % of the animal/object
3. 20.	Automatic Printout of Results
3. 20. 1 *	Off
3. 20. 2	On
3. 21.	Static Display After Load Removed
3. 21. 1 *	Display is static until unload threshold reached
3. 21. 2	Display is static until [CF] is pressed

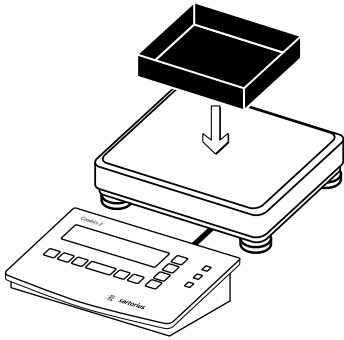Example:



Place empty container on the weighing instrument

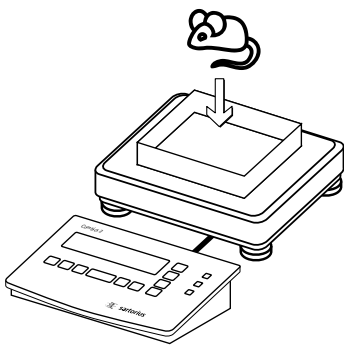


Tare the weighing instrument



REF (repeatedly, if necessary)

Set the desired number of measurements (in this example, 20)




Place 1st animal in bowl



Start measurement



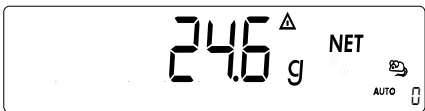Measurement starts as soon as all criteria are met. The number of subweighing operations remaining is shown in the numeric display.




At the conclusion of the series, the calculated mean value is displayed with the \triangle symbol.

Operating the Combics 2




Print results

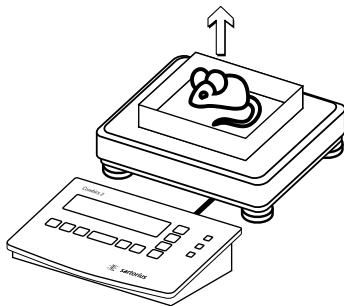x-Net + 24.6 g



Unload weighing instrument

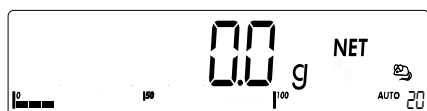



Once the weighing instrument is unloaded
(i.e., the load = less than 1/2 the minimum load), the current weight value is displayed.

Net-Total Formulation



When the Net-Total Formulation application is selected, you can weigh in different components up to a defined total. Each component is saved in the net-total memory.

You can print out both the total weight and the individual weights of the components.

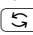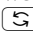Features:

- Weigh in up to 199 components in series
- Current component number shown in the numeric display (indicating the component to be added)
- Toggle the display between “component mode” and “additive mode” by pressing 
 - Component mode: Display the weight of the component currently on the instrument (for 1 second after it is saved; then the instrument is tared)
 - Additive mode: Display the weight of all components on the instrument (after it is saved, the net weight of the last component added is displayed briefly)
- Toggle to a second weighing instrument while weighing on the first
- Automatic printout of component or total weight when weight is saved
- Automatic taring of container weight
- Minimum load for automatic taring
- Minimum load for storing values
- Display the value in totalizing memory in the Info mode, accurately calculated in the active weight unit, by pressing and holding 

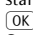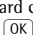Application Parameters: Net-Total Formulation



3. 5.	Minimum Load for Automatic Taring and Automatic Printing
3. 5. 1 *	1 digit
3. 5. 2	2 digits
3. 5. 3	5 digits
3. 5. 4	10 digits
3. 5. 5	20 digits
3. 5. 6	50 digits
3. 5. 7	100 digits
3. 5. 8	200 digits
3. 5. 9	500 digits
3. 5.10	1000 digits
3. 6.	Minimum Load for Automatic Storage/Transfer of Values
3. 6. 1 *	1 digit
3. 6. 2	2 digits
3. 6. 3	5 digits
3. 6. 4	10 digits
3. 6. 5	20 digits
3. 6. 6	50 digits
3. 6. 7	100 digits
3. 6. 8	200 digits
3. 6. 9	500 digits
3. 6.10	1000 digits
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3.17.	Automatic Individual or Component Printout when Value Stored
3.17. 1	Off
3.17. 2*	Generate printout with complete standard configuration each time  is pressed
3.17. 3	Generate printout with complete standard configuration only once when  is pressed

Minimum Load for Automatic Storage/Transfer of Values

The minimum load for automatic storage or transfer of weight values is configured under menu item 3.6. If the load exceeds this limit, the weight currently displayed can be stored in net-total memory.

Net-Total Formulation with Two Weighing Instruments

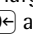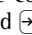
This mode is used for weighing large and small components at the same time.

In this mode, you can toggle from the small-component instrument to the large-component instrument during measurement. Once you toggle to the large-component instrument, the  and  keys are available until a component is value is saved. For example, you can tare a partially-filled container taken from the small-component instrument on the large component instrument.

The value in component memory on the small-component instrument is transferred to the large-component instrument and the weight unit is converted, if necessary. The Component and Additive display modes are both available on the large-component instrument.

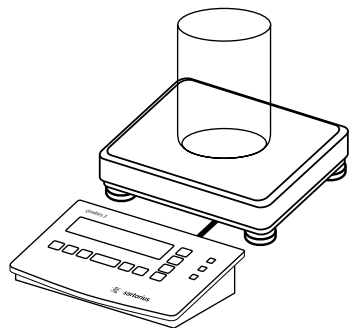
The value read by the active instrument is saved in component memory. The displayed result is accurately calculated in the active weight unit, for example:

1243 g
(from an instrument with 1 interval)
+ 1.40 kg
(from an instrument with 5 intervals)
= 2.643 kg (calculated result)


When you press  to stop a measurement series, the tare memories for both weighing instruments are cleared, unless the large-component instrument is in SBI mode, in which case the instrument is only tared.

Operating the Combics 2

Example:
Weigh in 3 components, display total weight after 2nd component (Additive mode)




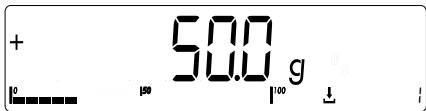
Place empty container on the weighing instrument




Container weight a prompt for first component are shown



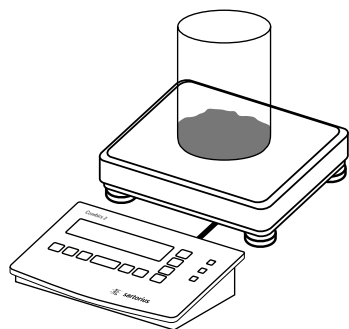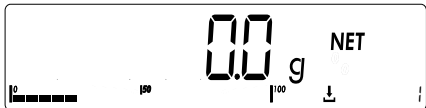Tare the weighing instrument



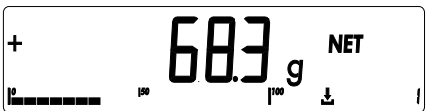Add the first component (in this example, 68.3 g)



The weight of the first component is displayed



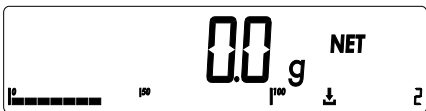
Store the weight of the first component

ACE HARDWARE
GOETTINGEN
20.02.2002 15:10


Cmp001 68.3 g

Total data record is printed, as configured under menu item 7.7. or 7.8.
The header is printed only once

Component weight



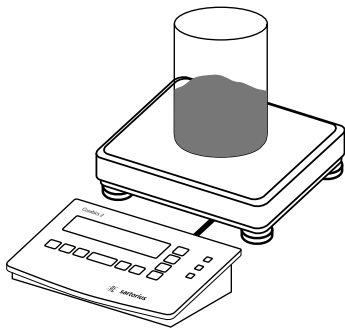The weighing instrument is tared automatically and the component counter value is increased by one (prompting the next component)



Add the second component (in this example, 82.7 g)



The weight of the second component is displayed



Store the weight of the second component

Cmp002 82.7 g

Component weight is printed



The weighing instrument is tared automatically and the component counter value is increased by one (prompting the next component)



Toggle to the Additive mode



The value displayed equals the weight of all components added up to now plus the current weight



Add the third component (bringing the total to the desired target)



Store the weight of the third component

Cmp003 49.1 g

Component weight is printed



End the weighing-in operation

Total data record is printed, as configured under 7.9. or 7.10.

```
Ser.no      12345678
n           3
S-Comp+    200.1 g
Cont.T+     50.0 g
-----
```

Serial number of the weighing instrument
 Number of components
 Content of component memory
 Content of tare memory (container weight)
 Dotted line

Operating the Combics 2

Totalizing Σ

With the Totalizing application, you can add weight values to the totalizing memory. The number of values added to the memory is also stored ("transaction counter").

Features:

- Weigh in up to 199 items
- Current transaction number displayed in the text line (indicating the items already added)
- Save weight values manually or automatically
- Accurate calculation of total of weight values from two weighing instruments
- Unfinished totalization stored in battery-backed memory
- Toggle the display between current weight and total value by pressing (Σ):
 - Totalizing mode: shows the value in totalizing memory plus the current weight; for filling to a defined total
- Automatic printout of component or total weight when weight is saved
- Automatic taring of container weight
- Minimum load for automatic tare
- Minimum load for storing values
- Press and hold (Σ) to display the value in totalizing memory (net, gross) in the Info mode, accurately calculated in the active weight unit.

The Combics has a totalizing memory for adding individual net and gross values. There are two ways to store weight values in the totalizing memory:

- Manually, by pressing the (OK) key. The value is added to the value already stored, and the transaction counter value increases by one. When a value is added manually, the program does not check whether the weighing instrument has been unloaded since the last time the (OK) key was pressed.
- Automatically, when the weighing instrument has stabilized and the minimum load is exceeded (menu item 3.16.2). The weighing instrument is considered to be unloaded when the load is less than 50% of the minimum load (menu item 3.6).

The numeric display shows the transaction counter.

Press the (CF) key to clear the totalizing memory. A printout is generated automatically.

With the weighing instruments connected, you can add values from both instruments to the totalizing memory. The displayed result is accurately calculated in the active weight unit. Example: 1.243 kg (measured on a weighing instrument using 1 scale interval) added to 1.40 kg (measured on a weighing instrument using 5 scale intervals) is displayed as 2.643 kg.

Application Parameters: Totalizing

- 3. 5. Minimum Load for Automatic Taring and Automatic Printing
 - 3. 5. 1 * 1 digit
 - 3. 5. 2 2 digits
 - 3. 5. 3 5 digits
 - 3. 5. 4 10 digits
 - 3. 5. 5 20 digits
 - 3. 5. 6 50 digits
 - 3. 5. 7 100 digits
 - 3. 5. 8 200 digits
 - 3. 5. 9 500 digits
 - 3. 5. 10 1000 digits
- 3. 6. Minimum Load for Automatic Storage/Transfer of Values
 - 3. 6. 1 * 1 digit
 - 3. 6. 2 2 digits
 - 3. 6. 3 5 digits
 - 3. 6. 4 10 digits
 - 3. 6. 5 20 digits
 - 3. 6. 6 50 digits
 - 3. 6. 7 100 digits
 - 3. 6. 8 200 digits
 - 3. 6. 9 500 digits
 - 3. 6. 10 1000 digits
- 3. 7. Automatic Taring: 1st Weight Tared
 - 3. 7. 1 * Off
 - 3. 7. 2 On
- 3. 8. Start Application with Most Recent Application Data when Combics is Switched On
 - 3. 8. 1 Automatic (on)
 - 3. 8. 2 * Manual (off)
- 3.16. Automatic Storage/Transfer of Value
 - 3.16. 1* Off
 - 3.16. 2 On
- 3.17. Automatic Individual or Component Printout when Value Stored
 - 3.17. 1 Off
 - 3.17. 2* Generate printout with complete standard configuration each time (OK) is pressed
 - 3.17. 3 Generate printout with complete standard configuration only once when (OK) is pressed

Minimum Load for Automatic Storage/Transfer of Values

The minimum load for automatic storage or transfer of weight values is configured under menu item 3.6. Once the load on the weighing instrument exceeds the specified minimum, the measured value is stored or transferred automatically, if the operating menu is configured for automatic storage (menu item 3.16.2).

Display

Press **[S]** to toggle the display between the weight currently on the weighing instrument and the current weight value plus the value in totalizing memory.

Press and hold **[S]** (> 2 sec) to toggle to the Info mode. In the "Info" mode, the total of gross weight values is shown first. Press **[S]** again to show the net value.

Press the **[S]** key again to show the transaction counter in the main display. Press **[S]** one more time to exit the Info mode.

To exit the "Info" mode earlier, press and hold **[S]** (> 2 sec).

Printout

Under menu item 3.17 you can configure whether a printout is generated manually, by pressing **[E]** (3.17.3), or automatically when a weight value is stored in the totalizing memory (3.17.2 or 3.17.3).

When you generate a printout manually by pressing **[E]**, the transaction counter value is not printed.

For automatic printouts, you can define whether a printout is generated after each individual transaction (3.17.2) or includes all totalized components (3.17.3).

Each printout starts with the header lines and the date and time, followed by a dotted line, before printing the gross, tare, net and transaction counter values.

The printout of components includes the header lines and date and time only on the first printout. Subsequent printouts include only the gross, tare, net and transaction counter values (and serial number, if configured).

The individual and component printouts are configured under menu items 7.7 and 7.8.

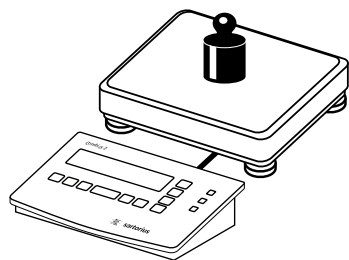
The total data record is printed when you clear the totalizing memory (by pressing the **[CF]** key).

The total data record includes data from the gross totalizing memory (*** G**), the net totalizing memory (*** N**), the total transaction count (**n**) and a dotted line.

The fields selected under menu items 7.9. and 7.10. are printed every time the **[CF]** key is pressed, regardless of whether individual or component printing is configured.

Operating the Combics 2

Example:
Totalizing weight values with printout of components



1. Place the sample on the weighing instrument



Weight value is displayed




Store weight value in totalizing memory




Main display: Current weight on the weighing instrument
Numeric display: Transaction counter (in this example: 1)

SARTORIUS			
GOETTINGEN			
06.02.2002		12:15	


G#	+	102.9	g
T	+	0.0	g
N	+	102.9	g
n		1	

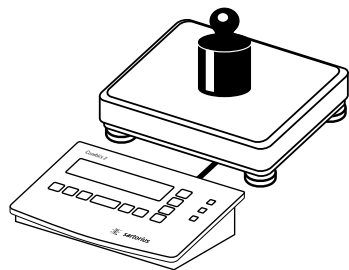

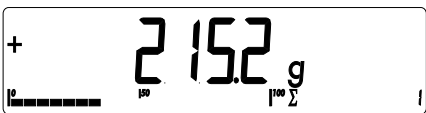
Unload the weighing instrument and then place the second sample on the instrument



Weight value is displayed

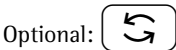

Toggle to Totalizing mode

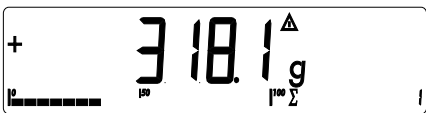


Main display: Totalizing memory content plus current weight
Numeric display: Transaction counter (in this example: 1)



Store weight value in totalizing memory

SARTORIUS
GOETTINGEN
06.02.2002 12:19

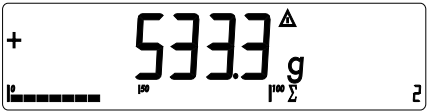G# + 215.2 g
T + 0.0 g
N + 215.2 g
n 2



Main display: Totalizing memory content plus current weight
Numeric display: Transaction counter (in this example: 2)

CF

End totalizing

Total data record is printed, as configured under 7.9. or 7.10.

*G + 318.1 g
*N + 318.1 g
n 2

Number of components
Contents of component memory
Contents of tare memory (container weight)
Dotted line

Configuration

Operating Menu Overview

You can configure the Combics to meet individual requirements by entering user data and setting selected parameters in the operating menu.

The operating menu is a combination of text levels and numeric codes.

- ////////// = Setting/function available on Combics 1 plus only
- ===== = Setting/function available on Combics 2 only
- ===== = Setting/function available on Combics 1 plus and Combics 2 only

First level display

2nd level display

Function

Menu

APPL			Select and configure application programs
		ΔΔ	Basic Weighing Function
		•••	Counting
		••• nπ	Neutral Measurement
		☞	Averaging (Animal Weighing)
		+/-	Checkweighing
		+/- CL	Classification
		%	Weighing in Percent
		↓	Net-Total Formulation
		Σ	Totalizing
Fn-KEY			Define the function of Fn
		oFF	No function
		Gross net	Gross/net toggling
		2. unit	Show 2nd weight unit
SETUP			Adapt Combics to user requirements
		WP1	Settings for weighing instrument on WP1
		COM1	Settings for the RS-232 interface
		unICoM	Settings for the optional second interface
		Ctrl In	Set the function of the universal input (control line)
		barCodE	Set the bar code scanner function
		PrntPro	Settings for the printout
		ut IL It	Settings for additional functions
		tIME	Set the time
		dAtE	Set the date
Info		Code	Enter a password to protect menu settings
			View device-specific information (service date, serial number, etc.)
LANG			Select language for calibration, adjustment and GMP printouts
		dEutSCH	German
		EngLISH	English
		US ModE	English with U.S. date/time format
		FRANc	French
		ItAL	Italian
		ESPAñOL	Spanish

Example: Change language from “English” to “English with U.S. date/time format”



Switch on the Combics



While all segments are lit, press the key

First level of menu for selecting the application program *APPL* is shown

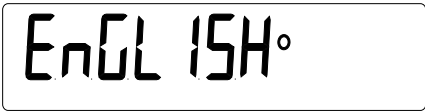


Select the *LANG* menu item (press repeatedly until *LANG* is displayed)



Select *LANG* device parameter

The currently active language is shown



Select the *U.S. Mode* menu item

Press repeatedly until *U.S. Mode* is displayed)



Confirm this menu item

Return to menu level 2 to configure other menu settings, or



return to menu level 1 to configure other menu settings, or



exit the operating menu

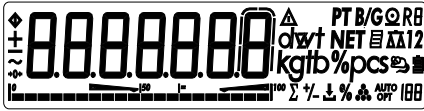


Configuration


Example: Assign a password to protect the application program settings “APPL” and the device parameters “SEtUP” from unauthorized changes (in this example: AB2)



Switch on the Combics



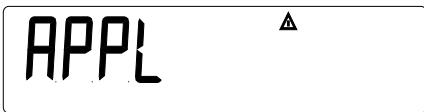
While all segments are lit, press the key




First level of menu for selecting the application program *APPL* is shown



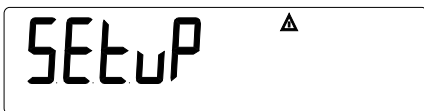
Select the *SEtUP* menu item
(press repeatedly until *SEtUP* is displayed)





Open the *SEtUP* menu

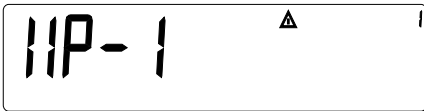

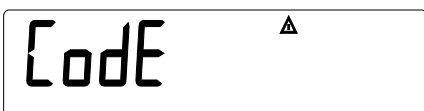Select the *Code* menu item
(Press repeatedly until *Code* is displayed)




Open the password input menu item

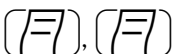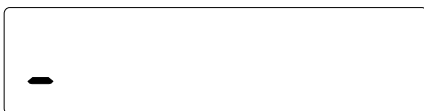

Enter the first character using the and keys (in this example: *A*)



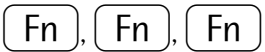Store character



Enter the second character using the and keys (in this example: b)

Store character

Enter the third character using the and keys (in this example: 2)

Store character

To enter other characters, use the and keys and press to confirm

Store the password

Exit this menu level to configure other menu settings, or

Exit the operating menu

To delete a password:
Overwrite the old password with the new password, or enter a space as the password and press to confirm.

Configuration

Operating Menu Overview (Parameters)



- ////// = Setting/function available on Combics 1 plus only
- ===== = Setting/function available on Combics 2 only
- ===== = Setting/function available on Combics 1 plus and Combics 2 only

* Factory setting

Menu

<div>APPL</div> <div><div><div><div>3. 5.</div><div>3. 5. 1 *</div><div>3. 5. 2</div><div>3. 5. 3</div><div>3. 5. 4</div><div>3. 5. 5</div><div>3. 5. 6</div><div>3. 5. 7</div><div>3. 5. 8</div><div>3. 5. 9</div><div>3. 5. 10</div></div><div><div>3. 7.</div><div>3. 7. 1 *</div><div>3. 7. 2</div></div><div><div>9. 1.</div><div>9. 1. 1</div><div>9. 1. 2 *</div></div></div><div><div><div>3. 5.</div><div>Numeric menu as for Weighing</div></div><div><div>3. 6.</div><div>3. 6. 1 *</div><div>3. 6. 2</div><div>3. 6. 3</div><div>3. 6. 4</div><div>3. 6. 5</div><div>3. 6. 6</div><div>3. 6. 7</div><div>3. 6. 8</div><div>3. 6. 9</div><div>3. 6. 10</div></div><div><div>3. 7.</div><div>3. 7. 1 *</div><div>3. 7. 2</div></div><div><div>3. 8.</div><div>3. 8. 1</div><div>3. 8. 2 *</div></div><div><div>3. 9.</div><div>3. 9. 1 *</div><div>3. 9. 2</div><div>3. 9. 3</div><div>3. 9. 4</div></div><div><div>3. 11.</div><div>3. 11. 1 *</div><div>3. 11. 2</div></div><div><div>3. 12.</div><div>3. 12. 1</div><div>3. 12. 3 *</div></div><div><div>3. 13.</div><div>3. 13. 1 *</div><div>3. 13. 2</div><div>3. 13. 3</div></div><div><div>9. 1.</div><div>9. 1. 1</div><div>9. 1. 2 *</div></div></div></div>	<div><h3>Application Programs</h3><h4>Basic Weighing Function</h4><p>Minimum Load for Automatic Taring and Automatic Printing</p><p>1 digit</p><p>2 digits</p><p>5 digits</p><p>10 digits</p><p>20 digits</p><p>50 digits</p><p>100 digits</p><p>200 digits</p><p>500 digits</p><p>1000 digits</p><p>Automatic Taring: 1st Weight Tared</p><p>Off</p><p>On</p><p>Factory Settings for All Application Programs</p><p>Yes</p><p>No</p></div> <div><h3>Counting</h3><p>Minimum Load for Automatic Taring and Automatic Printing</p><p>Minimum Load for Initialization</p><p>1 digit</p><p>2 digits</p><p>5 digits</p><p>10 digits</p><p>20 digits</p><p>50 digits</p><p>100 digits</p><p>200 digits</p><p>500 digits</p><p>1000 digits</p><p>Automatic Taring: 1st Weight Tared</p><p>Off</p><p>On</p><p>Start Application with Most Recent Application Data when Combics is Switched On</p><p>Automatic (on)</p><p>Manual (off)</p><p>Resolution for Calculation of Reference Value</p><p>Display resolution</p><p>Display resolution + 1 decimal place</p><p>Display resolution + 2 decimal places</p><p>Internal resolution</p><p>Storage Parameter</p><p>At stability</p><p>At increased stability</p><p>Reference Sample Updating</p><p>Off</p><p>Automatic</p><p>Reference Weighing Instrument</p><p>No reference instrument selected</p><p>WP1</p><p>WP2</p><p>Factory Settings for All Application Programs</p><p>Yes</p><p>No</p></div>
--	--

Configuration

APPL	
	Neutral Measurement
3. 5.	Minimum Load for Automatic Taring and Automatic Printing
3. 6.	Minimum Load for Initialization
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3. 8.	Start Application with Most Recent Application Data when Combics is Switched On
3. 8. 1	Automatic (on)
3. 8. 2 *	Manual (off)
3. 9.	Resolution for Calculation of Reference Value
3. 9. 1 *	Display resolution
3. 9. 2	Display resolution + 1 decimal place
3. 9. 3	Display resolution + 2 decimal places
3. 9. 4	Internal resolution
3.10.	Decimal Places for Display of Results
3.10. 1 *	None
3.10. 2	1 decimal place
3.10. 3	2 decimal places
3.10. 4	3 decimal places
3.11.	Storage Parameter
3.11. 1 *	At stability
3.11. 2	At increased stability
3.13.	Reference Weighing Instrument
3.13. 1 *	Off
3.13. 2	To weighing instrument WP1
3.13. 3	To weighing instrument WP2
9. 1.	Factory Settings for All Application Programs
9. 1. 1	Yes
9. 1. 2 *	No
	Averaging (Animal Weighing)
3. 5.	Minimum Load for Automatic Taring and Automatic Printing
3. 6.	Minimum Load for Automatic Start
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3. 8.	Start Application with Most Recent Application Data when Combics is Switched On
3. 8. 1	Automatic (on)
3. 8. 2 *	Manual (off)
3.18.	Start of Averaging Routine
3.18. 1 *	Manual
3.18. 2	Automatic
3.19.	Animal Activity
3.19. 1	0.1 % of the animal/object
3.19. 2 *	0.2% of the animal/object
3.19. 3	0.5 % of the animal/object
3.19. 4	1% of the animal/object
3.19. 5	2% of the animal/object
3.19. 6	5% of the animal/object
3.19. 7	10% of the animal/object
3.19. 8	20% of the animal/object
3.19. 9	50% of the animal/object
3.19.10	100 % of the animal/object
3.20.	Automatic Printout of Results
3.20. 1 *	Off
3.20. 2	On

APPL



- 3.21.
 - 3.21. 1 *
 - 3.21. 2
- 9. 1.
 - 9. 1. 1
 - 9. 1. 2 *

Static Display After Load Removed
 Display is static until unload threshold reached
 Display is static until **CF** is pressed

Factory Settings for All Application Programs

Yes

No



- 3. 5.
 - Numeric menu as for Weighing
- 3. 7.
 - 3. 7. 1 *
 - 3. 7. 2
- 3. 8.
 - 3. 8. 1
 - 3. 8. 2 *
- 4. 2.
 - 4. 2. 1 *
 - 4. 2. 2
- 4. 3.
 - 4. 3. 1 *
 - 4. 3. 2
- 4. 4.
 - 4. 4. 1
 - 4. 4. 2
 - 4. 4. 3
 - 4. 4. 4 *
 - 4. 4. 5
- 4. 5.
 - 4. 5. 1 *
 - 4. 5. 2
- 4. 6.
 - 4. 6. 1 *
 - 4. 6. 2
 - 4. 6. 3
 - 4. 6. 4
- 9. 1.
 - 9. 1. 1
 - 9. 1. 2 *

Checkweighing

Minimum Load for Automatic Taring and Automatic Printing

Automatic Taring: 1st Weight Tared

Off

On

Start Application with Most Recent Application Data when Combics is Switched On

Automatic (on)

Manual (off)

Check Range

30% to 170%

10% to infinity

Activate Control Line for "Set" as:

"Set" output

Ready to operate (for process control systems)

Activation of Outputs

Off

Always active

Active at stability

Active within check range

Active at stability within the check range

Parameter Input

Min, max, target

Only target with percent limits

Automatic Printing

Off

On

Only values within tolerance

Only values outside tolerance

Factory Settings for All Application Programs

Yes

No



- 3. 5.
 - Numeric menu as for weighing
- 3. 6.
 - Numeric menu as for Counting
- 3. 7.
 - 3. 7. 1 *
 - 3. 7. 2
- 3. 8.
 - 3. 8. 1
 - 3. 8. 2 *
- 4. 3.
 - 4. 3. 1 *
 - 4. 3. 2
- 4. 7.
 - 4. 7. 1
 - 4. 7. 2
 - 4. 7. 3 *

Classification

Minimum Load for Automatic Taring and Automatic Printing

Minimum Load for Initialization and Definition of "Class 1" Lower Limit

Automatic Taring: 1st Weight Tared

Off

On

Start Application with Most Recent Application Data when Combics is Switched On

Automatic (on)

Manual (off)

Activate Control Line for "Set" as:

"Set" output

Ready to operate

Activation of Outputs

Off


Always active

Active at stability

Configuration

APPL	
% CL	
4. 8.	Number of Classes
4. 8. 1 *	3 classes
4. 8. 2	5 classes
4. 9.	Parameter Input
4. 9. 1 *	Weight values
4. 9. 2	Percentage
4.10.	Automatic Printing
4.10. 1 *	Off
4.10. 2	On
9. 1.	Factory Settings for All Application Programs
9. 1. 1	Yes
9. 1. 2 *	No
%	Weighing in Percent
3. 5.	Minimum Load Automatic Taring and Automatic Printing
Numeric menu as for Weighing	
3. 6.	Minimum Load for Automatic Initialization
Numeric menu as for Counting	
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On
3. 8.	Start Application with Most Recent Application Data when Combics is Switched On
3. 8. 1	Automatic (on)
3. 8. 2 *	Manual (off)
3. 9.	Resolution for Calculation of Reference Value
3. 9. 1 *	Display resolution
3. 9. 2	Display resolution + 1 decimal place
3. 9. 3	Display resolution + 2 decimal places
3. 9. 4	Internal resolution
3.10.	Decimal Places for Display of Results
3.10. 1 *	None
3.10. 2	1 decimal place
3.10. 3	2 decimal places
3.10. 4	3 decimal places
3.11.	Storage Parameter
3.11. 1 *	At stability
3.11. 2	At increased stability
3.13.	Reference Weighing Instrument
3.13. 1 *	Off
3.13. 2	To weighing instrument WP1
3.13. 3	To weighing instrument WP2
3.15.	Display of Calculated Value
3.15. 1 *	Residue
3.15. 2	Loss
9. 1.	Factory Settings for All Application Programs
9. 1. 1	Yes
9. 1. 2 *	No
↓	Net-Total Formulation (2nd Tare Memory)
3. 5.	Minimum Load Automatic Taring and Automatic Printing
Numeric menu as for Weighing	
3. 6.	Minimum Load for Automatic Storage/Transfer of Values
Numeric menu as for Counting	
3. 7.	Automatic Taring: 1st Weight Tared
3. 7. 1 *	Off
3. 7. 2	On

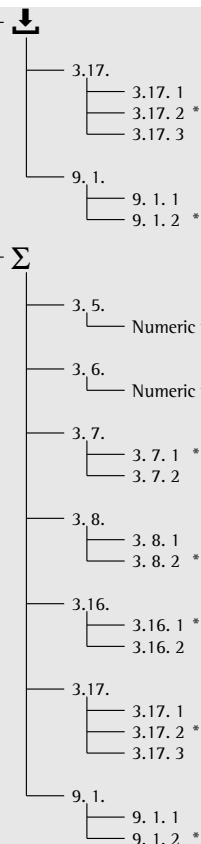APPL



Printout when Value is Stored in Totalizing Memory
Automatic printout of results: Off
Generate printout with complete standard configuration each time **OK** is pressed
Generate printout with complete standard configuration only once when **OK** is pressed

Factory Settings for All Application Programs

Yes
No

Totalizing

Minimum Load Automatic Taring and Automatic Printing

Minimum Load for Automatic Storage/Transfer of Values

Automatic Taring: 1st Weight Tared

Off
On

Start Application with Most Recent Application Data when Combits is Switched On

Automatic (on)
Manual (off)

Automatic Storage

Off
On

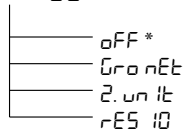
Printout when Value is Stored in Totalizing Memory

Automatic printout of results: Off
Individual of transaction by pressing **OK**
Print components of transaction by pressing **OK**

Factory Settings for All Application Programs

Yes
No

Fn-Fn



Key Assignment: **Fn**

No **Fn** key function

Gross/net toggling

Show 2nd unit

10-fold increased resolution. Display: max. 10 seconds

Configuration

SEtUP		Device Parameters	
11P-1		Password prompt displayed, if a password is configured	
OFF		Weighing Instrument 1	
ON		(Display designation of this menu level: 1)	
		Off	
		On	
1. 1.		Adapting Weighing Instrument to Ambient Conditions (Adapt Filter)	
1. 1. 1		Very stable conditions	
1. 1. 2		Stable conditions	
1. 1. 3		Unstable conditions	
1. 1. 4		Very unstable conditions	
1. 2.		Application Filter	
1. 2. 1		Final readout	
1. 2. 2		Filling mode	
1. 2. 3		Low filtering	
1. 2. 4		Without filtering	
1. 3.		Stability Range	
1. 3. 1		¼ digit	
1. 3. 2		½ digit	
1. 3. 3		1 digit ¹⁾	
1. 3. 4		2 digits ¹⁾	
1. 3. 5		4 digits ¹⁾	
1. 3. 6		8 digits ¹⁾	
1. 4.		Stability delay	
1. 4. 1		No delay	
1. 4. 2		Short delay	
1. 4. 3		Average delay	
1. 4. 4		Long delay	
1. 5.		Taring ¹⁾	
1. 5. 1		Without stability	
1. 5. 2		After stability	
1. 6.		Auto Zero	
1. 6. 1		On	
1. 6. 2		Off	
1. 7.		Weight Unit 1 ²⁾	
1. 7. 2		Grams /g	
1. 7. 3		Kilograms /kg	
1. 7. 4		Carats /ct ¹⁾	
1. 7. 5		Pounds /lb ¹⁾	
1. 7. 6		Ounces /oz ¹⁾	
1. 7. 7		Troy ounces /ozt ¹⁾	
1. 7. 8		Hong Kong taels /tlh ¹⁾	
1. 7. 9		Singapore taels /tls ¹⁾	
1. 7.10		Taiwanese taels /tlt ¹⁾	
1. 7.11		Grains /GN ¹⁾	
1. 7.12		Pennyweights /dwt ¹⁾	
1. 7.13		Milligrams /mg ¹⁾	
1. 7.14		Parts per pound /lb ¹⁾	
1. 7.15		Chinese taels /tlc ¹⁾	
1. 7.16		Mommies /mom ¹⁾	
1. 7.17		Austrian carats /k ¹⁾	
1. 7.18		Tola /tol ¹⁾	
1. 7.19		Baht /bat ¹⁾	
1. 7.20		Mesghal /MS ¹⁾	
1. 7.21		Tons /t	
1. 8.		Display Accuracy 1	
1. 8. 1		All digits	
1. 8. 2		Reduced by 1 decimal place for load change	
1. 8. 3		Reduced by 1 decimal place ¹⁾	
1. 8. 4		Reduced by 2 decimal places ¹⁾	
1. 8. 5		Reduced by 3 decimal places ¹⁾	

¹⁾ = Not available on instruments verified for use in legal metrology
²⁾ = Depends on weighing instrument model


SETUP		WP-1	
	on		
1. 9.			Calibration, Adjustment
1. 9. 1	*		Ext. calibration/adjustment; default weight
1. 9. 3			Ext. calibration/adjustment; weight can be selected under menu item 1.18.1
1. 9.10			Press $\rightarrow T \leftarrow$ -> 2 sec: blocked
1.10.			Calibration/adjustment sequence
1.10. 1			Calibration with automatic adjustment
1.10. 2	*		Calibration with adjustment triggered manually
1.11.			Zero-setting Range
1.11. 1			1 percent/max. cap.
1.11. 2	*		2 percent/max.cap.
1.12.			Initial Zero-setting Range
1.12. 2			2 percent/max. cap.
1.12. 3	*		5 percent/max.cap.
1.13.			Tare/Zero at Power On
1.13. 1	*		On
1.13. 2			Off, load previous tare value
1.13. 3			Only zero at power on
1.15.			Calibration Prompt
1.15. 1	*		Off
1.15. 2			Calibration prompt "ΔΔ" flashes in the display
1.16.			External Calibration/Adjustment ¹⁾
1.16. 1	*		Accessible
1.16. 2	²⁾		Blocked
1.18.			Enter Calibration Weight
1.18. 1			External user-defined weight; enter value; e.g.: 10.000 kg)
3. 1.			Weight Unit ^{2 3)}
3. 1. 2			Grams /g
3. 1. 3			Kilograms /kg
3. 1. 4			Carats /ct ¹⁾
3. 1. 5			Pounds /lb ¹⁾
3. 1. 6			Ounces /oz ¹⁾
3. 1. 7			Troy ounces /ozt ¹⁾
3. 1. 8			Hong Kong tael /tlh ¹⁾
3. 1. 9			Singapore tael /tls ¹⁾
3. 1.10			Taiwanese tael /tlt ¹⁾
3. 1.11			Grains /GN ¹⁾
3. 1.12			Pennyweights /dwt ¹⁾
3. 1.13			Milligrams /mg ¹⁾
3. 1.14			Parts per pound /lb ¹⁾
3. 1.15			Chinese tael /tlc ¹⁾
3. 1.16			Mommes /mom ¹⁾
3. 1.17			Austrian carats /k ¹⁾
3. 1.18			Tola /tol ¹⁾
3. 1.19			Baht /bat ¹⁾
3. 1.20			Mesghal /MS ¹⁾
3. 1.21			Tons /t
3. 2.			Display Accuracy 2
3. 2. 1	*		All digits
3. 2. 2			Reduced by 1 decimal place for load change
3. 2. 3			Reduced by 1 decimal place ¹⁾
3. 2. 4			Reduced by 2 decimal places ¹⁾
3. 2. 5			Reduced by 3 decimal places ¹⁾
9. 1.			Restore Factory Settings for WP1 Numeric Menu
9. 1. 1			Yes
9. 1. 2	*		No

¹⁾ = Not available on instruments verified for use in legal metrology

²⁾ = Factory setting on instrument verified for use in legal metrology

³⁾ = Menu depends on weighing instrument model

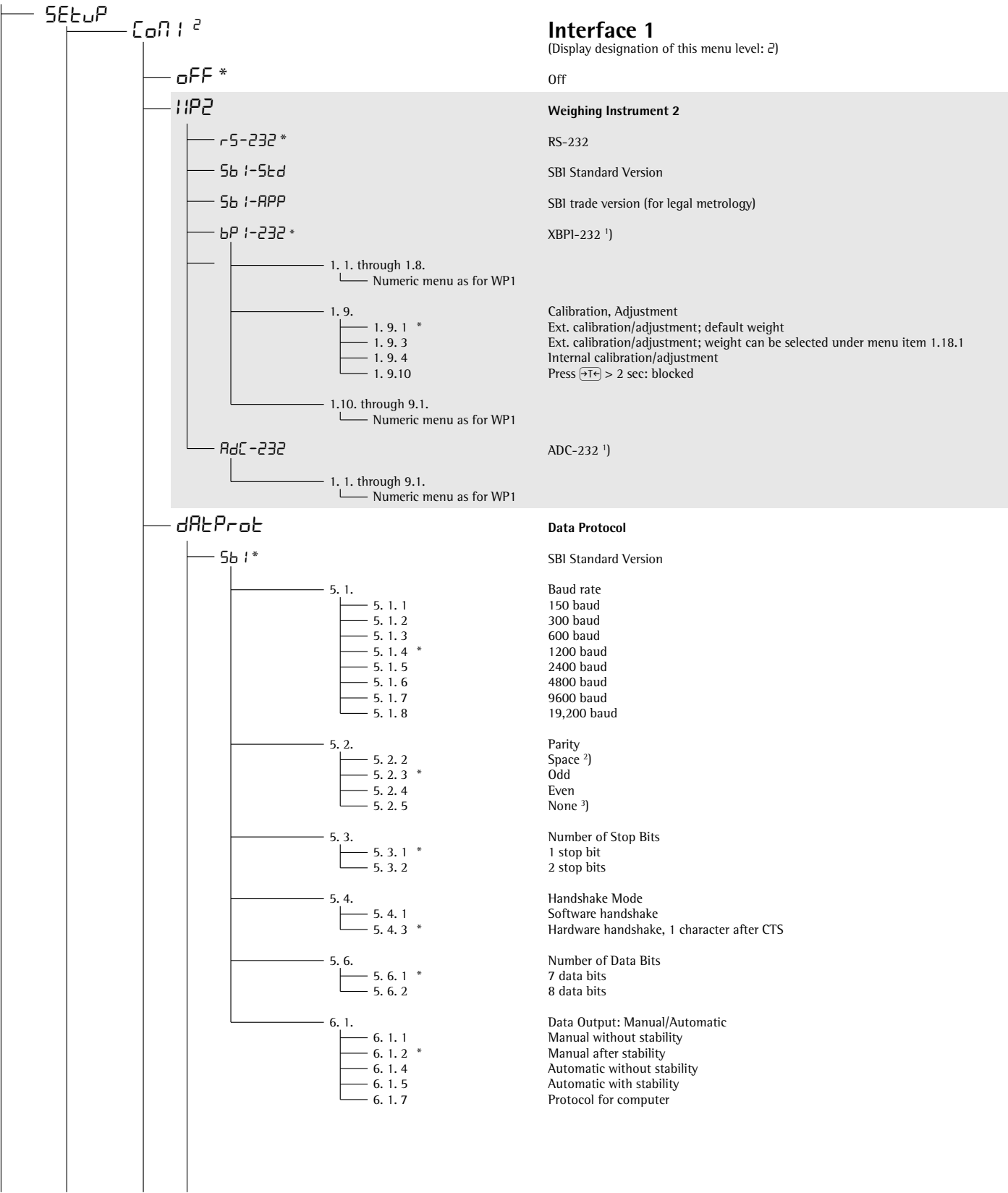
Configuration




¹⁾ = Menu depends on weighing instrument model
²⁾ = not with setting 5.6.2 (8 bits)
³⁾ = Not with setting 5.6.1 (7 bits)

SEtUP	
Con1 2	
dAtP rot	
Sb1 *	
6. 3.	Time-dependent automatic data output
6. 3. 1 *	1 display update
6. 3. 2	2 display updates
6. 3. 4	10 display updates
6. 3. 7	100 display updates
7. 2.	Data Output: Line Format
7. 2. 1	For raw data: 16 characters
7. 2. 2 *	For other applications: 22 characters
9. 1.	Restore Factory Settings for Numeric Menu COM1: SBI
9. 1. 1	Yes
9. 1. 2 *	No
bP1-232	XBPI-232
MP8	MP8 Interface Emulation
ProG	Program selection
P 3-1-1	MP8: 3-1-1
...	
P 3-1-9	MP8: 3-1-9
P 3-2-1	MP8: 3-2-1
...	
P 3-2-9	MP8: 3-2-9
P 3-3-1	MP8: 3-3-1
...	
P 3-3-9	MP8: 3-3-9
Ind. 2	Index 2
1 2. 1	Ind. 2.1
1 2. 2	Ind. 2.2
1 2. 3	Ind. 2.3
1 2. 4	Ind. 2.4
bAud	Baud Rate
5. 1. 1	150 baud
5. 1. 2	300 baud
5. 1. 3	600 baud
5. 1. 4 *	1200 baud
5. 1. 5	2400 baud
5. 1. 6	4800 baud
5. 1. 7	9600 baud
PAR 1	Parity
5. 2. 1	Mark
5. 2. 2	Space
5. 2. 3 *	Odd
5. 2. 4	Even
Pr int	Manual/Auto Print Mode
6. 1. 1	Manual without stability
6. 1. 2 *	Manual after stability
6. 1. 4	Automatic without stability
6. 1. 5	Automatic at stability

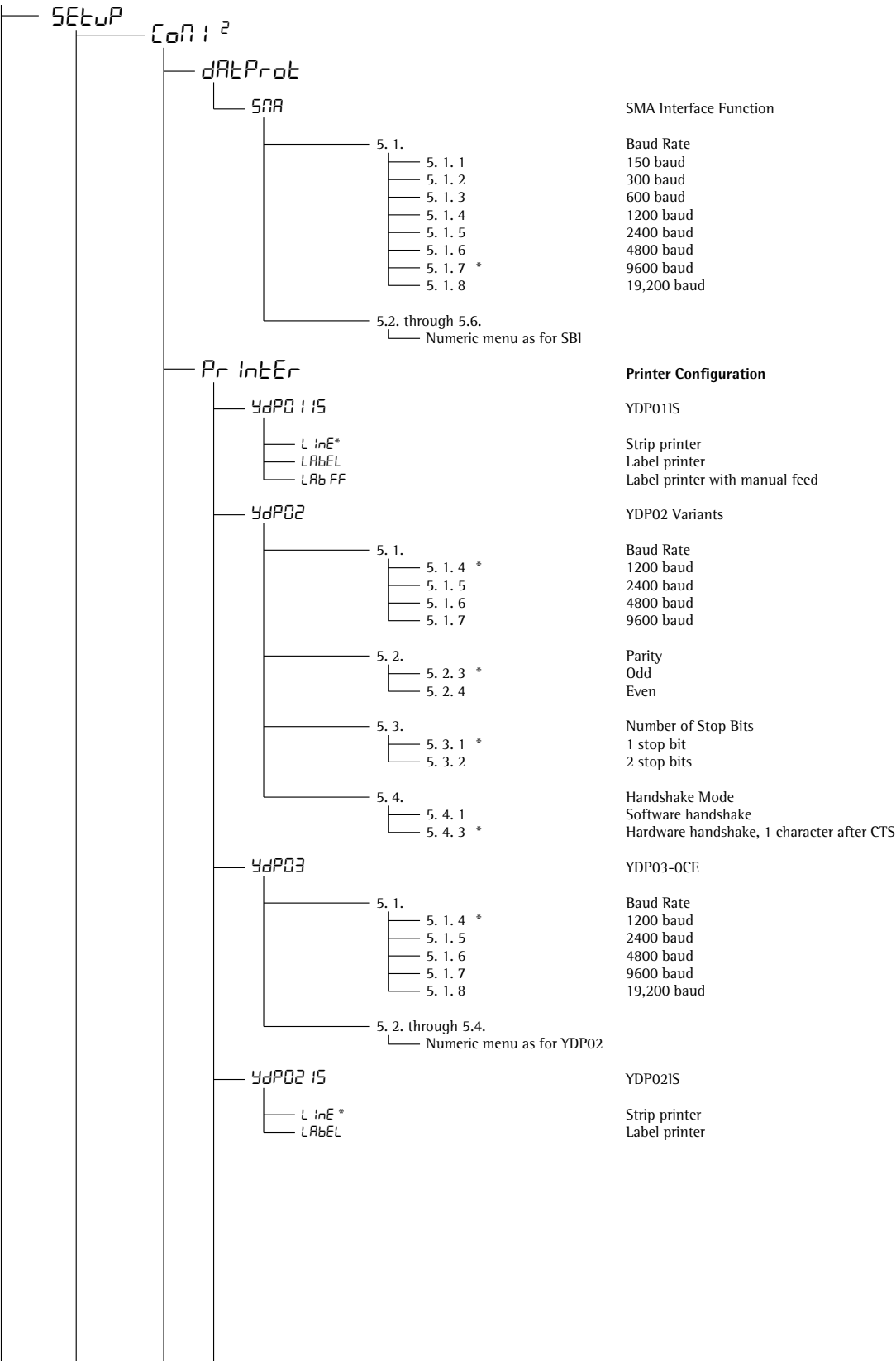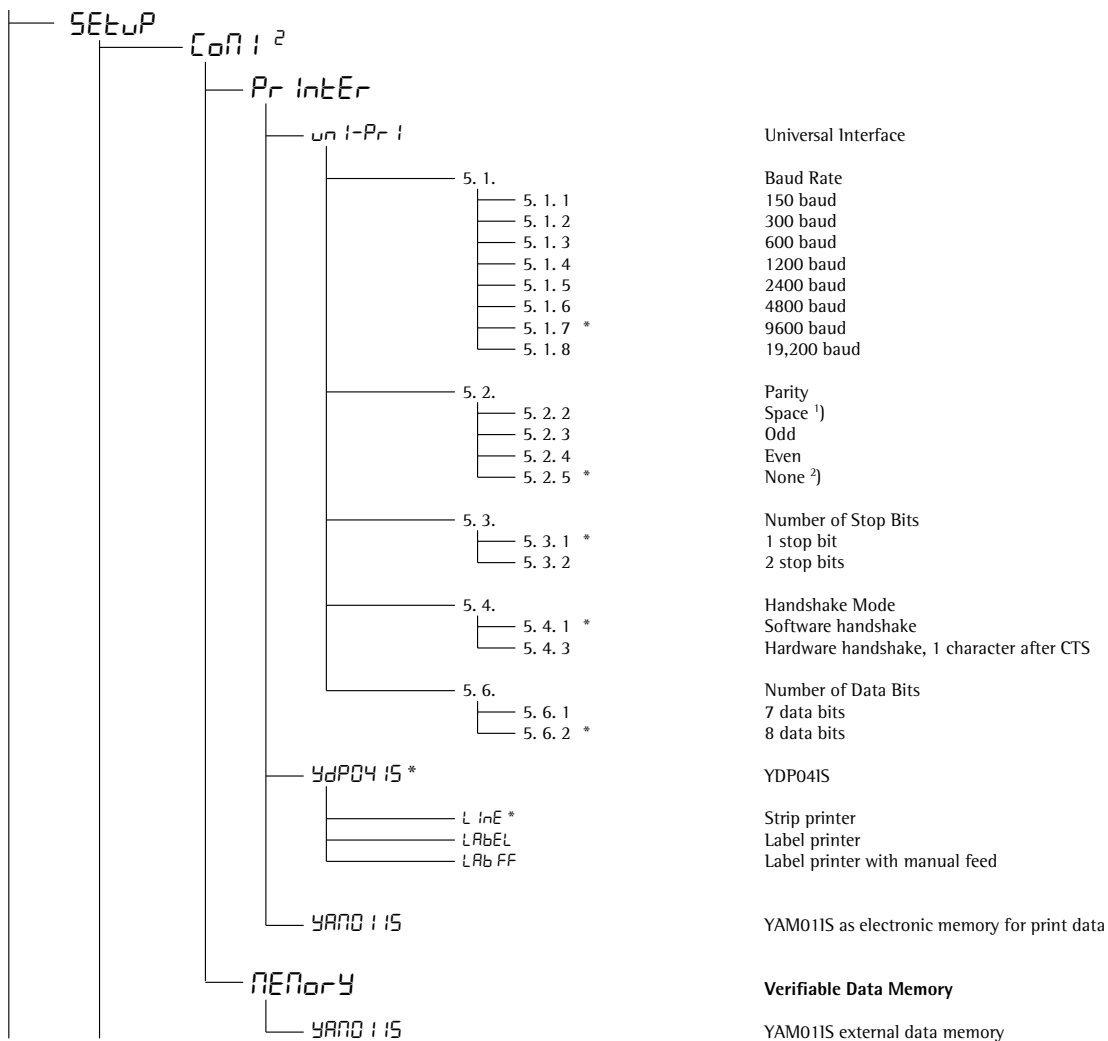
Configuration






¹⁾ = Not with setting 5.6.2 (8 bits)
²⁾ = Not with setting 5.6.1 (7 bits)

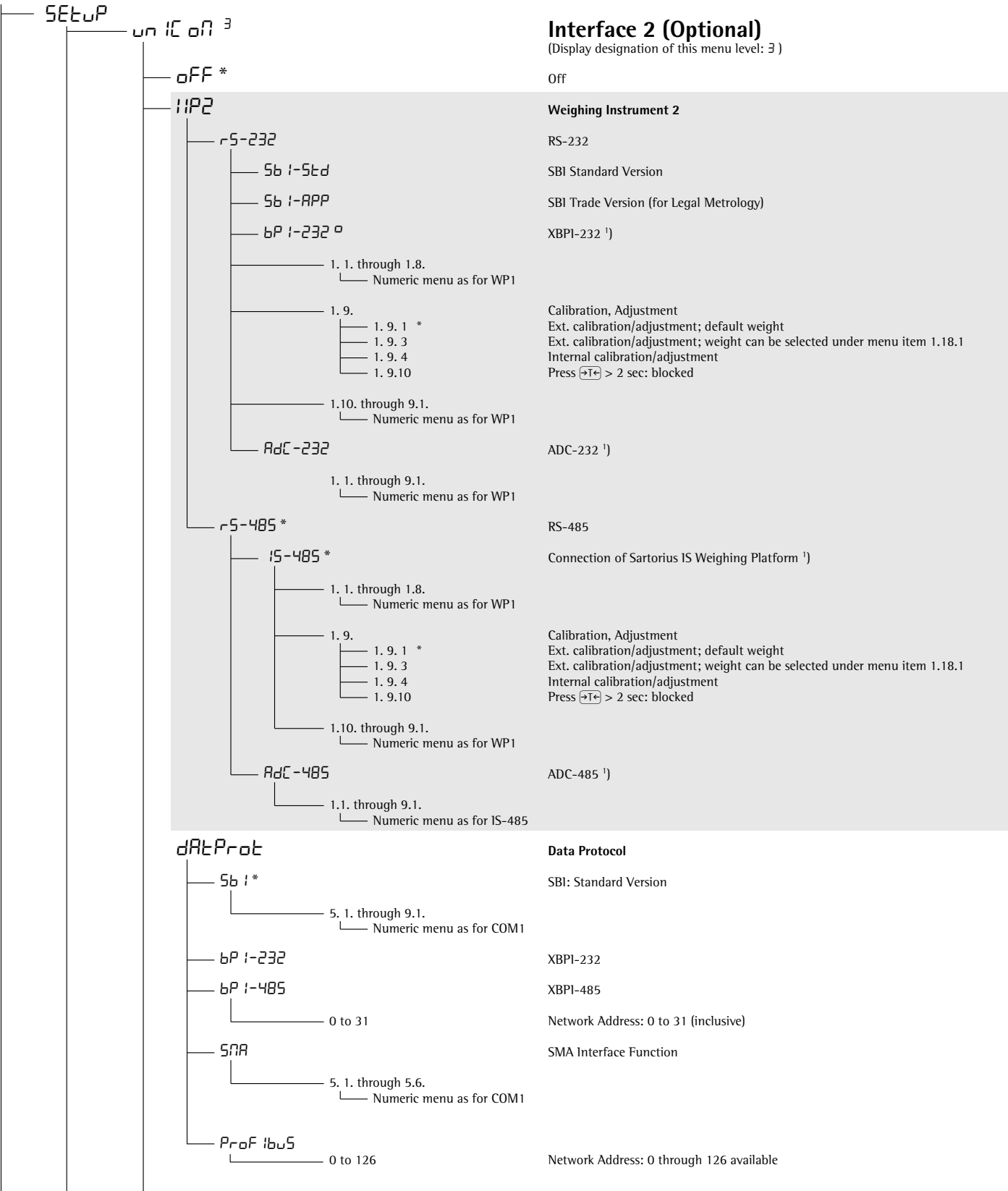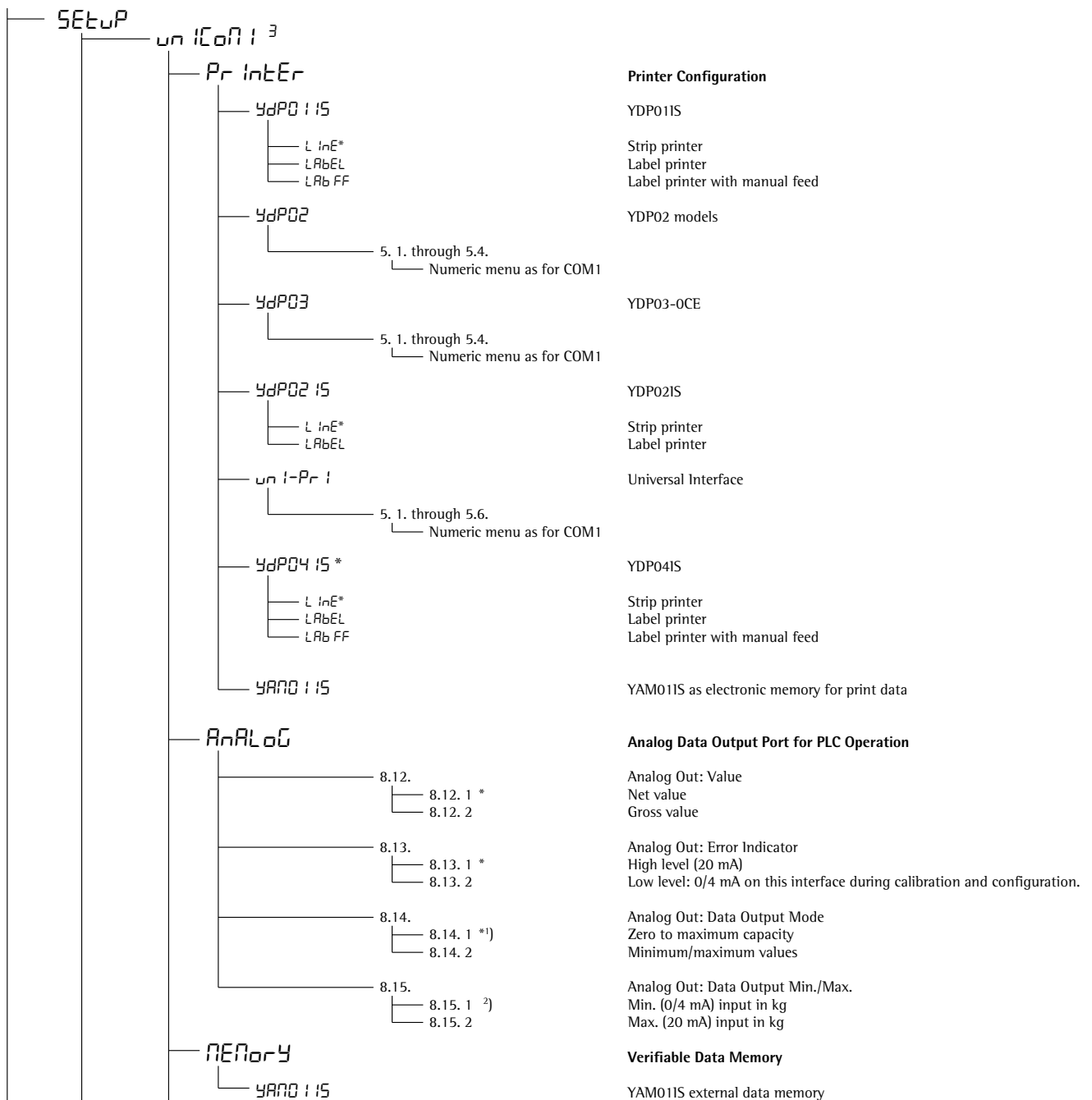
Configuration





¹⁾ = When setting 8.14.1 is active, the analog data output only works for XBPI weighing instruments
²⁾ = Not with setting 8.14.1

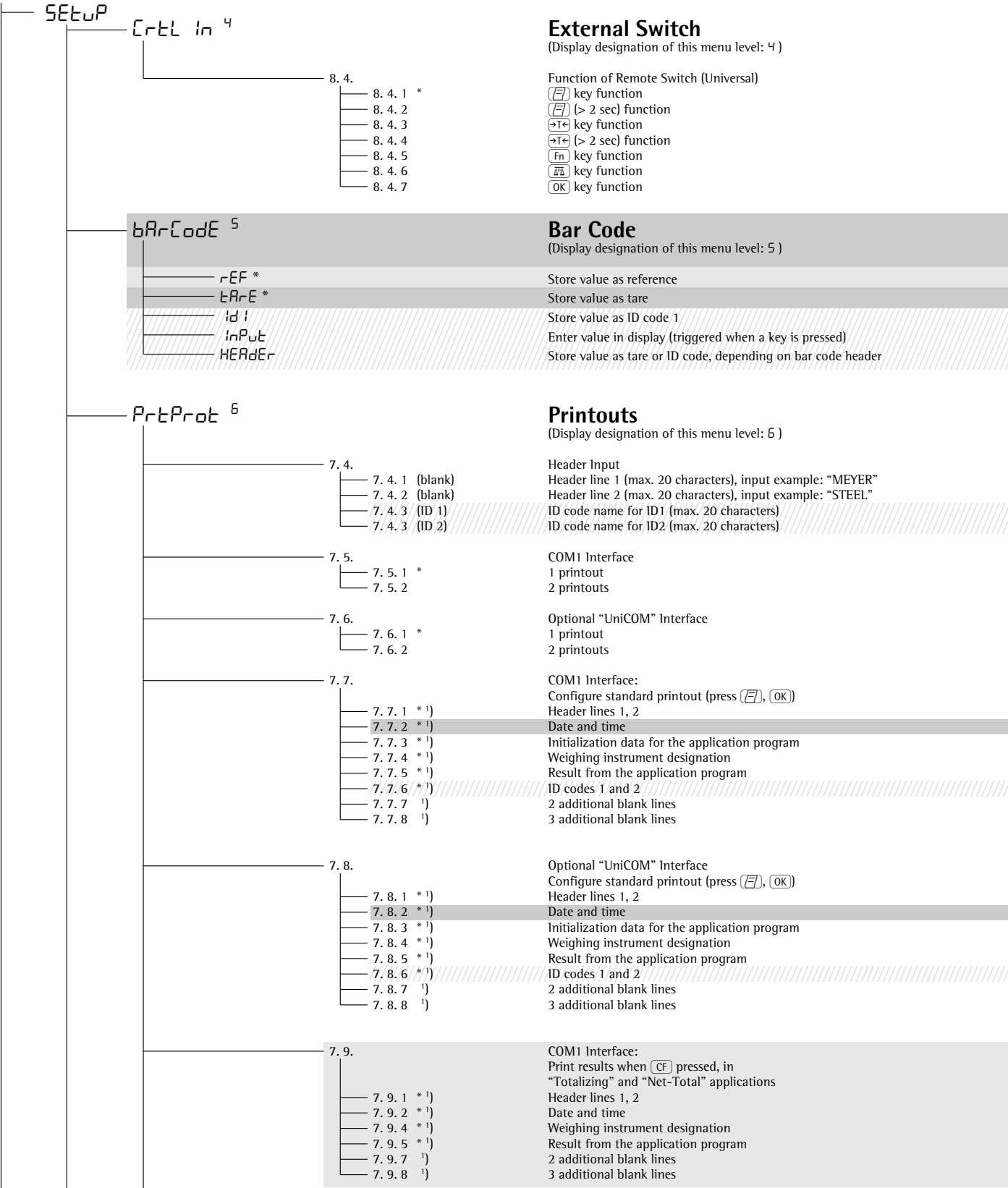Configuration



*) = More than one can be selected

SEtUP	
PrtProt ⁶	
7.10.	Optional "UniCOM" Interface Print results when [CF] pressed, in "Totalizing" and "Net-Total" applications
7.10. 1 * ¹⁾	Header lines 1, 2
7.10. 2 * ¹⁾	Date and time
7.10. 4 * ¹⁾	Weighing instrument designation
7.10. 5 * ¹⁾	Result from the application program
7.10. 7 ¹⁾	2 additional blank lines
7.10. 8 ¹⁾	3 additional blank lines
7.11.	ISO/GMP-compliant Printout
7.11. 1 *	Off
7.11. 2	Always for one measurement result
7.11. 3	Always for multiple measurement results
7.12.	Date/time printout line: Time not printed
7.12. 1 *	Off
7.12. 2	On
7.13.	One-time automatic printout at stability
7.13. 1 *	Off
7.13. 2	On
7.14.	FlexPrint
7.14. 1 *	Off
7.14. 2	On
9. 1.	Restore Factory Settings of the Numeric Menu for Data Protocol
9. 1. 1	Yes
9. 1. 2 *	No
ut IL It ⁷	
8. 2.	Acoustic Signal
8. 2. 1 *	On
8. 2. 2	Off
8. 3.	Keys
8. 3. 1 *	Accessible
8. 3. 2	Blocked
8. 7.	Automatic Shutoff, Indicator
8. 7. 1	Automatic shutoff acc. to menu item 8. 9.
8. 7. 2 *	No automatic shutoff
8. 8.	Display Lighting
8. 8. 1 *	On
8. 8. 2	Off
8. 8. 3	Automatic shutoff acc. to menu item 8.9.
8. 9.	Timer
8. 9. 1 *	After 1 + 1 minute not in use (after 1 min.: warning ²⁾ displayed for 2 minute)
8. 9. 2	After 2 + 2 minutes not in use (after 2 min.: warning displayed ²⁾ for 2 minutes)
8. 9. 3	After 5 + 5 minutes not in use (after 5 min.: warning ²⁾ displayed for 5 minutes)
8.11.	First Instrument Displayed
8.11. 1 *	Weighing instrument WP1
8.11. 2	Weighing instrument WP2
8.12.	Show the Geographical Data before Calibration
8.12. 1 *	No
8.12. 2	Yes
9. 1.	Restore Factory Settings of the Numeric Operating Menu
9. 1. 1	Yes
9. 1. 2 *	No

Operation

(Display designation of this menu level: ⁷)

¹⁾ More than one can be selected

²⁾ Warning: "ΔΔ" and 12 flash simultaneously

SETUP	TIME	Time Input example: 10.07.41 (hours.minutes.seconds)
	DATE	Date Input example: 01.05.02 (day.month.year) U.S. Mode: (month.day.year)
	CODE	Password Set, change and delete (max. 8 characters); example: 12345678
INFO	SERVICE	Device Information Service Information Service Date
	TERM	Indicator Model Serial number Version number of the indicator Software version
	11P-1	Weighing Instrument 1 Software version: 1st weighing instrument Latitude (in degrees) ¹⁾ Altitude (in meters) ¹⁾ Acceleration of gravity m/s ² ¹⁾
	11P-2	Weighing Instrument 2 (e.g., IS Weighing Platform) Model: 2nd weighing instrument Software version: 2nd weighing instrument Serial number Latitude (in degrees) ¹⁾ Altitude (in meters) ¹⁾ Acceleration of gravity m/s ² ¹⁾
	FLEX-Inf	FlexPrint File name ²⁾ ID ²⁾ Version ²⁾
LANG.	DEUTSCH	German
	ENGLISH *	English
	U.S. Mode	English with U.S. date/time format
	FRANC.	French
	ITAL.	Italian
	ESPAÑOL	Spanish

¹⁾ Output: either latitude and altitude or acceleration of gravity (depends on the input before verification)

²⁾ These three parameters are shown for each file loaded

Data Interfaces

The indicator is equipped with the following data interfaces:

- COM1: Standard interface (RS-232)
- UniCOM: Universal data interface (optional)

Both of these interfaces can be configured in the Setup menu for various input/output functions (e.g., Printer, second weighing instrument, PC, checkweighing/classification display). The optional UniCOM interface can be used for RS-232, RS-485 or RS-422 communication, or as voltage-/current (analog) interface. A bar code scanner (Combits 2 only) or an external rechargeable battery pack can be connected to the female UniCOM port (on CW1S, CW1NS and CW2S models, use the corresponding terminal screws).

Features

- Indicator models CW1P, CW1NP and CW2P (IP44 protection):
Connect via a 25-contact D-Sub female connector.
Use a T-connector (see “Accessories”) to connect a second device to the same interface.
- Indicator models CW1S, CW1NS and CW2S (IP67 protection):
Route connecting cable from the peripheral device to the indicator via a cable gland. Then connect the free ends of the cable using the terminal screws.
If you wish to connect a second peripheral device to the same interface port, use a separate cable gland to route the connecting cable of this device into the indicator.

⚠ Warning When Using RS-232 Connecting Cables Not Supplied by Sartorius:

The pin assignments in the cable might not be compatible with Sartorius equipment. Check all pin assignments against the cabling diagrams and disconnect any lines that are not assigned. Failure to do so may damage or even completely ruin your indicator and/or peripheral device.

Specifications

Serial interface:

Operating mode:	Full duplex	
Standard:	COM1: UniCOM ¹⁾ :	RS-232, RS-232 or RS-422/RS-485
Interface connector:	CW1P, CW1NP and CW2P (IP44 protection): 25-contact D-Sub female connector CW1S, CW1NS and CW2S (IP67 protection): The cable is connected to terminal screws inside the housing and routed into the housing via a cable gland.	
Transmission rates:	150, 300, 600, 1200, 2400, 4800, 9600 and 19,200 baud (depending on the operating mode)	
Number of data bits:	7 or 8 bits	
Parity:	Space, odd, even, none (depending on the operating mode)	
Number of stop bits:	1 or 2 stop bits	
Handshake mode:	Software (XON/XOFF) or hardware (1 character after CTS)	
Communication mode:	SBI, XBPI-232 ²⁾ , XBPI-485 ¹⁾²⁾ , MP8 binary ³⁾ , SMA Available printers: – YDP01IS – YDP02IS-Label – YDP01IS-Label – Universal – YDP02 – YDP04IS – YDP03 – YDP04IS-Label – YDP02IS – YAM01IS Alibi memory	
Network address ⁴⁾ :	0, 1, 2, (...), 31	
SBI: Manual data output:	Without stability, after stability, configurable printout	
SBI: Automatic data output:	Without stability, at stability, at user-defined intervals	
SBI: Output format:	16 or 22 characters	
Printout of application data:	Output of a configurable printout	

Analog UniCOM interface (optional)

Standard:	4 to 20 mA, 0 to 20 mA, 0 to 5V
Power supply:	Internal or external
Factory setting:	4 to 20 mA, internal power supply
Interface connector:	CW1P, CW1NP and CW2P indicators (IP44 protection): 25-contact D-Sub female connector CW1S, CW1NS and CW2S indicators (IP67 protection): The free ends of the cable are connected to terminal screws inside the housing; the cable is routed into the housing via a cable gland.

¹⁾ Optional UniCOM universal data interface

²⁾ XBPI operating mode: 9600 baud, 8 data bits, parity: odd, 1 stop bit

³⁾ Only with the standard COM1 interface

⁴⁾ Network address is valid only in the XBPI mode

Data Interfaces

Options for Connecting Peripherals

You can connect the following printers to the COM1 and UniCOM interfaces:

- YDP02 (user-definable interface parameters)
- YDP03 (user-definable interface parameters)
- YDP011S (strip or label printer)
- YDP021S (strip or label printer)
- YDP041S (strip or label printer)
- Universal printer (user-definable transmission parameters)
- YAM011S Alibi memory

The following devices can also be connected to the COM1 interface:

- Foot switch / hand switch
- PC (RS-232 interface)
- Second weighing instrument (Combics 2 only, RS-232 interface)
- External checkweighing display (red/yellow/green) over the digital I/O (Sartorius standard)

The connector of the UniCOM universal port enables the following devices to be connected:

- External rechargeable battery pack
- Bar code scanner (Combics 1 plus and Combics 2 only)

The following devices can also be connected to the UniCOM universal interface:

- PC (RS-232 interface)
- Second weighing instrument (Combics 2 only; can be switched from RS-232 to RS-485 operating mode)
- Second printer (external power source required)
- Remote display
- Current interface (0/4 - 20 mA)

⚠ If necessary, use an external power source to power peripheral devices.

Connecting a second weighing instrument:

Combics 2 enables connection of a second weighing instrument to either the COM1 or the UniCOM port.

COM1 operates in RS232 mode.

A second weighing instrument on this port can use the following operating modes:

- SBI
- XBPI-232 (factory setting)
- ADC-232

UniCOM can operate in either the RS-232 mode or in RS-485 mode.

A second weighing instrument on this port can use the following operating modes:

- SBI (RS-232 mode)
- XBPI-232 (RS-232 mode)
- ADC-232 (RS-232 mode)
- IS-485 (RS-485 mode, XBPI mode; factory setting)
- ADC-485 (RS-485 mode)

Connecting a Printer

The standard COM1 port or the optional universal UniCOM interface (or both) can be used as a printer interface.

Operation as a COM Port

For operation as a COM port, you can adapt the data protocol to the following operating modes:

- SBI (factory setting)
- XBPI-232
- XBPI-485 (only UniCOM)
- MP8-binary (only COM1)
- SMA

You can operate the COM1 port and the UniCOM port independently of one another (i.e., for transferring data and controlling equipment via a PC while simultaneously outputting data to the COM1 printer port).

In the SBI communication mode, you can control a display unit and a connected weighing instrument by sending ESC commands from a PC to the communications port (COM1 or UniCOM) (see page 96).

⚠ When using two interfaces (COM1 and UniCOM) each device can be configured only once. Devices are as follows:

- WP1-2
- Data communication
- External Alibi memory
- Analog output port

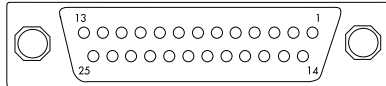
If you attempt to configure a device a second time (e.g., for the UniCOM port) after it has been configured another port (e.g., "Data communication" on COM1), error code INF 74 is displayed.

Pin Assignment Charts

**Models CW1P, CW1NP, CW2P
(IP44-protected):**

Female Connectors COM1 and UniCOM:

25-contact D-Submini DB25S with
screw lock hardware



Front view

Male interface connector used
(please use connectors with the same
specifications):

25-pin D-Submini DB25, with inte-
grated shielded cable clamp assembly
(Amp type 826 985-1C) and fastening
screws (Amp type 164868-1)

Pin assignments, COM1:

Pin 1: Shield
Pin 2: Data output (TxD)
Pin 3: Data input (RxD)
Pin 4: Not connected
Pin 5: Clear to send (CTS)
Pin 6: Internally connected
Pin 7: Internal ground (GND)
Pin 8: Internal ground (GND)
Pin 9: Not connected
Pin 10: Not connected
Pin 11: +12 V for printer
Pin 12: RES_OUT\
Pin 13: +5V
Pin 14: Internal ground (GND)
Pin 15: Universal remote switch
Pin 16: Control output "lighter"
Pin 17: Control output "equal"
Pin 18: Control output "heavier"
Pin 19: Control output "set"
Pin 20: Data terminal ready (DTR)
Pin 21: Supply ground (GND)
Pin 22: Not connected
Pin 23: Not connected
Pin 24: Power supply +15...25 V
Pin 25: +5V

Pin Assignment Chart: 2nd Interface:
Connection of external rechargeable
battery and bar code scanner¹⁾ (optional
UniCOM interface not installed)

Pin 1: Shield
Pin 2: Not connected / *
Pin 3: Not connected / *
Pin 4: Internal ground (GND)
Pin 5: Not connected / *
Pin 6: Not connected / *
Pin 7: Not connected / *
Pin 8: Not connected / *
Pin 9: Not connected / *
Pin 10: Not connected / *
Pin 11: +12 V for printer
Pin 12: RES_OUT\
Pin 13: +5 V switch
Pin 14: Internal ground (GND)
Pin 15: Keyboard data
Pin 16: Not connected / *
Pin 17: Not connected / *
Pin 18: Not connected / *
Pin 19: Keyboard clock
Pin 20: Not connected / *
Pin 21: LINE_1 _GND
Pin 22: LOW_BATT²⁾
Pin 23: BATT_ON_OFF³⁾
Pin 24: LINE_1_B
Pin 25: +5V

* Pin assignments depending on the UniCOM used

1) Combics 2 only

2) Signal from battery pack: battery drained

3) Switch off battery pack when weighing instrument
switched off

Data Interfaces

Models CIS1, CISN1 and CIS2 (IP67-protected):

Connection of open cable ends to terminal screws inside the indicator

COM1 interface connections:


□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□
11	12	13	14	15	16	17	18	19	20											
⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘
1	2	3	4	5	6	7	8	9	10											
□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□

Top view

Terminal assignments

- No. 1: Universal switch
- No. 2: Control output “set”
- No. 3: Control output “heavier”
- No. 4: Control output “equal”
- No. 5: Control output “lighter”
- No. 6: Clear to send (CTS)
- No. 7: Data output (TxD)
- No. 8: Data input (RxD)
- No. 9: Data terminal ready (DTR)
- No. 10: Internal ground (GND)
- No. 11: LINE_A
- No. 12: LINE_A
- No. 13: GND_LINE_A
- No. 14: GND_LINE_A
- No. 15: +12 V for printer
- No. 16: Reset output
- No. 17: +5 V
- No. 18: +5 V
- No. 19: Ground (GND)
- No. 20: Ground (GND)

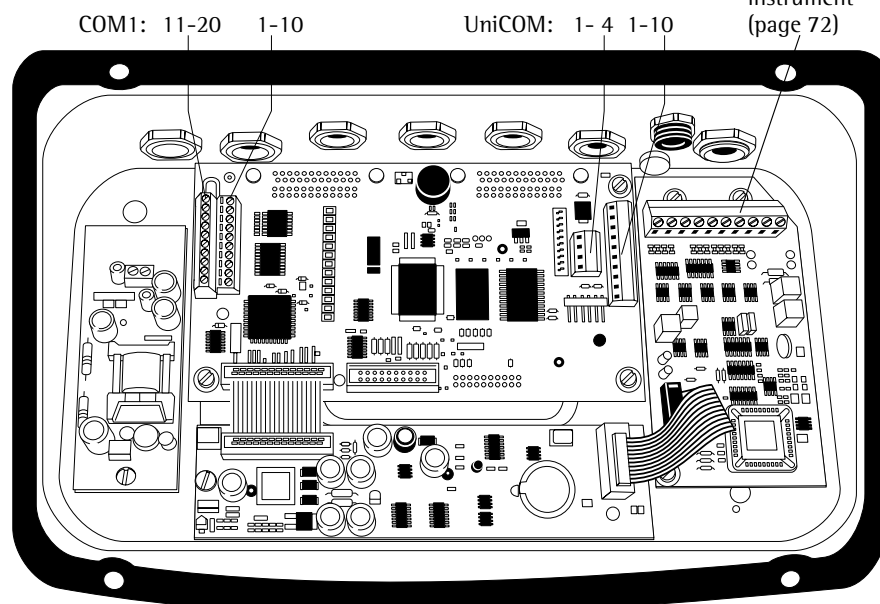Models CIS1, CIS1N and CIS2:
Terminals on the PCB



Second connection: For connecting an external battery pack and a bar code scanner ¹⁾
(optional UNICOM interface not installed)

⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘	⌘
1	2	3	4	5	6	7	8	9	10											
□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□	□
⌘	⌘	⌘	⌘																	
1	2	3	4																	
□	□	□	□																	

Terminal assignments in the
10-terminal strip

- No. 1: Not connected*
- No. 2: GND
- No. 3: GND
- No. 4: +5V Switch
- No. 5: Not connected *
- No. 6: Keyboard clock
- No. 7: Keyboard data
- No. 8: Not connected *
- No. 9: Not connected *
- No. 10: Not connected *

Terminal assignments in the
4-terminal strip

- No. 1: Supply, ground (GND_LINE_B)
- No. 2: Supply, battery pack
- No. 3: LOW_BATT ²⁾
- No. 4: BATT_ON_OFF ³⁾

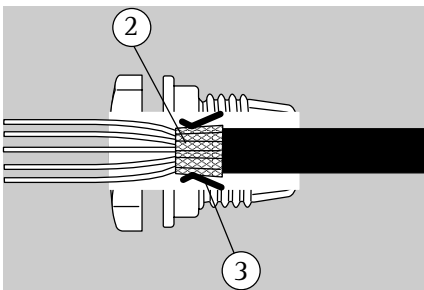
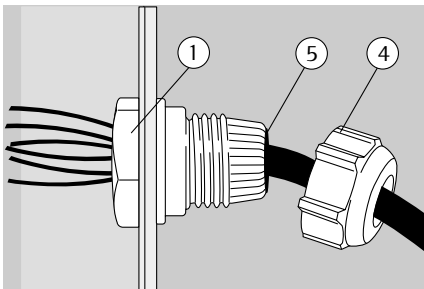
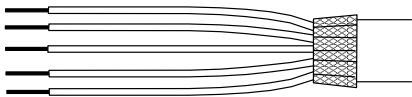
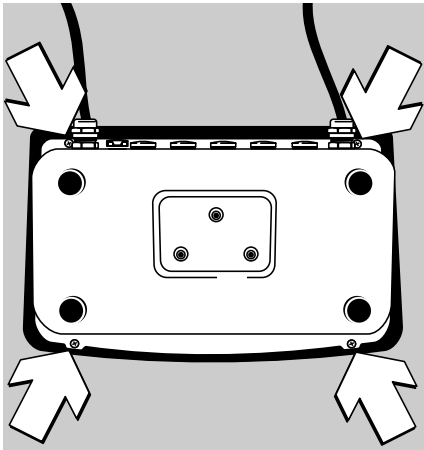
* Pin assignments depending on the UniCOM used

¹⁾ Combics 2 only

²⁾ Signal from battery pack: battery drained

³⁾ Switch off battery pack when weighing instrument switched off

Installing the Interface Cable



⚠ Peripheral devices should be connected to the interfaces only by a trained and authorized Sartorius technician. Any installation work that does not conform to the instructions in this manual will result in forfeiture of all claims under the manufacturer's warranty.

⚠ Disconnect the equipment from power (unplug from the wall socket) before beginning any installation work.

- Open the Combics indicator:
remove the four cap nuts from the front panel and remove the panel.

- Prepare the cable:

- Expose approx. 10 cm (4 in) of the cable end for installation
- Remove all but approx. 1 cm (½ in) of the shielding and fold it back over the casing
- Strip the casing from approximately 1 cm (½ in) of the wires and attach ferrules to the wire ends.

- Attach the cable gland:

⚠ Please use extreme caution when performing any work on the equipment that affects this cable gland. Use a torque wrench and tighten the cable gland to 5 Nm.

- Remove the protective cap from the bore hole on the indicator.
- Guide the enclosed cable gland through the bore hole and secure it inside the housing with the nut (1).

- Slide the cable gland over the cable until the clamps (3) are in contact with the shield (2). Tighten the nut (4) until the sealing clamp (5) forms a slight ridge between nut and cable.
- Make sure the shield is in contact with the clamps.


- Connect the wires securely in accordance with the terminal assignments.

- After you close the housing again, use a pressure gauge to check the integrity of the IP67-protection. For details, contact the Sartorius Service Center.

Data Interfaces



Cabling Diagram (Adapter Cable for PC)

(Model CW1P | CW1NP | CW2P: adapter cable 7357312; model CW1S | CW1NS | CW2S: connecting cable YCC02-D9F6).
Diagram for connecting a computer or other peripheral device to the indicator using the RS-232-C/V24 standard and cables up to 15 m (50 ft.) long:

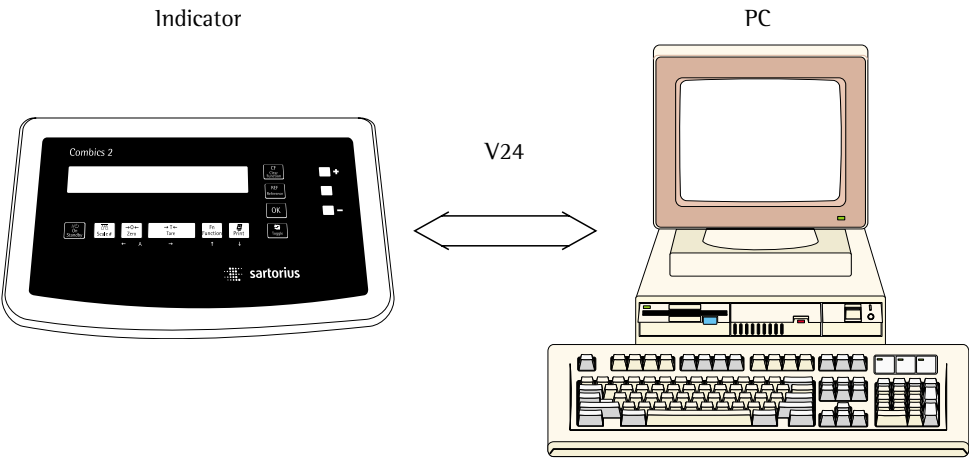

Cabling diagram

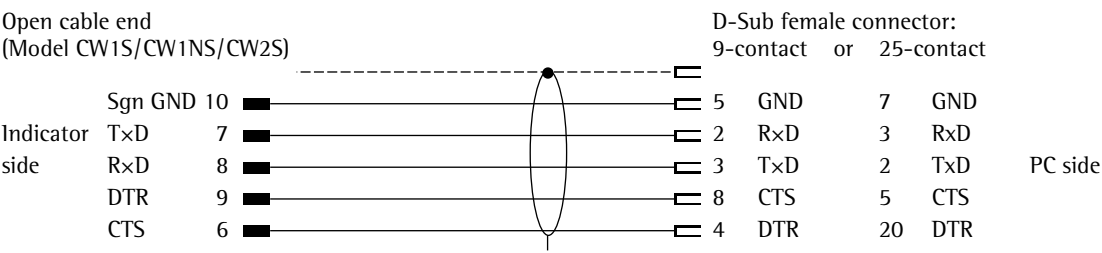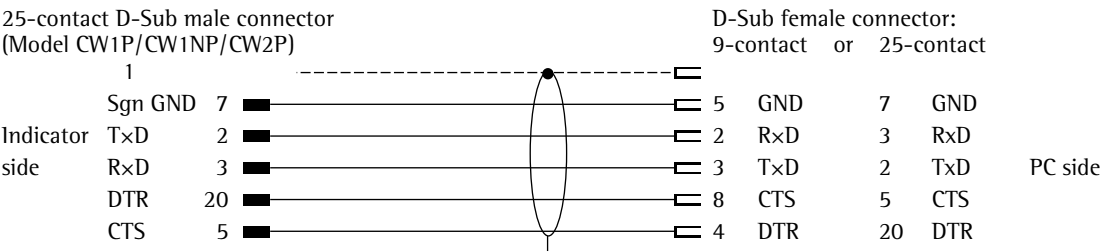Connection assignments for the cable from the indicator to an RS-232 PC interface



Synchronization

Data communication between the indicator and a computer takes the form of messages (“telegrams”) made up of ASCII code. For error-free data communication, the settings for baud rate, parity, handshake mode and character format must be the same at both ends.

You can configure the interface settings in the Setup menu so that they match those of the computer. You can also define parameters in the indicator to make data output dependent on various conditions. The conditions that can be configured are listed in the descriptions of the application programs.

If you do not connect a peripheral device to the indicator’s interface port, this will not generate an error message.

Handshake

The weighing instrument interface (Sartorius Balance Interface = SBI) has transmit and receive buffers. You can define the handshake parameter in the indicator’s Setup menu:

- Hardware handshake (CTS/DTR)
- Software handshake (XON, XOFF)

Hardware Handshake

Hardware handshake with a 4-wire interface: 1 more character can be transmitted after CTS (clear-to-send).

Software Handshake

The software handshake is controlled via XON and XOFF. When a device switched on, XON must be transmitted to enable a connected device to communicate.

When the software handshake is configured in the Setup menu, the hardware handshake becomes active after the software handshake.

The data transmission sequence is as follows:

Scale	---	byte	---	Computer
(trans-	---	byte	---	(receiving
mitting	---	byte	---	device)
device)	---	byte	---	
	<--	XOFF	---	
	---	byte	---	
	---	byte	---	
	...			
	(Pause)			
	...			
	<--	XON	---	
	---	byte	---	
	---	byte	---	
	---	byte	---	
	---	byte	---	

Transmitting Device

Once XOFF has been received, it prevents further transmission of characters. When XON is received, it re-enables the transmitting device to send data.

Receiving Device

To prevent too many control commands from being received at one time, XON is not transmitted until the buffer is almost empty.

Data Interfaces

Configuring the Data Interface as a COM Port (dAtEPrOt)

Configure the interface as a COM port in the Setup menu as a COM1 or UniCOM port under the “Data Protocol” (dAtEPrOt) menu item.

SBI Communication

This is a simple ASCII interface. Data output is configured under menu items 6.1 and 6.3:

- Manual output of displayed value with or without stability (menu items 6.1.1 and 6.1.2)
- Automatic output of displayed value with or without stability (menu items 6.1.4 and 6.1.5) at intervals defined in display updates. The number of display updates comprising an output interval is configured under menu item 6.3.
- Output of a printout as configured in the Setup program (menu item 6.1.7). Output is linked to the “Data Protocol” menu item (dAtEPrOt) (see page 100, “Configuring Printouts”).

If you do not activate and configure a user-definable data record, the printout simply contains the current value displayed on the indicator (weight with unit, calculated value, alphanumeric display).

SMA Communication

Standardized communications protocol of the Scale Manufacturers Association

MP8 Binary Purpose

With the MP8 interface you can connect MP8-generation peripheral devices with separate power supplies to the Combics indicator.

Features

- The weighing instrument is used only for determining weight values.
- The data interface transmits only in MP8 binary protocol
- The application program for MP8 can be selected under menu item 3
- The program index 2 for MP8 can be selected under menu item 4.
- “MP8 interface emulation” is not permitted in legal metrology.

Data Input Format

You can connect a computer to your indicator to send commands controlling weighing instrument functions and applications via the interface port. All commands use the same format (data input format) starting with the ESC character (ASCII 27) and ending with a carriage return (CR; ASCII 13) and a line feed (LF; ASCII 10). The total length of a command is anywhere from 4 characters (1 command character between the start and end described above) to 7 characters (4 command characters).

The table below shows the available command characters; each command must be flanked by the start and end characters as described above.

Example: The command character for output is “P” (“output to Port”). The string “ESC P CR LF” triggers this command.

Command	Meaning
K	Weighing mode 1
L	Weighing mode 2
M	Weighing mode 3
N	Weighing mode 4
O	Block all keys
P	Output readout to data interface
Q	Output acoustic signal
R	Release (unblock) keys
T	Tare and zero (combination tare function)
f3_	Zero (see also the “kZE_” command)
f4_	Tare without zeroing (see also the “kT_” command)
i_	Information about the indicator Example of output: “C12/012502/1” Meaning: Indicator: Combics 2, software version: 012502, active weighing instrument: 1
kF1_	F1: Trigger (Fn) key function
kF2_	F2: (CF) key function (Combics 2 only)
kF3_	F3: (REF) key function (Combics 2 only)
kF4_	F4: (OK) key function (Combics 2 only)
kF5_	F5: (S) key function (Combics 2 only)

Command	Meaning
kF6_	F6: Trigger (Info) key function (Combics 1 plus only)
kF7_	F7: (ID1) key function (Combics 1 plus only)
kF8_	F8: (ID2) key function (Combics 1 plus only)
kCF_	CF: (CF) key function (Combics 2 only)
kP_	(E) key function Output to printer port
kT_	(T) key function (tare)
kNW_	(T/A) key function (Toggle the weighing instrument)
kZE_	(+0<) key function (zero the instrument)
x1_	Output model designation of active weighing instrument. Example: “LP6200S-0C”
x2_	Output serial number of active weighing instrument; example: “0012345678 ”
x3_	Output software version of active weighing instrument; example: “ 00-20-04 ”
x4_	Output software version of indicator; example: “ 01-25-02 ”
x9_	Output serial number of indicator; example: “0012345678 ”
x10_	Output model of indicator; example: “CW2P4-1500RR-LCE”
z1_	Activate input for printout header 1
z2_	Activate input for printout header 2
txx...x_	xx...x: Input text Length acc. to input (Combics 1 plus only)

The ASCII code for the “underline” character is 95.

Format for entering printout header lines: ESC z x a ... a _ CR LF where x=(header line) 1 or 2; a ... a: a...a= up to 20 characters of text, followed by the “underline” character, carriage return and line feed.

Data Output Format

You can output the values displayed in the line for measured values and the weight unit with or without a data ID code. The first 6 characters, called the data header, show the data ID code identifying the subsequent value. Select menu item 7.2 to define whether the data ID code is included in output or not.

Examples:

+	235	pcs	Without data ID code
Qnt +	235	pcs	With data ID code

Display segments that are not activated are output as spaces. Values with no decimal point are output without a decimal point.

Data Output Format with 16 Characters (without Data Header)

Normal Operation

Position 1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	+	*	D	D	D	D	D	D	D	*	U	U	U	CRLF	
or	-	*	D	D	D	D	D	D	D	*	U	U	U	CRLF	
or	*	*	*	*	*	*	*	*	*	*	*	*	*	*	CRLF

+ -: Plus or minus sign
 *: Space
 D: Digit or letter (max. 7 characters plus decimal point)
 U: Unit symbol (1 - 3 letters, followed by 0 - 2 spaces)
 CR: Carriage return
 LF: Line feed

Special Codes

Position 1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	*	*	*	*	*	*	-	-	*	*	*	*	*	*	CRLF
or	*	*	*	*	*	*	H	*	*	*	*	*	*	*	CRLF
or	*	*	*	*	*	*	H	H	*	*	*	*	*	*	CRLF
or	*	*	*	*	*	*	L	*	*	*	*	*	*	*	CRLF
or	*	*	*	*	*	*	L	L	*	*	*	*	*	*	CRLF
or	*	*	*	*	*	*	C	*	*	*	*	*	*	*	CRLF

*: Space
 - -: Final readout mode
 H: Overload
 HH: Overload in Checkweighing
 L: Underload
 LL: Underload in Checkweighing
 C: Calibration/adjustment

Error Codes

Position 1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	*	*	*	E	r	r	*	*	#	#	*	*	*	*	CRLF
or	*	*	*	E	r	r	*	*	#	#	*	*	*	*	CRLF

*: Space
 #: Error code number (2 or 3 digits)

Example (output of value: +1255.7 g):

Position 1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	+	*	*	*	1	2	5	5	.	7	*	g	*	*	CRLF

Position 1: Plus or minus sign or space
 Position 2: Space
 Positions 3-10: Weight value with decimal point; leading zeros are output as spaces.
 Position 11: Space
 Positions 12-14: Unit symbol or space
 Position 15: Carriage return
 Position 16: Line feed

Data Output Format with 22 Characters

Normal Operation

Position 1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22
	1	1	1	1	1	1	+	*	D	D	D	D	D	D	D	*	U	U	U	CRLF	
or	1	1	1	1	1	1	-	*	D	D	D	D	D	D	D	*	U	U	U	CRLF	
or	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	*	CRLF

11: ID code character, right-justified with spaces
 + -: Plus or minus sign
 *: Space
 D: Digit or letter (max. 7 characters plus decimal point)
 U: Unit symbol (1 - 3 letters, followed by 0 - 2 spaces)
 CR: Carriage return
 LF: Line feed

Special Codes

Position 1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	
	S	t	a	t	*	*	*	*	*	*	*	*	*	*	-	-	*	*	*	*	*	CRLF
or	S	t	a	t	*	*	*	*	*	*	*	*	*	*	H	*	*	*	*	*	*	CRLF
or	S	t	a	t	*	*	*	*	*	*	*	*	*	*	H	H	*	*	*	*	*	CRLF
or	S	t	a	t	*	*	*	*	*	*	*	*	*	*	L	*	*	*	*	*	*	CRLF
or	S	t	a	t	*	*	*	*	*	*	*	*	*	*	L	L	*	*	*	*	*	CRLF
or	S	t	a	t	*	*	*	*	*	*	*	*	*	*	C	*	*	*	*	*	*	CRLF

** : Space
 H: Overload
 L: Underload
 C: Calibration/adjustment
 - -: Final readout mode
 HH: Overload in Checkweighing
 LL: Underload in Checkweighing

Data Interfaces

Error Codes

Position	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22
		S	t	a	t	*	*	*	*	*	E	r	r	*	*	#	#	*	*	*	*	CRLF
or		S	t	a	t	*	*	*	*	*	E	r	r	*	*	#	#	*	*	*	*	CRLF

*: Space

#: Error code number (2 or 3 digits)

ID Code Characters

ID character	Meaning
G #	Gross value
N	Net value
T	Application tare memory 1
T 2	Application tare memory 2
D i f f	Difference from calibration value
N o m .	Exact calibration weight
N o m .	Exact calibration weight using SBI output
n R e f	Reference sample quantity
p R e f	Reference percentage
w R e f	Reference sample weight
Q n t	Result from Counting application Result from Counting (piece count) and Neutral Measurement applications
m D e f	Target value for Animal weighing
x - N e t	Result from Animal Weighing
S e t p	Target value for Checkweighing
W . D i f f	Absolute difference (e.g., in kg) in Checkweighing
L i m	Deviation in % in Checkweighing
M a x	Upper limit for Checkweighing
M i n	Lower limit for Checkweighing
S t a t	Status
C l a s s x	Classification
L i m x	Class limit
D	Percentage (as loss)
P r c	Percentage (as residue)
W x x %	Reference percentage weight
C o m p x x x	Component xxx
C o n t . T	Contents of the tare memory in Net-total Formulation
T o t . c p	Total weight in Net-Total Formulation
P T 2	Preset tare
n	Transaction counter
* G	Sum of gross weights in Totalizing
* N	Sum of net weights in Totalizing
S e r . n o	Serial number of the weighing instrument or indicator

Example (output of value: +1255.7 g):

Position	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22
		G	#	*	*	*	*	+	*	*	*	1	2	5	5	.	7	*	g	*	*	CRLF

Positions 1-6: ID code character, right-justified with spaces

Position 7: Plus or minus sign or space

Position 8: Space

Positions 9-16: Weight value with decimal point;
leading zeros are output as spaces.

Position 17: Space

Positions 18-20: Unit symbol or space

Position 21: Carriage return

Position 22: Line feed

⚠ If the weight value is output with 10-fold increased resolution, this value is not permitted to be printed or saved in a weighing instrument operated in legal metrology in the SBI mode. In this case, the unit symbol is not included with output.

Configuring the Data Interface as a Printer Port (PrInTEr)

You can connect one or two strip printers or one or two label printers to the Combics. Configure the COM1 and UniCOM interfaces as printer ports under the “PrInTEr” menu item.

There are several actions that generate the command for outputting data to the printer port:

- Pressing the **[F]** key. If the operating menu is active, all menu settings under the active menu level are printed.
- On receipt of the SBI command “Esc k P _”. For details, see “Data Input Format” in this chapter.
- In some applications, pressing a given key (e.g., to save a value or start a routine) also generates a print command. In this case, a configurable printout is generated with application-specific data.

The **☉** **↕** symbols are displayed when data is being output to the printer port.

Configuring Printouts

Configure printouts in the Setup menu is under the “Printouts” menu item (PrTEPrTE). The printout should be formatted only after the desired application has been configured, as some of the positions are application-dependent. You can configure a different printout for each interface. Each printout contains your choice of the information blocks described in the following; to enable or disable a block in the printout, select it or deselect it in the Setup menu. Combics 2 only: For the “Totalizing” and “Net-Total” applications, you can also configure summarized printouts (results) independent of individual component value printouts.

Block 1: Headers

You can define 2 headers, each with 20 characters per line (e.g., for printing your company’s name). Enter the header lines under menu items 7.4.1 and 7.4.2. Blank header lines are not printed.
Example: format of Block 1:

ACE HARDWARE
GOETTINGEN

In this example, the company name is centered on the printout. This was achieved by entering blank spaces at the beginning of each line.

Block 2: Date/Time

(not on Combics 1)
Example: format of Block 2:

21.01.2001 16:02

To obtain a standardized time stamp (e.g., for documentation within a completely automated system), you can suppress the printout of the time stamp in the “Date / time” information block. This function is controlled under menu code 7 -12, “Time not printed”. The factory setting is “Off” (i.e., the time is included on the printout). If you select “On” for this menu item, the time stamp can be inserted by a higher-level controller or central computer to maintain consistent time stamping. This setting is especially important for communication with a PC.

Separating Block:

Dotted line, blank line (for the Weighing application).
This block is automatically inserted before further information blocks are printed.

Block 3: Initialization Data

Which data is included in this block depends on the active application. It can include, for example, reference sample quantity, reference piece weight, target weight, etc. The block is terminated with a blank line.
This block can only be activated for the standard printout. It cannot be selected for the printout of results.
Example: format of Block 3 (Counting application)

nRef 10 pcs
wRef + 0.035 kg

Block 4: Serial Number

Example: format of Block 4:

Ser.no. 1234567890

Block 5: Results

Which data is included in this block depends on the active application. Normally it includes gross, net and tare weights, followed by a blank line. The result is printed after this data; for example, the piece count. The block is terminated with a blank line.
Example: format of Block 3 (Counting application)


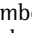
G# + 1.402 kg
T + 0.200 kg
N + 1.202 kg


Qnt 34 pcs

Data Interfaces

GMP-compliant Printouts

When the corresponding menu item is active, the measured result is bracketed on the printout by a GMP header and a GMP footer (GMP = “Good Manufacturing Practice”).

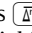The GMP header precedes the first measured result. The GMP footer is printed either after each measured result (“ISO/GLP/GMP: For 1 application result,” menu item 7.11.2), or after the last result in a series of measurements (“ISO/GLP/GLP: For several application results,” menu item 7.11.3). To end a series of measured results, press and hold the  key (> 2 sec). In this case, the  symbol is displayed after the GMP header is printed and remains in the display until the GMP footer is printed.

If you toggle to a different weighing instrument (Combics 2 only) while a GMP printout of several measured results is being generated, the GMP footer for the weighing instrument used up to that point is generated when you press . The GMP header for the other weighing instrument is included on the next printout generated.

A GMP-compliant printout is generated automatically at the conclusion of calibration/adjustment routines, as well as when you set or clear a preload.

If you use a label printer for GMP-compliant printouts and menu item 7.11.3 is active, the header and footer are printed on two different labels. To generate GMP-compliant printouts on labels, select menu item 7.11.2.

Three examples of GMP headers and one example of a footer are shown in the following. On Combics 1 models, the “date and time” line is not included.

Weighing platform WP1:	
-----	Dotted line
14.01.2002 09:43	Date/time ¹⁾
Model CW2P1-30ED-LCE	Combics model
Ser.no. 12345678	Combics serial no.
Vers. 1.1007.12.1	Software release for applications
BVers. 01-25-01	Software release for basic version
-----	Dotted line
Weighing platform WP2 (xBPI protocol): ²⁾	
-----	Dotted line
14.01.2002 09:45	Date/time ¹⁾
Model CW2P1-30ED-LCE	Combics model,
Ser.no. 12345678	Combics serial no.
Vers. 1.1007.12.1	Software release for applications
BVers. 01-25-01	Software release for basic version
Type IS12000S	Platform model
Ser.No. 12345678	Platform serial no.
-----	Dotted line
Weighing platform WP2 (SBI protocol): ²⁾	
-----	Dotted line
14.01.2002 09:45	Date/time ¹⁾
Model CW2P1-30ED-LCE	Combics model
Ser.no. 12345678	Combics serial no.
Vers. 1.1007.12.1	Software release for applications
BVers. 01-25-01	Software release for basic version
Type SBI	(Platform model)
-----	Dotted line
GMP footer:	
-----	Dotted line
14.01.2002 09:45	Date/time ¹⁾
Name :	Field for operator signature
-----	Blank line
-----	Dotted line

¹⁾ Not applicable for Combics 1 indicators
²⁾ Combics 2 indicators only

Sample Printouts

For details on the individual information blocks, see “Configuring Printouts” above. For details on configuring the header lines, refer to the chapter describing the particular application.

Weighing Application

There is no data for the “initialization data” block. If this block is enabled for the printout, a blank line is output.

```
      HEADER LINE 1
      HEADER LINE 2

14.01.2002      09:43
-----

G#      +      1.402 kg
T        +      0.200 kg
N        +      1.202 kg
-----
```

With weighing instrument serial number:

```
Ser.no.      80705337

G#      +      1.402 kg
T        +      0.200 kg
N        +      1.202 kg
-----
```

Counting Application

The “Initialization data” block contains the reference sample quantity and reference piece weight. The “Results” block contains gross, net and tare weights, as well as the calculated piece count.

```
-----
nRef      10 pcs
wRef      +      0.035 kg

G#      +      1.402 kg
T        +      0.212 kg
N        +      1.190 kg

Qnt      34 pcs
-----
```

Neutral Measurement Application

The “Initialization data” block contains the reference sample quantity and reference weight. The “Results” block contains gross, net and tare weights, as well as the calculated piece count.

```
Ref      2 o
wRef      +      1.200 kg

G#      +      14.700 kg
T        +      0.300 kg
N        +      14.400 kg

Qnt      12 o
-----
```

Weighing in Percent Application

The “Initialization data” block contains the reference percentage and reference weight. The results block shows gross, net and tare weights, as well as the percentage, which is shown as either the loss or the residual amount.

Percentage = residue:

```
-----
pRef      100 %
Wxx%      +      2.100 kg

G#      +      1.859 kg
T        +      0.200 kg
N        +      1.659 kg

Prc      79 %
-----
```

Percentage = loss:

```
-----
pRef      100 %
Wxx%      +      2.100 kg

G#      +      0.641 kg
T        +      0.200 kg
N        +      0.441 kg

D        21 %
-----
```

Checkweighing Application

The “Initialization data” block contains the nominal, minimum and maximum weights. The “Results” block always contains the gross, net and tare weights. The other results can be displayed in one of two ways:

- Result = Weight:
The deviation from the nominal weight is given both as a percentage and as an absolute (weight) value, whether the result lies within the tolerance limits or not.
- Result = Threshold status:
If the result lies within the tolerance limits, the printout shows the deviation from the nominal weight both as a percentage and as an absolute (weight) value, just as in the “Weight” printout mode described above.
If the result is outside the tolerance limits, the last line of the printout indicates the status as follows:

Result in “OK” range; “Weight” or “Threshold” printout:

```
-----
Setp      +      1.300 kg
Min        +      1.235 kg
Max        +      1.365 kg

G#      +      1.312 kg
T        +      0.000 kg
N        +      1.312 kg

Lim        +      0.92 %
W.Diff+    0.012 kg
-----
```

Result outside (over) the “OK” range; “Threshold” printout:

```
-----
Setp      +      1.300 kg
Min        +      1.235 kg
Max        +      1.365 kg

G#      +      1.400 kg
T        +      0.000 kg
N        +      1.400 kg

Stat      HH
-----
```

Data Interfaces

Classification Application

The "Initialization data" block contains the upper limits of Classes 1 through 4. The "Results" block contains gross, net and tare weights, as well as the class that the sample belongs to (1 through 5, where Class 5 means that the upper limit of Class 4 was exceeded).

Lim1	+	10.000	kg	
Lim2	+	11.000	kg	
Lim3	+	12.000	kg	
Lim4	+	13.000	kg	
G#	+	9.700	kg	
T	+	0.000	kg	
N	+	9.700	kg	
Class		1		

Animal Weighing Application

The "Initialization data" block contains the number of measured values that averaging is based on. The "Results" block contains the tare weight and the mean value.

mDef		8		
T	+	0.000	kg	
x-Net	+	4.202	kg	

Net-Total Formulation Application

The "Initialization data" block is empty. If this block is enabled for the printout, a blank line is output.

Which data is contained in the "Results" block value depends on the program operating status at the time of printing. The following options are available:

- Total/results printout
After **[CF]** is pressed
(tare memory is cleared)
- Individual/components printout
After **[OK]** is pressed
(component is stored in tare memory)
- Standard
After **[F]** is pressed (component is not stored in tare memory)

"Total" printout:

n		3		
S-Comp+		3.400	kg	
Cont.T+		0.200	kg	

Component printout (menu item 3.17.3)

When the components printout is configured, the header is printed only once, followed by all components.

If you are using a label printer, make sure a single label is large enough for the list of all components. For printer models YDP01IS and YDP04IS, you can configure manual form feed in the operating menu. If the corresponding setting is active, you can activate "form feed" manually. With the YDP02IS printer, form feed is automatic after each print command (fixed setting). Example with 2 components

HEADER LINE 1				
HEADER LINE 2				
14.01.2002		09:43		

Cmp001+		1.200	kg	
Cmp002+		2.000	kg	

Component printout

(menu item 3.17.2)

The entire standard printout is generated for each component.

Example for the second component:

HEADER LINE 1				
HEADER LINE 2				
14.01.2002		09:46		

Cmp002+		2.000	kg	

Standard printout

Example before the 2nd component is stored:

G#	+	4.400	kg	
T	+	0.200	kg	
T2	+	4.200	kg	
N	+	0.000	kg	

Totalizing Application

The "Initialization data" block is empty. If this block is enabled for the printout, a blank line is output.

Which data is contained in the "Results" block value depends on the program operating status at the time of printing. The following options are available:

- Printout of results
After **[CF]** is pressed
(totalizing memory is cleared)
- Individual/component printout of one transaction
After **[OK]** is pressed
(component is stored in tare memory)
- Standard printout
After **[F]** is pressed (component is not stored in tare memory)

"Total" printout:

*G		9.200	kg	
*N	+	8.600	kg	
n		3		

Component printout (menu item 3.17.3)

The header is printed only once, all transaction are printed one after the other. For printing on a label printer, see also "Component printout, Net-Total."

Example with 2 transactions:

HEADER LINE 1				
HEADER LINE 2				
14.01.2002		09:43		

G#	+	1.400	kg	
T	+	0.200	kg	
N	+	1.200	kg	
n		1		
G#	+	3.400	kg	
T	+	0.200	kg	
N	+	3.200	kg	
n		2		

Component printout
(menu item 3.17.2)

The entire standard printout is generated for each component.

Example: 2. Print second transaction:

```

      HEADER LINE 1
      HEADER LINE 2
14.01.2002      09:43
-----
G#      +      2.400 kg
T      +      0.200 kg
N      +      2.200 kg
n              2

```

Standard printout

The transaction counter value is not printed.

Example: 2. Print second transaction:

```

G#      +      2.400 kg
T      +      0.200 kg
N      +      2.200 kg

```

Print menu parameters:

All active menu item settings below the active menu level are printed.

```

-----
MENU
      SETUP.
WP1
-----
  1
    1.1
      1.1.2
      1.2.1
      1.3.2
    ...
  1.18
    1.18.1
      CAL
        10.000 kg

```

etc.

GMP-compliant Printouts

Linearization record:

```

-----
14.01.2002      13:00
Model CW2P1-30ED-LCE
Ser.no.      12345678
Vers.      1.1007.12.1
BVers.      01-25-01
-----

```

Linearization

```

Wt.1 +      7.00 kg
Wt.2 +     15.00 kg
Wt.3 +     22.00 kg
Wt.4 +     30.00 kg
      completed
-----

```

```

14.01.2002      13:02
Name:
-----

```

Calibration/adjustment record:

```

-----
14.01.2002      13:50
Model CW2P1-30ED-LCE
Ser.no.      12345678
Vers.      1.1007.12.1
BVers.      01-25-01
-----

```

```

External calibration
Nom. +     30.000 kg
Diff. -      0.003 kg
External adjustment
Diff. +      0.000 kg
-----

```

```

14.01.2002      13:52
Name:
-----

```

“Set preload” record:

```

-----
14.01.2002      13:50
Model CW2P1-30ED-LCE
Ser.no.      12345678
Vers.      1.1007.12.1
BVers.      01-25-01
-----

```

Set preload

completed

```

-----
14.01.2002      13:52
Name:
-----

```

“Clear preload” record:

```

-----
14.01.2002      13:50
Model CW2P1-30ED-LCE
Ser.no.      12345678
Vers.      1.1007.12.1
BVers.      01-25-01
-----

```

Clear preload

completed

```

-----
14.01.2002      13:52
Name:
-----

```

Weighing printout with multiple results;
example with 2 results:

```

-----
14.01.2002      09:43
Model CW2P1-30ED-LCE
Ser.no.      12345678
Vers.      1.1007.12.1
BVers.      01-25-01
-----

```

```

      HEADER LINE 1
      HEADER LINE 2
14.01.2002      09:43
-----

```

```

G#      +      2.40 kg
T      +      0.20 kg
N      +      2.20 kg
-----

```

```

      HEADER LINE 1
      HEADER LINE 2
14.01.2002      09:44
-----

```

```

G#      +      3.40 kg
T      +      0.30 kg
N      +      3.10 kg
-----

```

```

14.01.2002      09:45
Name:
-----

```

Error Codes

Error codes are shown on the main display. There are three types of error:

- Dynamic errors are indicated until the error is corrected, by an error code (e.g. *lnF 01*).
- Temporary errors are indicated for 2 seconds (e.g., *lnF 07*)
- "Fatal" errors are displayed continuously (e.g., *Err 101*) until the indicator is switched off and back on again.

Error Code	Cause	Solution
No display segments shown	No power connection	Check power supply
-----	Key pressed is not currently available	
<i>H</i>	Weighing capacity exceeded	Unload the weighing instrument
<i>L</i> or <i>Err 54</i>	Load plate/weighing pan is not on the weighing instrument	Place the load plate/weighing pan on the weighing instrument
<i>Err 101 - 104</i>	Key is stuck Key pressed when switching on the Combics	Release key or Contact your local Sartorius Service Center
<i>Err 320</i>	Operating program memory (EEPROM) defective	Contact your local Sartorius Service Center
<i>Err 335</i>	Verified platform not compatible with terminal	Connect a compatible weighing platform
<i>Err 340</i>	Operating parameter (EEPROM) error	Turn the weighing instrument off, then back on again. If this error code remains displayed, please contact your local Sartorius Service Center
<i>Err 341</i>	Data lost from RAM, battery needs to be recharged	Leave the weighing instrument power on for at least 10 hrs.
<i>Err 343</i>	Loss of data in the memory area for transaction numbers in external Alibi memory modules	Contact your local Sartorius Service Center
<i>lnF 01</i>	Data output not compatible with output format	Change the configuration in the Setup menu
<i>lnF 02</i>	Calibration/adjustment condition not met, e.g., – the weighing instrument was not tared – the weighing instrument is loaded	Calibrate only when zero is displayed Press  to tare Unload the weighing instrument
<i>lnF 03</i>	Calibration/adjustment could not be completed within a certain time	Allow the weighing instrument to warm up again and repeat the adjustment process
<i>lnF 06</i>	Internal calibration weight defective	Contact your local Sartorius Service Center
<i>lnF 07</i>	Function not allowed in weighing instruments verified for use in legal metrology	Contact your local Sartorius Service Center for information on having the settings changed
<i>lnF 08</i>	The load on the weighing instrument is too heavy to zero the readout	Check whether "tare/zero with power on" is set (1.12)
<i>lnF 09</i>	Taring is not possible when the gross weight is < zero	Zero the weighing instrument
<i>lnF 10</i>	Tare key is blocked when there is data in the tare memory	The data stored for the application program (Combics 2 only) must be deleted (clear the memory) before taring.
<i>lnF 18</i>	Preload is too light	
<i>lnF 19</i>	Preload is too heavy	
<i>lnF 29</i>	Minimum load not reached	Define a lower value for the minimum load (in the Application settings, under 3.6)
<i>lnF 30</i>	BPI ID (BPI byte) in active weighing instrument not deleted (cannot deactivate XBPI mode in COM1)	Restore factory settings for device parameters in the active weighing instrument
<i>lnF 31</i>	Interface handshake not completed	Send XON, CTS
<i>lnF 71</i>	Cannot store the current weight value (e.g., if control limits are too low or too high)	None
<i>lnF 72</i>	Cannot store the current weight value (e.g., the transaction counter has reached its limit)	None
<i>lnF 73</i>	Data not found or unreadable	Contact your local Sartorius Service Center
<i>lnF 74</i>	Function is blocked (e.g., menu is locked)	None
<i>no !IP</i>	No weighing platform connected	Contact your local Sartorius Service Center

Care and Maintenance

Service

Regular servicing by a Sartorius technician will extend the service life of your Combics indicator and ensure its continued weighing accuracy. Sartorius can offer you service contracts, with your choice of regular maintenance intervals ranging from 1 month to 2 years.

The optimum maintenance interval depends on the operating conditions at the place of installation and on the individual tolerance requirements.

Repairs

- ⚠ Disconnect defective equipment from power immediately (unplug the equipment from the wall outlet (mains supply)). Repairs may be performed only by authorized Sartorius service technicians using original Sartorius parts. Repairs performed by untrained persons may result in considerable hazards for the user. Important Note: If the equipment is still under warranty, send the entire indicator to the factory for repairs.

- ⚠ If a cable or cable gland is damaged or defective, replace the cable as a complete unit with all its connectors.

- ⚠ Do not open the indicator while it is carrying current. Allow at least 10 seconds to elapse after disconnecting the equipment from power before opening the equipment housing. Proper fitting of all surfaces is essential for the IP rating of the housing; for this reason the device must be opened and closed by a certified technician.

Cleaning

- ⚠ Disconnect the indicator from power (unplug the from the wall outlet (mains supply) and disconnect any data cables.
- ⚠ Make sure that no liquid enters the indicator housing.
- ⚠ Do not use any aggressive cleaning agents (solvents or similar agents).
- ⚠ Do not wash down the equipment with water or dry it with compressed air; this is not permitted.
- Clean the indicator using a piece of cloth which has been wet with a mild detergent (soap).
- If used in the food industry, use a cleaning agent suitable for the particular working environment.
- After cleaning, wipe down the indicator with a soft, dry cloth.

- If the weighing platform is installed in a pit, make sure that no dirt accumulates between the edge of the pit and the weighing platform. This will prevent measuring errors.
- Observe your company's internal regulations and standard industry guidelines with regard to cleaning intervals and cleaning agents.
- Regularly remove all dirt from the floor of the pit.

Cleaning Stainless Steel Surfaces

Clean all stainless steel parts regularly. Use a damp cloth or sponge to clean stainless steel parts on the weighing instrument. You can use any household cleaning agent that is suitable for use on stainless steel. Clean stainless steel surfaces by wiping them down. Then rinse the equipment thoroughly, making sure to remove all residues. Afterwards, allow the equipment to dry. If desired, you can apply oil to the cleaned surfaces as additional protection.

Solvents are permitted for use only on stainless steel parts.

Cleaning the Interior of the Weighing Platform

- If dirt enters the interior of the weighing platform, the load plate has to be removed. Particular caution is advised when working with models larger than 1000 × 1000 mm.
- Blow out the interior with pressurized air or rinse with a weak water jet (max. 60°C). Make sure no dirt enters the gap in the overload safety mechanism (the gap between the load cell and the screwed-on plate).

Corrosive Environment

- > Remove all traces of corrosive substances from the weighing platform on a regular basis.

Replacing the Dust Cover

- > Replace damaged dust covers.
- Place the new dust cover on the indicator and press down on the front and back along the edges until the cover is firmly seated.

Safety Inspection

Safe operation of the device is no longer ensured when:

- there is visible damage to the device or power cord
- the built-in power supply no longer functions properly
- the device has been stored for a relatively long period under unfavorable conditions (e.g., extreme moisture)

If there is any indication that safe operation of the device is no longer warranted:

- Disconnect the equipment from power (unplug the equipment from the wall outlet (mains supply) and lock it in a safe place to ensure that it cannot be used.
- Notify your nearest Sartorius Service Center or the International Technical Support Unit based in Goettingen, Germany.

Maintenance and repair work may be performed only by authorized Sartorius service technicians who:

- have access to the required service and maintenance manuals, and
- have attended the relevant service training courses.

- ⚠ The seals affixed to this equipment indicate that the equipment may be opened only by authorized service technicians, to ensure safe and trouble-free operation of the equipment and to maintain the conditions for warranty coverage.

Recycling

Sartorius products are packaged to ensure safe shipment using environmentally friendly materials. After successful installation of the indicator, you should return this packaging for recycling because it is a valuable source of secondary raw material.

For information on recycling options, including recycling of old weighing equipment, contact your municipal waste disposal center or local recycling depot.

If the equipment contains batteries, make sure to remove them before disposal. Batteries are hazardous waste and must be disposed of separately. Contact your municipal waste disposal center or local recycling depot for details on the proper disposal of batteries.

Overview

Common Specifications

Maximum readability	31,250 scale intervals (not in legal metrology)
Accuracy class	Ⓜ (on models ...-CE)
Verification scale intervals	≤3000e, (single-range scale) or 2 × 3000e (multiple range scale acc. to EN45501)
Digital protective interface	acc. to EN45501
Data interface	Bi-directional RS-232C interface with control outputs (standard equipment)
Additional data interface	optional
Display	20 mm LCD, 7-segment plus status symbols, backlit
Housing: Material	AISI 304 stainless steel
Dust and water protection acc. to EN60529	CW1P, CW1NP, CW2P: IP44 (optional IP65) CW1S, CW1NS, CW2S: IP67
Operating temperature range	-10°C to +40°C
Power supply	100–240 VAC (–15/+10 %), 50–60 Hz, max. 17 W/23 VA
DC Supply	optional 15.5–24 VDC (±10%), max. 12 W
AC Supply	optional 13–17 VAC (±10%), 50–60 Hz, max. 12 W
Battery operations	External rechargeable battery pack YBR10Z
Emissions	Acc. to EN61326+A1 Class B (IEC 61326+A1)
Immunity to interference	Acc. to EN61326+A1, industrial environment (IEC 61326+A1)
Electrical safety	Acc. EN61010-1 (IEC 1010-1), EN60950 (IEC 950)

Model-specific Specifications (Platform Specifications)

Model code:	CW...-L	CW...-I	CW...-LCE	CW....-NCE (2 x 3000e)			
	Readability	Readability	Readability	Weighing range 1		Weighing range 2	
Weighing capacity	15000d	30000d	1x3000e	Maximum capacity	Readability	Maximum capacity	Readability
3 kg	0.2 kg	0.1 g	1 g	1.5 kg	0.5 g	3 kg	1 g
6 kg	0.5 kg	0.2 g	2 g	3 kg	1 g	6 kg	2 g
15 kg	1 kg	0.5 g	5 g	6 kg	2 g	15 kg	5 g
30 kg	2 kg	1 g	10 g	15 kg	5 g	30 kg	10 g
60 kg	5 kg	2 g	20 g	30 kg	10 g	60 kg	20 g
150 kg	10 kg	5 g	50 g	60 kg	20 g	150 kg	50 g
300 kg	20 kg	10 g	100 g	150 kg	50 g	300 kg	100 g
600 kg	50 kg	20 g	200 g	300 kg	100 g	600 kg	200 g
1500 kg	100 kg	50 g	500 g	600 kg	200 g	1500 kg	500 g
3000 kg	200 kg	100 g	1000 g	1500 kg	500 g	3000 kg	1000 g

Type Designation

Here, we'll show you an example of how to put together order numbers. CW1P1-60 FE-LCE means the following:

Complete Combics scale with indicator 1
With one load cell
A maximum, single-range weighing capacity of 60 kg
A platform size of 500 x 400 mm
A resolution of 3,000 e for accuracy class III (CW1P...);
(...1...);
(...60...);
(...FE...); and
(...LCE)

Complete Combics scale

CW1P



Indicator with 20 mm LCD, backlit;
RS-232C interface port as a standard feature;
port for optional battery operation.
Indicator material: stainless steel
Type of protection: IP44

CW1NP



With selectable application programs.
Indicator with 20 mm LCD, backlit; integrated LEDs (red-green-yellow) for checkweighing or classification; RS-232C interface port as a standard feature; port for bar code scanner or optional battery pack.

Applications: weighing; counting; checkweighing; classification into 3 or 5 classes; totalizing; net-total formulation; filling; weighing in percent; neutral measurement; animal weighing.
Indicator material: stainless steel.
Type of protection: IP44.

CW2P



With selectable application programs.
Dot-matrix display, backlit. Integrated LEDs (red-green-yellow) for checkweighing or classification; RS-232C interface port as a standard feature; port for bar code scanner or optional battery pack.

Applications: weighing; counting; checkweighing; classification into 3 or 5 classes; totalizing; net-total formulation; filling; weighing in percent; neutral measurement; animal weighing.
Indicator material: stainless steel.
Type of protection: IP44.

Applications	Material Design	Number of load cells	-	Capacity (kg)	Platform size mm order code	-	Resolution
CW1	P	1	-	60	FE	-	LCE



1	3 kg	300×240	(DC)	L
1	6 kg	300×240	(DC)	15,000 d
1	15 kg	300×240	(DC)	
1	30 kg	400×300	(ED)	
1		500×400	(FE)	I
1	60 kg	400×300	(ED)	30,000 d
1		500×400	(FE)	
1		650×500	(GF)	
1		800×600	(IG)	LCE
1	150 kg	500×400	(FE)	3,000 e
1		650×500	(GF)	(verification scale intervals)
1		800×600	(IG)	
1	300 kg	650×500	(GF)	
1		800×600	(IG)	
4		1,000×1,000	(LL)	NCE
4		1,250×1,000	(NL)	2×3,000 e
4		1,500×1,250	(RN)	(verification scale intervals)
4		1,500×1,500	(RR)	
4		2,000×1,500	(WR)	
4	600 kg	1,000×1,000	(LL)	
4		1,250×1,000	(NL)	
4		1,500×1,250	(RN)	
4		1,500×1,500	(RR)	
4		2,000×1,500	(WR)	
4	1,500 kg	1,000×1,000	(LL)	
4		1,250×1,000	(NL)	
4		1,500×1,250	(RN)	
4		1,500×1,500	(RR)	
4		2,000×1,500	(WR)	
4	3,000 kg	1,000×1,000	(LL)	
4		1,250×1,000	(NL)	
4		1,500×1,250	(RN)	
4		1,500×1,500	(RR)	
4		2,000×1,500	(WR)	

CW2S4-1500RR-L, our example of a complete, stainless steel scale, means the following:

Complete Combics stainless steel scale with indicator 2 (CW2S...);
 With four load cells (...4...),
 A maximum, single-range weighing capacity of 1,500 kg (...1500...);
 A platform size of 1,500 × 1,500 mm (...RR...); and
 A resolution of 15,000 digits (...L)

Complete Combics stainless steel scale

CW1S



Indicator with 20 mm LCD, backlit;
 RS-232C interface port as a standard feature;
 port for optional battery operation.
 Indicator material: completely made of stainless steel Type of protection: IP67

CW1NS



With selectable application programs.
 Indicator with 20 mm LCD, backlit; integrated LEDs (red-green-yellow) for checkweighing or classification; RS-232C interface port as a standard feature; port for bar code scanner or optional battery pack.

Applications: weighing; counting; checkweighing; classification into 3 or 5 classes; totalizing; net-total formulation; filling; weighing in percent; neutral measurement; animal weighing. Indicator material: completely made of stainless steel. Type of protection: IP67.

CW2S



With selectable application programs.
 Dot-matrix display, backlit. Integrated LEDs (red-green-yellow) for checkweighing or classification; RS-232C interface port as a standard feature; port for bar code scanner or optional battery pack.

Applications: weighing; counting; checkweighing; classification into 3 or 5 classes; totalizing; net-total formulation; filling; weighing in percent; neutral measurement; animal weighing. Indicator material: completely made of stainless steel. Type of protection: IP67.

Applications	Material Design	Number of fond cells		Capacity (kg)	Platform size mm order code		Resolution
CW2	S	4	-	1500	RR	-	L




1	3 kg	300 × 240	(DC)	L
1	6 kg	300 × 240	(DC)	15,000 d
1	15 kg	300 × 240	(DC)	
1	30 kg	400 × 300	(ED)	
1		500 × 400	(FE)	I
1	60 kg	400 × 300	(ED)	30,000 d
1		500 × 400	(FE)	
1		650 × 500	(GF)	
1		800 × 600	(IG)	LCE
1	150 kg	500 × 400	(FE)	3,000 e (verification scale intervals)
4		650 × 500	(GF)	
4		800 × 600	(IG)	
4	300 kg	650 × 500	(GF)	
4		800 × 600	(IG)	
4		1,000 × 1,000	(LL)	
4		1,250 × 1,000	(NL)	NCE
4		1,500 × 1,250	(RN)	2 × 3,000 e (verification scale intervals)
4		1,500 × 1,500	(RR)	
4		2,000 × 1,500	(WR)	
4	600 kg	800 × 600	(IG)	
4		1,000 × 1,000	(LL)	
4		1,250 × 1,000	(NL)	
4		1,500 × 1,250	(RN)	
4		1,500 × 1,500	(RR)	
4		2,000 × 1,500	(WR)	
4	1,500 kg	1,000 × 1,000	(LL)	
4		1,250 × 1,000	(NL)	
4		1,500 × 1,250	(RN)	
4		1,500 × 1,500	(RR)	
4		2,000 × 1,500	(WR)	
4	3,000 kg	1,000 × 1,000	(LL)	
4		1,250 × 1,000	(NL)	
4		1,500 × 1,250	(RN)	
4		1,500 × 1,500	(RR)	
4		2,000 × 1,500	(WR)	

Overview


Dimensions (Scale Drawings)

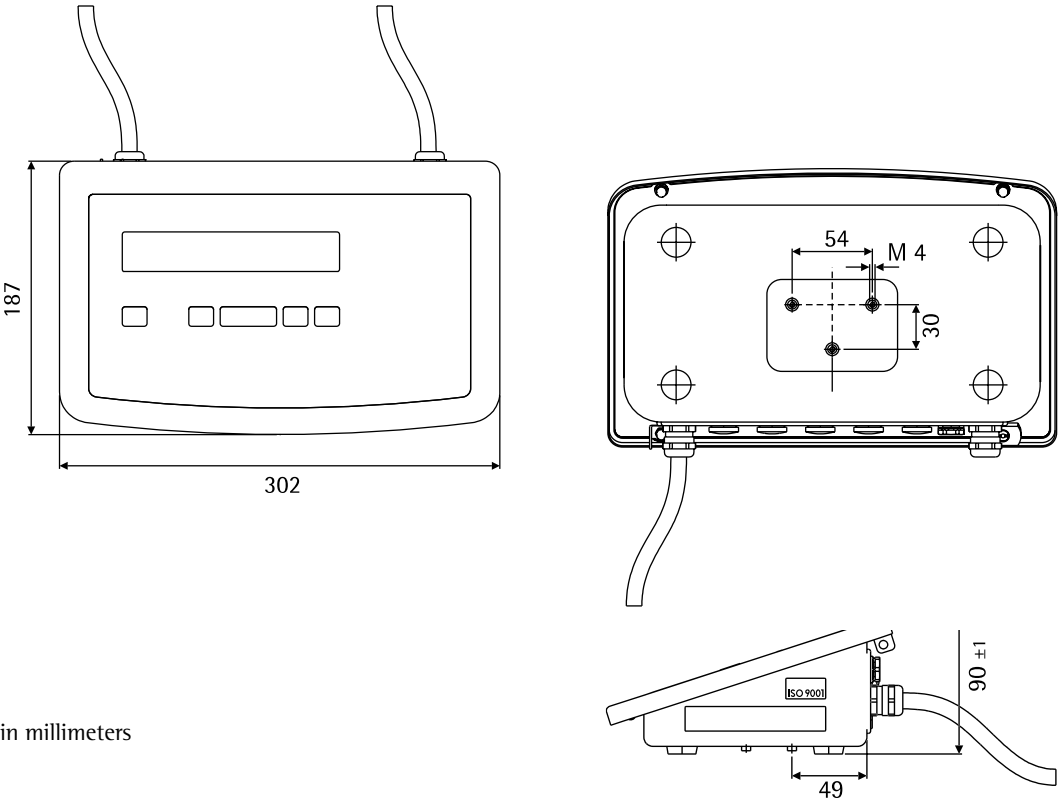Combics:



in millimeters

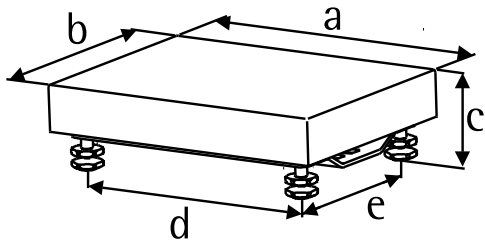Weighing Platforms:



Standard and Stainless Steel Versions			Height leveling feet	Distance between		Cable length
Length	Width Standard Version	Height Stainless Steel Version		d	e	
a	b	c	c	d	e	(m) ca.
(mm)	(mm)	(mm)	(mm)	(mm)	(mm)	
320	240	72	82	264	184	1.5
400	300	94	101	344	244	1.5
500	400	96	104	443	343	1.5 (painted) 3.0 (stainless steel)
650	500	145	90	530	434	3.0
800	600	145	90	680	534	3.0
1000	1000	90	90	804	804	6.0
1250	1000	90	90	1054	804	6.0
1500	1250	90	90	1304	1054	6.0
1500	1500	90	90	1304	1304	6.0
2000	1500	100	100	1804	1304	6.0


Accessories

Combiacs full-range scales:


Product	Order No.
Optional Interface	
Interface module (RS-232C) for UniCOM data interface	YD001C-232
Interface module (RS-485/422) for UniCOM data interface	YD001C-485
Profibus-DP module for UniCOM interface	YD001C-DP
Bluetooth® wireless interface module with antenna outside the housing	YD001C-BT
Analog current output, 0-20 mA, 4-20 mA, 0-10 V, 16-bit	YDA01C-20MA

Printers and Printer Accessories


Verifiable printer with functions for date, time and statistical evaluations	YDP03-OCE
Printer paper (5 rolls; length per roll: 50 m) for data printer	6906937
Replacement ink ribbon cartridge for printer	6906918

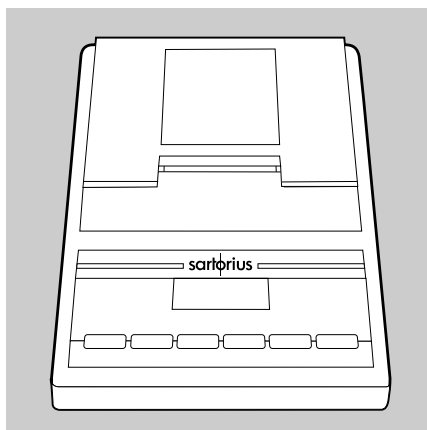

Verifiable strip and label printer with thermal print head, paper width 101 mm, with adapter cable (12-pin round male connector) and external power supply Adapter cable YCC01-01CISLM3 required for CW.P models. Adapter cable YCC02-R12F6 required for CW.S models.	YDP12IS-OCE-UV
Printer paper (1 roll) for YDP02IS -OCE printer, 101 mm × 75 m, thermo paper	69Y03196
Labels for YDP02IS-OCE printer, extra large, 101 × 127 m, 305 labels	69Y03195

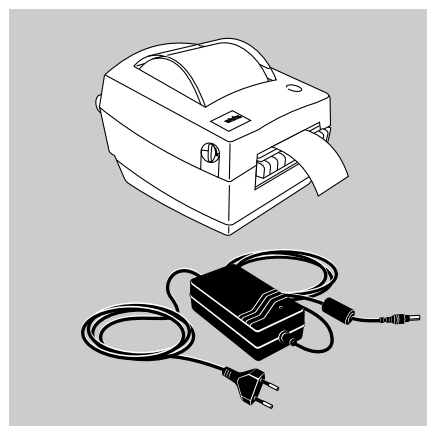

Verifiable strip and label printer with thermal print head, Paper width 60 mm, with adapter cable (12-contact round male connector) and external power supply. Adapter cable YCC01-01CISLM3 required for CW..P models. Adapter cable YCC02-R12F6 required for CW..S models.	YDP04IS-OCE-UV
Printer paper (3 rolls) for YDP02/04IS-OCE, 60 mm + 75 m, thermo paper	69Y03090
Labels for YDP02/04IS-OCE small, 58 + 30 mm, 1000 labels	69Y03092
Labels for YDP02/04IS-OCE medium, 58 + 76 mm, 500 labels	69Y03093
Labels for YDP02/04IS-OCE large, 58 + 100 mm, 350 labels	69Y03094
Cable for direct connection of YDP02/04IS-OCE to Combiacs CW.P models CISL indicator	YCC01-01CISLM3

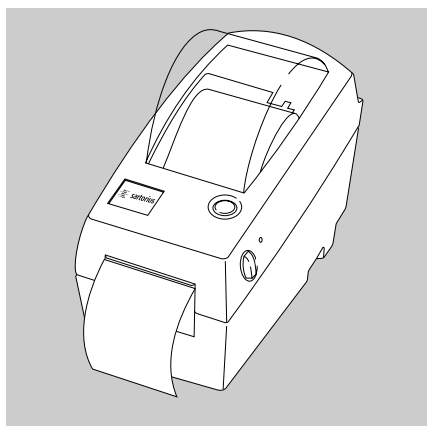

* not for use in legal metrology

Product	Order No.
Electrical Accessories	
External red/green/red display for CW.P models	YRD11Z
Second display* for CW.P models	YRD02Z
Remote display*, 7-segment, e. g. up to 45 mm characters	On request
Bar code scanner, 120 mm scanning width, with cable for connection to Combics CW.P indicators	YBR02CISL
Bar code scanner for CW.S models, adapter cable YCC02-R12F6	YBR02FC
Foot switch, incl. D-SUB 25-pin T-connector	YFS01
Hand switch, incl. D-SUB 25-pin T-connector	YHS02
External Alibi memory for electronic storage of weighing data	YAM01IS
Scanner for loading weighing data from PC, from YAM13IS	YAM02IS
Power supply for YAM01IS or YAM02IS	YAM11IS
Memory card for YAM01IS Alibi memory	YAM13IS
Cable for connecting Combics indicator to YAM01IS Alibi memory	YCC01-10CIM3
Cable for connecting (D-SUB 9-pin) YAM01IS Alibi memory to a PC (2 m)	69EM0012
Flow control for pumps with analog or digital pulse interface	YFC02Z-V2
Additional Accessories	
Dust covers (2 pcs)	YDC01CI
IP65 Kit for IP44 Combics	YAS01CISL
Anti-theft device	YTP01CI
Cable gland for Combics CW1S/CW4S models, IP67 protected	YAS04CIS
Mechanical Accessories	
Retainer plate for mounting indicator on platform (front-mounted indicator), made of AISI 304 stainless steel, for platform dimensions 240×300 mm	YDH11CWS
Installation kit for installing the Combics in a pit (with connecting hardware that lets you disconnect the indicator as desired)	YAS99I
Retainer for wall mounting, stainless steel	YDH02CIS
Floor-mounted column	YDH03CIP
Floor-mounted column, stainless steel	YDH03CIS
Base for installing floor-mounted column	YBP03CIP
Base for installing floor-mounted column, stainless steel	YBP03CIS
Retainer for bar code scanner, for attachment to floor-mounted column, bench stand or complete scale retainer	YBH01CWS
Plates for attaching printer to floor-mounted column or bench stand	YPP01CWS
Software	
Flexible formatting options for printouts (e.g., for bar codes with variable font size, graphics, etc.) – Ask your sales representative	Information available on request!
Sartorius WinScale scale driver software for Windows 95/98/2000/NT. Displays the scale readout on your PC monitor and provides secure memory for storing data that is subject to legal control.	
YCC01-09ISM5 RS-232 connecting cable required	YSW03
SartoConnect data transfer software for connecting your Sartorius scale to a computer running Windows 95/98/NT. This software lets you load the data recorded by your scale in a PC application program such as MS Excel or Access. Includes a cable for connecting the scale to a PC	YSC011

* not for use in legal metrology

Power Supply

24V industrial power supply module*	YAS02CI
External rechargeable battery pack, up to 40 h operation, incl. charger	YRB10Z
Connecting cable (25-pin D-SUB) for YRB10Z battery pack (2 m)	YCC02-RB01
Connecting cable with cable gland for YRB10Z battery pack (2 m)	YCC02-RB02
Connecting cable with cable gland for car battery (2 m)	YCC02-CB02

Connecting Cables

Connecting cable with cable gland for YBR02FC bar code scanner ¹⁾	YCC02-BR02
Connecting cable with cable gland 9-pin D-SUB male connector, 6 m ¹⁾	YCC02-D09M6
Connecting cable with cable gland, 9-pin D-SUB female connector, 6 m ¹⁾	YCC02-D09F6
Connecting cable with cable gland, 25-pin D-SUB female connector, 1.5 m ¹⁾	YCC02-D25F6
Connecting cable with cable gland, 12-pin round male connector, 6 m ¹⁾	YCC02-R12M6
Connecting cable with cable gland, 12-pin round female connector, 1.5 m ¹⁾	YCC02-R12F6
Cable for YDA01C-20MA power interface, with open cable ends e.g., 5 × 5 m*	6906926
Cable for connecting a PC; 25-contact D-Sub, approx. 1.5 m	7358312
Cable for connecting a PC; 9-contact D-Sub, approx. 1.5 m	7358314
Cable for connecting an isi terminal or QA/QC or FB/FC scale; 25-pin D-Sub to 12-pin round connector	YCC01-02ISM3
Cable for connecting a scale with 25-contact D-Sub female connector (25-pin D-Sub male connectors on both ends)	YCCD1-01M3
Cable for connecting an IS weighing platform; 25-pin D-Sub male connector to 25-contact D-Sub female connector, 3 m	YCC01.03CISLM3

* not for use in legal metrology

¹⁾ for stainless steel Combics models only (CW1S, CW1NS, CW2S)

Platform Accessories:

Dimensions Length × Width	800 × 600 mm	1000×1000 mm	1250×1000 mm	1500×1250 mm	1500×1500 mm	2000×1500 mm
Drive-on ramp, painted						
Order No.	YAR01CWP	YAR02CWP	YAR02CWP	YAR03CWP	YAR04CWP	YAR05CWP
Drive-on ramp, painted (tread plate)						
Order No.	YAR01CWPT	YAR02CWPT	YAR02CWPT	YAR03CWPT	YAR04CWPT	YAR05CWPT
Drive-on ramp, AISI 304 stainless steel						
Order No.	YAR01CWS	YAR02CWS	YAR02CWS	YAR03CWS	YAR04CWS	YAR05CWS
Drive-on ramp, AISI 304 stainless steel (tread plate)						
Order No.	YAR01CWST	YAR02CWST	YAR02CWST	YAR03CWST	YAR04CWST	YAR05CWST
Drive-on ramp, AISI 316 Ti stainless steel						
Order No.	YAR01CWS4	YAR02CWS4	YAR02CWS4	YAR03CWS4	YAR04CWS4	YAR05CWS4
Drive-on ramp, AISI 316 Ti stainless steel (tread plate)						
Order No.	YAR01CWST4	YAR02CWST4	YAR02CWST4	YAR03CWST4	YAR04CWST4	YAR05CWST4
Frame for pit installation, painted						
Order No.	YEG01CWP	YEG02CWP	YEG03CWP	YEG04CWP	YEG05CWP	YEG06CWP
Frame for pit installation, stainless steel						
Order No.	YEG01CWS	YEG02CWS	YEG03CWS	YEG04CWS	YEG05CWS	YEG06CWS
Length × Width	320 × 240 mm	400 × 300 mm	500 × 400 mm	650 × 500 mm	800 × 600 mm	
Roller conveyor, painted:						
Order No.	YRC01DCA	YRC01EDA	YRC01FEA	YRC01GFP	YRC01IGP	
Roller conveyor, AISI 304 stainless steel, for platform sizes:						
Order No.	YRC01DCS	YRC01EDS	YRC01FES	YRC01GFS	YRC01IGS	
Roller-ball load plate	–					
Order No.	–	YLP01CWS	YLP02CWS	YLP03CWS	YLP04CWS	
Column, painted, for attaching indicator to platform						
Order No.	YDH01CWP (Height 330 mm)	YDH02CWP (Height 500 mm) YDH03CWP (Height 750 mm) YDH03CWP (Height 750 mm)	YDH02CWP (Height 500 mm)	YDH03CWP (Height 750 mm)	–	
Column, AISI 304 stainless steel, for attaching indicator to platform						
Order No.	YDH01CWS (Height 330 mm)	YDH02CWS (Height 500 mm) YDH03CWS (Height 750 mm)	YDH02CWS (Height 500 mm)	–	–	
Bench, painted						
Order No.	–	YWT01CWP	YWT02CWP	YWT03CWP	YWT04CWP	
Bench, AISI 304 stainless steel						
Order No.	–	YWT01CWS	YWT02CWS	YWT03CWS	YWT04CWS	

Product	Order No.
Options	
Set of stainless steel floor fasteners (2 fastening plates, 4 special dowel screws), stainless steel	YFP01CWS
Column for bench, painted, for attaching indicator Adjustable height	YDH01WTCWP
Column for bench, stainless steel, for attaching indicator Adjustable height	YDH01WTCWS
Set of castors (2 guide castors, 2 lockable castors) for bench	YR001WTCW
Plate for bench column, for attaching indicator and printer	YPP01CWS
Retainer for bar code scanner, for attachment to bench column	YBH01CWS

Connecting an IS Weighing Platform to a Combics 2 Indicator

You can connect an IS weighing platform as WP2.

Characteristics:

- IS weighing platforms process weighing data independently of the indicator
- IS platforms can be internally adjusted
- ...-OCE models have a separate approval number, on a tag affixed to the cable
- Please observe the conditions described in the instruction manual for the platform connected

Declaration of Conformity

CE Marking on Sartorius Equipment

In 1985, the Council of the European Community approved a resolution concerning a new approach to the technical harmonization and standardization of national regulations. The organization for monitoring compliance with the directives and standards concerning the CE marking is governed in the individual EU Member States through the implementation of the EC Directives adopted by the respective national laws. As of December 1993, the scope of validity for all EC Directives has been extended to the Member States of the European Union and the Signatories of the Agreement on the European Economic Area.

Sartorius complies with the EC Directives and European Standards in order to supply its customers with weighing instruments and related equipment that feature the latest technology and provide many years of trouble-free service.

The CE marking may be affixed only to weighing instruments and associated equipment that comply with the following Directives:

Council Directive 89/336/EEC

"Electromagnetic compatibility (EMC)"

Applicable European Standards

1. Electromagnetic compatibility

1.1 Reference to 89/336/EEC:

Official Journal of the European Communities, No. 2001/C 105/03

EN 61326-1 Electrical equipment for measurement, control and laboratory use EMC requirements

Part 1: General requirements
Defined immunity to interference:
Industrial areas, continuous, un-monitored operation
Limitation of emissions:
Residential areas, Class B

Important Note:

The operator shall be responsible for any modifications to Sartorius equipment and for any connections of cables or equipment not supplied by Sartorius and must check and, if necessary, correct these modifications and connections. On request, Sartorius will provide information on the minimum operating specifications (in accordance with the Standards listed above for defined immunity to interference).

Council Directive 73/23/EEC "Electrical equipment designed for use within certain voltage limits"

Applicable European Standards:

EN 60950	Safety of information technology equipment including electrical business equipment
EN 61010	Safety requirements for electrical equipment for measurement, control and laboratory use
Part 1:	General requirements

If you use electrical equipment in installations and under ambient conditions requiring higher safety standards, you must comply with the provisions as specified in the applicable regulations for installation in your country.

Weighing Instruments for Use in

Legal Metrology:

Council Directive 90/384/EEC

"Non-automatic weighing instruments"

This Directive regulates the determination of mass in legal metrology.

For the respective Declaration of Type Conformity for Sartorius weighing instruments verified for use as legal measuring instruments that have an EC Type-approval Certificate, see the instruction manual(s) for each weighing instrument connected or refer to the enclosed "Guide to Verification."

This Directive also regulates the performance of the EC verification by the manufacturer, provided that an EC Type-approval Certificate has been issued and the manufacturer has been accredited by an officer of a Notified Body registered at the Commission of the European Community for performing such verification. Sartorius complies with EC Directive No. 90/384/EEC for non-automatic weighing instruments, which has been in effect since January 1, 1993, within the Single European Market, as well as the accreditation of the Quality Management System of Sartorius AG by Lower Saxony's Regional Administrative Department of Legal Metrology (Niedersächsische Landesverwaltungsamt – Eichwesen) from February 15, 1993.

For additional information on the CE mark on Sartorius equipment, see Sartorius Publication No. W--0052-e93081.

"EC Verification" - A Service Offered by Sartorius

Our service technicians authorized to perform the verification of your weighing instruments that are acceptable for legal metrological verification can inspect and verify the metrological specifications at the place of installation within the Member States of the European Union and the Signatories of the Agreement on the European Economic Area.

Subsequent Verifications within the European Countries

The validity of the verification will become void in accordance with the national regulations of the country in which the weighing instrument is used. For information on verification and legal regulations currently applicable in your country, and to obtain the names of the persons to contact, please contact your local Sartorius office, dealer or service center.

For more information on the verification of weighing instruments for use in legal metrology, contact the Sartorius Service Center.

Declaration of Type Conformity to Directive No. 90/384/EEC

This declaration is valid for non-automatic electromechanical weighing instruments for use in legal metrology. These weighing instruments accepted for legal metrological verification have an EC Type-Approval Certificate. The model(s) concerned is(are) listed below along with the respective type, accuracy class, and number of the EC Type-Approval Certificate:

Model	Type weighing instrument	Type indicator	Accuracy class	EC type-approval certificate no.	Indicator test certificate no.
CW...-CE	SARTICS	TN	III	D04-09-015	D09-03.13

SARTORIUS AG declares that its weighing instrument types comply with the requirements of the Council Directive on non-automatic weighing instruments, no. 90/384/EEC of 20 June 1990; the associated European Standard "Metrological aspects of non-automatic weighing instruments," No. EN 45501; the amended, currently valid versions of the national laws and decrees concerning legal metrology and verification in the Member States of the European Union, the EU, and the Signatories of the Agreement on the European Economic Area, which have adopted this Council Directive into their national laws; and with the requirements stipulated on the Type-Approval Certificate for verification. This Declaration of Type Conformity is valid only if the ID label on the weighing instrument has the CE mark of conformity and the green metrology

sticker with the stamped letter "M" (the two-digit number in large print stands for the year in which the mark has been affixed):



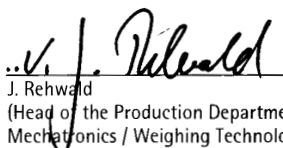
If these marks are not on the ID label, this Declaration of Type Conformity is not valid. Validity can be obtained, for example, by submitting the weighing instrument for final action to be taken by an authorized representative of SARTORIUS AG. The period of validity of this Declaration of Type Conformity shall expire upon any tampering with, repair or modification of this weighing instrument or, in some Member States, on the date of expiration.

The operator of this weighing instrument shall be responsible for obtaining an authorized renewal of the verification, such as subsequent or periodic verification, of the weighing instrument for use as a legal measuring instrument.

Sartorius AG
37070 Goettingen, Germany
Signed in Göttingen, 01.11.2004



Dr. G. Maaz
(President of the Mechatronics Division)



J. Rehwald
(Head of the Production Department
Mechatronics / Weighing Technology Division)

Physikalisch-Technische Bundesanstalt

Braunschweig und Berlin

PTB



EG-Bauartzulassung

EC type-approval certificate

Zulassungsinhaber:

Sartorius AG

Issued to:

Weender Landstr. 94-108
37075 Göttingen

Rechtsbezug:

In accordance with:

§ 13 des Gesetzes über das Mess- und Eichwesen (*verification act*)
vom/dated 23. März 1992 (BGBl. I S. 711), zuletzt geändert am (*last*
amended on) 25.11.2003 (BGBl. I S. 2304), in Verbindung mit Richtlinie
(*in connection with council directive*) 90/384/EWG, geändert durch (*amended*
by) 93/68/EWG

Bauart:

In respect of:

Nichtselbsttätige elektromechanische Waage mit oder ohne Hebelwerk
Nonautomatic electromechanical weighing instrument with or without
lever system

Typ / Type:

SARTICS

Max 0,5 kg ... 300 t

(III) $n \leq 6250$

(III) $n \leq 1000$

Option: Mehrbereichs- und Mehrteilungswaage
multi-interval and multiple range instrument

Zulassungsnummer:

Approval number:

D04-09-015 1. Revision

Gültig bis:

Valid until:

07.04.2014

Anzahl der Seiten:

Number of pages:

12

Geschäftszeichen:

Reference No.:

PTB-1.12-4014622

Benannte Stelle:

Notified Body:

0102

Im Auftrag

By order


Marcus Link



Braunschweig, 18.10.2004

Siegel
Seal

R3-0023

Die Hauptmerkmale, Zulassungsbedingungen und Auflagen sind in der Anlage enthalten, die Bestandteil der Revision der EG-Bauartzulassung ist. Hinweise und eine Rechtsbehelfsbelehrung befinden sich auf der ersten Seite der Anlage

The principal characteristics, approval conditions and special conditions, if any, are set out in the Annex which forms an integral part of this Revision of the EC type-approval certificate. For notes and information on legal remedies, see first page of the Annex.

Physikalisch-Technische Bundesanstalt

Braunschweig und Berlin

PTB



Prüfschein

Test certificate

Ausgestellt für:

Sartorius AG

Issued to:

Weender Landstr. 94-108
37075 Göttingen

Prüfgrundlage:

EN 45501 (1992), Nr. 8.1, WELMEC-Dokument 2.1 (1997) EWG
Richtlinie 90/384/EWG, OIML R 76-1

In accordance with:

Gegenstand:

Auswertegerät *Indicator*
und Terminal *and Terminal*
Typ: TN und / and TN-X

Object:

Typ:

Kennnummer:

Serial number:

Prüfscheinnummer:

D09-03.13 1. Revision

Test certificate number:

D09-03.13 Revision 1

Datum der Prüfung:

Date of Test:

Anzahl der Seiten:

9

Number of pages:

Geschäftszeichen:

PTB-1.12-4009190

Reference No.:

Benannte Stelle:

0102

Notified Body:

Im Auftrag

By order

Marcus Link



Braunschweig, 31.10.2003

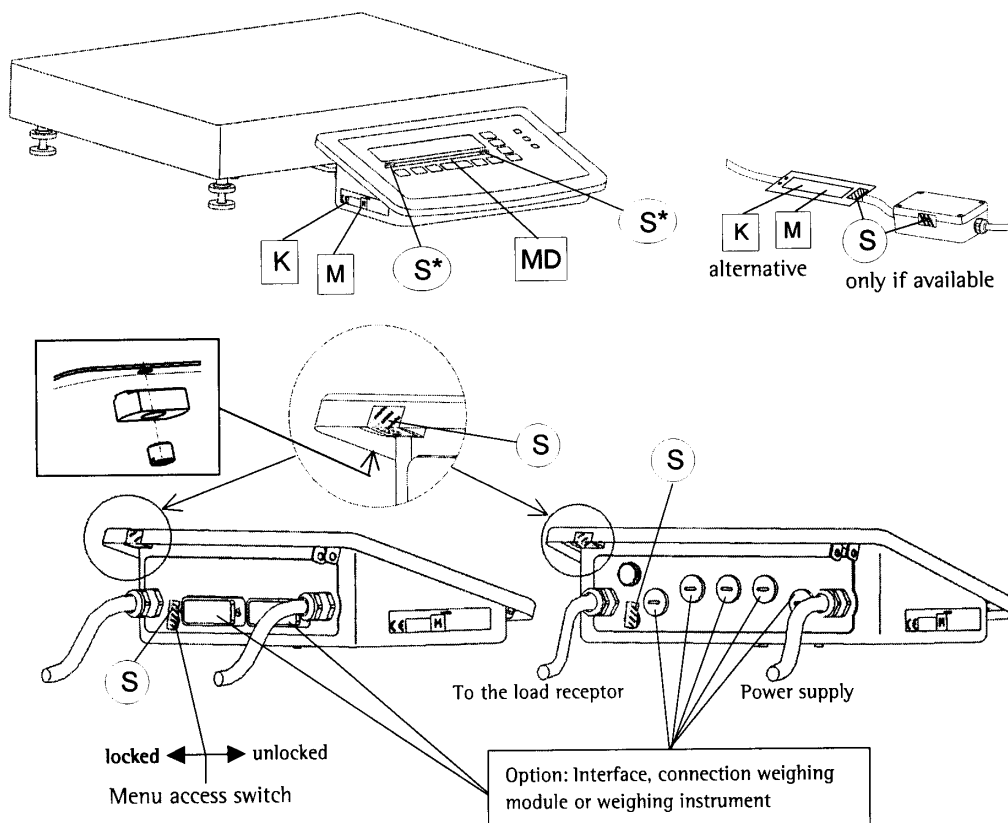
Siegel

Seal


R3-39200

Hinweise siehe erste Seite der Anlage, die Bestandteil des Prüfscheines ist.
For notes, see first page of the Annex which forms an integral part of the test certificate.


Plates and Markings



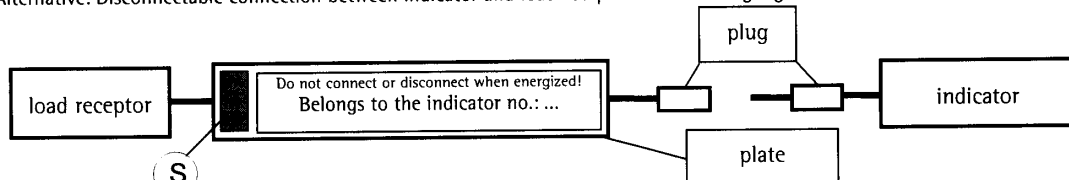Alternative: Disconnectable connection between indicator and load receptors with strain-gauge load cells:



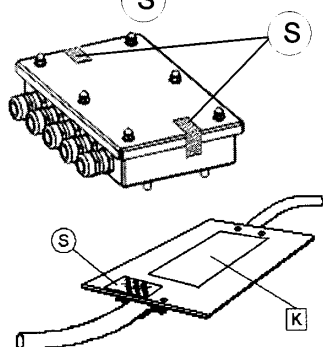If there is a junction box between load receptor and electronic evaluation unit the junction box has to be secured against inadmissible manipulation.




Alternative place for the Descriptive Plate of the weighing instrument
Handling in this case:

Affix the ID label of the weighing instrument to the delivered tag plate. Affix the ID tag plate to the data cable of the weighing module near the indicator. The verification officer or an authorized Sartorius representative must then seal the ID tag plate to the fastener.


- M** Mark for EC verification (green metrology sticker)
- K** Descriptive plate with CE-conformity
- MD** Metrological data: Max, Min and e

- S** Protective mark
- S*** Protective mark, only for detachable labels that remain intact

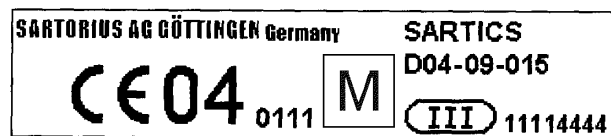
Example of descriptive plate of the already verified weighing instrument K



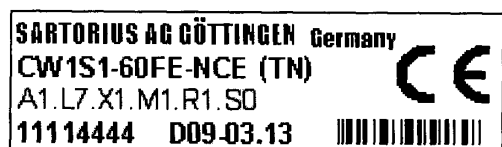
Example of plate with model designation of the terminal T



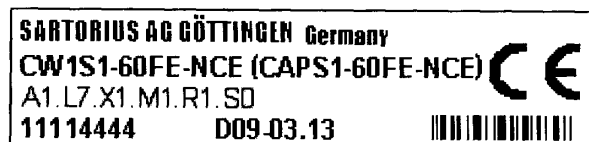Example of plate with model designation of the platform



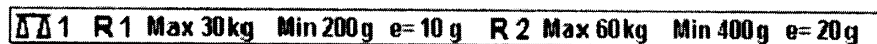Example of label with metrological data MD



In case of using as an indicating and operator terminal for weighing modules

- The same locations for protective marks are valid at the terminal
- A disconnectable connection between indicating and operator terminal and weighing module needs not to be secured with a protective mark
- For verification the regulations (Protective marks at the weighing module, fixing of the descriptive plate, etc.) of the weighing module have to be disregarded. To connect a weighing module the menu access switch of the Terminal must be set to the "unlocked" position. The serial number of the weighing module is stored in the Terminal and the Terminal is only working with this weighing module. After connecting the weighing module the menu access switch of the Terminal must be set to the "locked" position for verification. The weighing module has its own descriptive plate. The requirements of the approval no. D97-09-018 are valid.
- As well the metrological data of the weighing module have to be attached under the display and have to be secured against manipulation.

Example of label with metrological data for a connected weighing module



Index

	Page		Page
Access code	59	Net-total formulation	50
Accessories	98	Neutral measurement	29
Acoustic signal	14, 75		
Adjustment	23	Operating design	10
Animal weighing	47	Operating limits	9
Applications, consultation	2	Operating menu overview	57
Automatic data output	84	Operating menu parameters, overview	57
Averaging	47	Operation	13
Bar code scanner, connecting	83	Password	59
Battery, external rechargeable, connecting	83	Password, general	Appendix
		Peripherals, connecting	78
Cabling diagram	82	Pin assignment charts	79
Calibration/adjustment	23	Plates and markings	107
CE marking	103	Precautions	4
Checkweighing	33	Printer, configuring	70
Classification	38	Printing	87
Cleaning	93	Printout: samples	89
COM1 interface: settings	68	Printout: settings	74
Configuration	57		
Connecting cables to interfaces	81	Recycling	93
Connecting the equipment to AC power	7	Reference sample updating	26
Counting	25	Repairs	93
Data input format	84	Safety inspection	93
Data interfaces	77	Sample printouts	89
Data output format	85	SBI communication	84
Data output, configuring	87	Scale drawings	97
Data record: settings	68	Seal, verification	5
Data records: samples	89	Setting the time	76
Date and time, setting	76	Shipping conditions	5
Declaration of Conformity	103	Shutoff, automatic	14
Device information	76	Specifications	94
Device parameters	14, 66	Stainless steel surfaces, cleaning	93
Dimensions	97	Startup	5
Display elements	11	Storage and transport conditions	5
Display lighting	14, 75	Switch, remote (external)	74
EC type-approval certificate	105	Totalizing	53
Equipment supplied	5	Troubleshooting guide	92
General password	Appendix	Universal interface	71
General view of the equipment	6	Unpacking the equipment	5
Geographical data, checking	5		
GMP-compliant printouts	75, 88	Warmup time	7
		Warnings	4
ID codes	21	Weighing	13
Individual data ID codes	21	Weighing in percent	43
Installation instructions	5	Weighing platform WP1: settings	66
Installing the interface cable	81		
Intended use	2		
IP protection	4		
Key to model designations	95		
Keys, functions of	10, 12		
Language, setting	58, 76		
Legal metrology	5		
Leveling the weighing platform	9		
Maintenance	93		

General Password



When you select “Setup”, the password prompt is displayed for 2 seconds.



The first character of the password flashes.

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

;

Fn

repeatedly ,

→T←

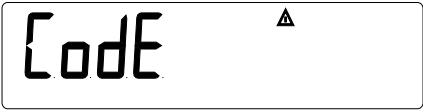
;

Entering the password:
– Press **Fn** to select numbers and press **→T←** to confirm. Press either **Fn** (to scroll through numbers from 0 up to 9) or **↵** (to scroll from 9 down to 0) repeatedly until the desired number is displayed.

The password is shown on the display.



Confirm the password entered.



Exit current menu level



Press and hold **→T←** (> 2 sec.)

Store input and exit the Setup menu



General password:
40414243

Service password:
202122

Sartorius AG
Weender Landstrasse 94–108
37075 Goettingen, Germany

Phone +49.551.308.0
Fax +49.551.308.3289
www.sartorius.com

Copyright by Sartorius AG,
Goettingen, Germany.
All rights reserved. No part
of this publication may
be reprinted or translated in
any form or by any means
without the prior written
permission of Sartorius AG.
The status of the information,
specifications and illustrations
in this manual is indicated
by the date given below.
Sartorius AG reserves the
right to make changes to
the technology, features,
specifications and design of
the equipment without notice.

Status:
December 2004, Sartorius AG,
Goettingen, Germany